



AMIS_Documents
Doc: ADOC_074

Originator:
Quality Specialist

Approver:
Managing
Director

Revision No: 006
Revision Date: 26.05.2019

Issued By: Management Rep

Certificate

AMIS0615

Certified Reference Material

Gold and Uranium Ore
Witwatersrand, South Africa

Certificate of Analysis

AMIS

A: 11 Avalon Road, West Lake View, Ext 11, Modderfontein, 1609, South Africa

P: PO Box 856, Isando, 1600, Gauteng, South Africa

T: +27 (0) 11 923 0800

W: www.amis.co.za



Table of Contents

Table of Contents	2
List of Tables	3
List of Appendices	3
Summary Statistics	4
1. Certified Concentrations and Uncertainties	6
2. Statistical Comparison of Means	8
3. Intended Use	8
4. Abbreviations and Symbols	9
5. Uncertified Concentration Values	10
6. Units	10
7. Analytical and Physical Methods	10
8. Origin of Material	10
9. Approximate Mineral and Chemical Composition	10
10. Quantitative Analysis by X-Ray Diffraction	11
11. Health and Safety	12
12. Method of Preparation	12
13. Particle Size Determination	12
14. Handling	13
15. Storage information	13
16. Methods of Analysis Requested	13
17. Information Requested of Participating Laboratories	13
18. Certification of Mean and Estimation of Measurement Uncertainty	14
19. Two Standard Deviations	14
20. Confidence Interval	14
21. Expanded Uncertainty	14
22. Confidence Interval and Expanded Uncertainty	14
23. Participating Laboratories	15
24. Accepted Assay Data	16
25. Reported Values	23
26. Validation of Accuracy (Trueness)	23
27. Limit of Detection and Limit of Quantitation in Gravimetric Fire Assay	23
28. Metrological Traceability	23
29. Period of Validity	23
30. Minimum Sample Size	23
31. Availability	23
32. Recommended use in Quality Control	23
33. Legal Notice	24
References	25
Appendices	26
End of certificate	35

List of Tables

Table 1. Certified concentrations, two standard deviations, combined and expanded uncertainty.	6
Table 2. Certified major oxides concentrations, two standard deviations, combined and expanded uncertainty.....	7
Table 3. The results of a two-sample equal or unequal variance t-test (two-tailed) data sets in which different analytical methods /instrumentation were used.....	8
Table 4. Abbreviations, symbols and descriptions.	9
Table 5. Results of XRD analysis.	12
Table 6. Particle Size Determination by laser diffraction.....	12
Table 7. Example of replicate assay data in which the <i>CI</i> , <i>2s</i> and <i>U</i> are shown.....	15
Table 8. Data used to calculate the certified values after removal of outliers.	16
Table 9. Uncertified element concentrations statistics.	26
Table 10. A single-factor ANOVA table showing key elements. Where <i>P</i> is the total number of groups, or laboratories. <i>P</i> -1 is 1 less than number of laboratories, <i>P</i> (<i>n</i> -1) is the number of data values minus number of groups (equals degrees of freedom for each group added together), and <i>P</i> -1 + <i>P</i> (<i>n</i> -1) is 1 less than the number of data points. <i>MS</i> is the mean squares of between laboratories and within laboratories. After Ellison <i>et al.</i> , (2009), Table 6.2, page 61.....	30
Table 11. CRM certified value quoted expanded uncertainty <i>U</i> , the coverage factor for the CRM, <i>k</i> =2.25 and mean for <i>n</i> =9 replicates and corresponding standard deviation for the replicate data.	31
Table 12. Mass of assay sample and corresponding limit of detection and limit of quantitation for an assay microbalance capability of smallest prill mass of 1µg or 0.001mg.	34
Table 13. Recommended reporting scheme for LOD and LOQ in fire assay.....	34
Table 14. T-distribution table for t-critical values (t crit.) for a two-tailed t-test at a 95% level of confidence.....	35

List of Appendices

Appendix 1: Uncertified Element Statistics.....	26
Appendix 2. Certification of Reference Material and Estimation of Measurement Uncertainty.....	27
Appendix 3. Example: Comparison of Mean and Certified Value for Validation of Accuracy	31
Appendix 4. Two-standard Deviations	32
Appendix 5. Confidence Interval.....	32
Appendix 6. Using the CRM in Quality Control.....	32
Appendix 7. Conversion to Air-dry Basis	33
Appendix 8. Example of Determination of LOD and LOQ in Fire Assay	34
Appendix 9. T-distribution table	35

Summary Statistics

Recommended Concentrations and Limits (at two Standard Deviations) Certified Concentrations

Analyte	Method	Certified (μ)⁷	(2s)⁹ \pm	Unit
Au	Pb Collection ¹	1.068	0.071	g/t
U	4A_MICP ²	15	2	ppm
C	Combustion/LECO ³	852	160	ppm
S	Combustion/LECO	0.44	0.04	%
S	4A_MICP	0.45	0.04	%
LOI	LOI ⁴	1.67	0.12	%
SG	SG ⁵	2.75	0.22	Dimensionless
Ag	4A_MICP	0.6	0.3	ppm
Al	4A_MICP	3.71	0.27	%
As	4A_MICP	70	8	ppm
Ba	4A_MICP	173	19	ppm
Be	4A_MICP	0.6	0.1	ppm
Bi	4A_MICP	0.8	0.09	ppm
Ca	4A_MICP	1080	134	ppm
Cd	4A_MICP	0.2	0.03	ppm
Ce	4A_MICP	36	4	ppm
Co	4A_MICP	24	2	ppm
Cs	4A_MICP	2	0.2	ppm
Cu	4A_MICP	42	4	ppm
Fe	4A_MICP	1.72	0.088	%
Ga	4A_MICP	9	0.5	ppm
Hf	4A_MICP	1	0.5	ppm
K	4A_MICP	9421	236	ppm
La	4A_MICP	19	2	ppm
Li	4A_MICP	23	1	ppm
Lu	4A_MICP	0.1	0.04	ppm
Mg	4A_MICP	4950	279	ppm
Mn	4A_MICP	180	13	ppm
Mo	4A_MICP	3	0.3	ppm
Na	4A_MICP	760	110	ppm
Nb	4A_MICP	2	0.8	ppm
Ni	4A_MICP	77	14	ppm
P	4A_MICP	90	16	ppm
Pb	4A_MICP	37	4	ppm
Rb	4A_MICP	37	1	ppm
Sb	4A_MICP	3	0.6	ppm
Sc	4A_MICP	6	0.4	ppm
Sn	4A_MICP	1	0.1	ppm
Sr	4A_MICP	17	2	ppm
Ta	4A_MICP	0.3	0.1	ppm
Tb	4A_MICP	0.2	0.08	ppm
Th	4A_MICP	5	0.5	ppm
Tl	4A_MICP	0.3	0.06	ppm
V	4A_MICP	35	4	ppm
W	4A_MICP	2	0.3	ppm
Y	4A_MICP	5	2	ppm
Yb	4A_MICP	0.6	0.3	ppm
Zn	4A_MICP	98	6	ppm
Zr	4A_MICP	51	12	ppm

Major Oxides
Certified Concentrations (at two Standard Deviations)

Analyte	Method	Certified (μ) ⁷	(2s) ⁹ \pm	Unit
U ₃ O ₈	4A MICP	0.002	0.0002	%
Al ₂ O ₃	XRF ⁶	7.34	0.055	%
CaO	XRF	0.14	0.01	%
Cr ₂ O ₃	XRF	0.040	0.002	%
Fe ₂ O ₃	XRF	2.45	0.057	%
K ₂ O	XRF	1.11	0.037	%
MgO	XRF	0.85	0.08	%
Na ₂ O	XRF	0.099	0.01	%
P ₂ O ₅	XRF	0.021	0.003	%
SiO ₂	XRF	85.76	1.2	%
SO ₃	XRF	1.07	0.094	%
TiO ₂	XRF	0.19	0.01	%

1. Certified Concentrations and Uncertainties

AMIS0615 is a new standard material, developed and certified in June 2019. Table 1 gives the certified concentrations, confidence interval, combined and expanded uncertainty for the certified reference material. Table 2 shows the certified major oxides concentrations, two standard deviations, confidence interval, combined and expanded uncertainty.

Table 1. Certified concentrations, two standard deviations, combined and expanded uncertainty.

Analyte	Method	Certified (μ) ⁷	N	n	k	% RSD	(u_c) ⁸	(2s) ⁹ \pm	(CI) ¹⁰ 95%	(U) ¹¹ \pm	Unit
Au	Pb Collection ¹	1.068	10	75	2.262	3	0.035	0.071	0.015	0.08	g/t
U	4A_MICP ²	15	5	40	2.776	6	0.8	2	1	2	ppm
C	Combustion/LECO ³	852	5	40	2.776	9	80	160	87	222	ppm
S	Combustion/LECO	0.44	6	47	2.571	5	0.02	0.04	0.02	0.05	%
S	4A_MICP	0.45	7	52	2.447	5	0.02	0.04	0.02	0.05	%
LOI	LOI ⁴	1.67	6	47	2.571	4	0.059	0.12	0.053	0.2	%
SG	SG ⁵	2.75	4	32	3.182	4	0.11	0.22	0.17	0.3	Dimensionless
Ag	4A_MICP	0.6	5	40	2.776	21	0.1	0.3	0.2	0.4	ppm
Al	4A_MICP	3.71	5	40	2.776	4	0.13	0.27	0.14	0.4	%
As	4A_MICP	70	7	54	2.447	6	4	8	4	10	ppm
Ba	4A_MICP	173	5	40	2.776	5	9	19	11	26	ppm
Be	4A_MICP	0.6	7	55	2.447	13	0.07	0.1	0.06	0.2	ppm
Bi	4A_MICP	0.8	5	40	2.776	5	0.04	0.09	0.04	0.1	ppm
Ca	4A_MICP	1080	7	56	2.447	6	67	134	55	163	ppm
Cd	4A_MICP	0.2	3	24	4.303	7	0.01	0.03	0.02	0.06	ppm
Ce	4A_MICP	36	5	39	2.776	5	2	4	2	5	ppm
Co	4A_MICP	24	6	48	2.571	4	1	2	0.8	2	ppm
Cs	4A_MICP	2	5	40	2.776	6	0.1	0.2	0.1	0.3	ppm
Cu	4A_MICP	42	6	48	2.571	4	2	4	2	5	ppm
Fe	4A_MICP	1.72	6	48	2.571	3	0.044	0.088	0.035	0.1	%
Ga	4A_MICP	9	5	40	2.776	3	0.3	0.5	0.3	0.8	ppm
Hf	4A_MICP	1	6	46	2.571	17	0.2	0.5	0.2	0.6	ppm
K	4A_MICP	9421	6	47	2.571	1	118	236	60	303	ppm
La	4A_MICP	19	5	39	2.776	5	0.9	2	1	2	ppm
Li	4A_MICP	23	5	40	2.776	3	0.6	1	0.5	2	ppm
Lu	4A_MICP	0.1	2	16	12.706	21	0.02	0.04	0.2	0.3	ppm
Mg	4A_MICP	4950	6	47	2.571	3	139	279	119	358	ppm
Mn	4A_MICP	180	5	40	2.776	4	6	13	7	18	ppm
Mo	4A_MICP	3	7	54	2.447	5	0.2	0.3	0.1	0.4	ppm
Na	4A_MICP	760	5	39	2.776	7	55	110	68	153	ppm
Nb	4A_MICP	2	4	32	3.182	19	0.4	0.8	0.6	1	ppm
Ni	4A_MICP	77	6	45	2.571	9	7	14	8	19	ppm
P	4A_MICP	90	5	37	2.776	9	8	16	8	22	ppm
Pb	4A_MICP	37	6	45	2.571	5	2	4	2	5	ppm
Rb	4A_MICP	37	5	38	2.776	2	0.7	1	0.5	2	ppm
Sb	4A_MICP	3	6	48	2.571	10	0.3	0.6	0.2	0.8	ppm
Sc	4A_MICP	6	6	45	2.571	4	0.2	0.4	0.2	0.5	ppm
Sn	4A_MICP	1	5	39	2.776	7	0.07	0.1	0.05	0.2	ppm
Sr	4A_MICP	17	6	47	2.571	5	0.8	2	0.8	2	ppm
Ta	4A_MICP	0.3	4	31	3.182	20	0.06	0.1	0.09	0.2	ppm
Tb	4A_MICP	0.2	3	24	4.303	18	0.04	0.08	0.1	0.2	ppm
Th	4A_MICP	5	5	40	2.776	5	0.3	0.5	0.3	0.7	ppm
Tl	4A_MICP	0.3	5	40	2.776	9	0.03	0.06	0.04	0.08	ppm
V	4A_MICP	35	6	48	2.571	5	2	4	2	5	ppm
W	4A_MICP	2	4	32	3.182	9	0.2	0.3	0.2	0.5	ppm
Y	4A_MICP	5	6	48	2.571	18	0.9	2	0.9	2	ppm
Yb	4A_MICP	0.6	3	24	4.303	21	0.1	0.3	0.3	0.6	ppm
Zn	4A_MICP	98	5	40	2.776	3	3	6	3	8	ppm
Zr	4A_MICP	51	5	40	2.776	12	6	12	7	16	ppm

Table 2. Certified major oxides concentrations, two standard deviations, combined and expanded uncertainty.

Analyte	Method	Certified (μ) ⁷	N	n	k	% RSD	(u_c) ⁸	(2s) ⁹ \pm	(CI) ¹⁰ 95%	(U) ¹¹ \pm	Unit
U ₃ O ₈	4A_MICP	0.002	5	40	2.776	6	0.0001	0.0002	0.0001	0.0003	%
Al ₂ O ₃	XRF ⁶	7.34	5	40	2.776	0.4	0.028	0.055	0.023	0.08	%
CaO	XRF	0.14	7	53	2.447	3	0.005	0.01	0.003	0.01	%
Cr ₂ O ₃	XRF	0.040	4	29	3.182	2	0.001	0.002	0.001	0.002	%
Fe ₂ O ₃	XRF	2.45	6	46	2.571	1	0.029	0.057	0.029	0.07	%
K ₂ O	XRF	1.11	6	44	2.571	2	0.019	0.037	0.023	0.05	%
MgO	XRF	0.85	4	30	3.182	5	0.04	0.08	0.07	0.1	%
Na ₂ O	XRF	0.099	4	32	3.182	6	0.006	0.01	0.007	0.02	%
P ₂ O ₅	XRF	0.021	5	39	2.776	7	0.001	0.003	0.002	0.004	%
SiO ₂	XRF	85.76	6	48	2.571	1	0.60	1.2	0.59	2	%
SO ₃	XRF	1.07	2	16	12.706	4	0.047	0.094	0.42	0.6	%
TiO ₂	XRF	0.19	6	48	2.571	4	0.007	0.01	0.006	0.02	%

1. Pb Collection
2. 4A_MICP is a Multi-acid digestion with either ICPOES/ICPMS/AAS finish
3. Combustion/LECO
4. LOI is Loss on Ignition
5. SG is Specific Gravity
6. XRF is X-ray Fluorescence
7. The certified value μ , is an unweighted grand mean of the means of N accepted sets of data from different laboratories and n number of test sample replicates. The certified value is traceable to SI units and is reported on a dry basis.
8. The combined uncertainty of the certified value is the within-laboratory reproducibility standard deviation derived from the analysis of variance of results from N number of laboratories and n number of sample replicates. (u_c)
9. Two standard deviations (2s)
10. Confidence interval at 95% level of confidence.
11. Expanded uncertainty (U) at a confidence level of 95% is determined by multiplication of the combined uncertainty (u_c) with a coverage factor (k) found from N-1 degrees of freedom (see Appendix 7 for t-distribution table). Example: $U = 2.36 \times 0.23 = 0.5\%$

2. Statistical Comparison of Means

A comparison of means for replicate data for the same element concentration determined by different analytical methods is done equating the variances between the two data sets; if the variances are found to be equal (F-test, p -value >0.05), then an equal variance t-test is applied. Should the variances be statistically significant, i.e. $p<0.05$, then an unequal variance t-test is performed. For either t-test, if the obtained p -value ≥ 0.05 , the null hypothesis that the means (certified values) are equal is accepted (Table 3). This gives the analyst confidence in the certified values reported by different analytical methods on the same analyte.

Table 3. The results of a two-sample equal or unequal variance t-test (two-tailed) data sets in which different analytical methods /instrumentation were used.

Method	Certified value	Method	Certified value	F-Test Outcome	p -value (t-test)	t-test Outcome
S Combustion/LECO	0.444 %	S 4A_MICP	0.453 %	Equal Variance ($p=0.309$)	0.485	Accept H_0 ; certified values are equal

3. Intended Use

AMIS0615 is a matrix matched Certified Reference Material, fit for use as a control sample in routine assay laboratory quality control when inserted within runs of test samples and measured in parallel to test samples. This material can also be used for method development, use as independent calibration verification check standard (i.e. if not used as a calibration standard in an instrument calibration), or for validation of accuracy in a method validation exercise (see Appendix 3). The recommend procedure for the use of this CRM as a control standard in laboratory quality control is to develop a Shewhart chart, where a mean value and corresponding 1, 2 and 3 standard deviations are derived from replicate measurements of the CRM (see Appendix 6). This CRM can also be used to assess inter-laboratory or instrument bias and establish within-laboratory precision and within-laboratory reproducibility. The certified concentrations and expanded uncertainty for this material are property values based on an inter-laboratory measurement campaign and reflect consensus results from the laboratories that took part in the exercise.

4. Abbreviations and Symbols

Abbreviations and symbols used in this document are shown in Table 4.

Table 4. Abbreviations, symbols and descriptions.

Abbreviation/Symbol	Description
Alpha (α)	Significance level (denoted by alpha, ' α ') of 0.05 or 5%
ANOVA	Analysis of variance by statistical means
Bq	The becquerel is the SI derived unit of radioactivity.
BIF	Banded iron formation
CRM	Certified reference material
df	Degrees of freedom, typically, $n-1$, or $N-1$
F_{calc}	Calculated F statistic from ANOVA or Fisher's test
F-critical or F_{crit}	F-critical value from F-distribution table
GOI	Gain on ignition
H_0	Null hypothesis
H_1	Alternate hypothesis
g/t	Grams per ton
k	Coverage factor, e.g. $k=2$ for 95% level of confidence
LOC	Level of confidence or confidence level
LOD	Limit of detection
LOQ	Limit of quantitation
LOI	Loss on ignition
MS	Mean squares (ANOVA)
MSb	Mean squares between (ANOVA)
MSw	Mean squares within (ANOVA)
N	Number of labs
n	Number of replicates
μ	Property or certified value of a CRM
p	' p -value' a measure of the strength of evidence against H_0
P	Total number of data points in ANOVA
ppm	Parts per million. Equivalent to g/t
RSD	Relative standard deviation usually expressed as % at a 68% LOC
Replicates	Replication is the repetition of an experimental condition so that the variability associated with an analysis can be estimated (ASTM E1847)
s	Standard deviation
s_r	Within laboratory repeatability as derived from ANOVA
s_s	Between laboratory standard deviation as derived from ANOVA
SS	Sum of squares in ANOVA
SST	Total variation in ANOVA
SSB	Between group (laboratory) variance
SSW	Within group (laboratory) variance

Abbreviation/Symbol	Description
2s	Two times standard deviation
SI	Standard International system of units
t_{calc}	Calculated t statistic from a one-sample, two-tailed t-test
t-critical or t_{crit}	t-critical value at given alpha and degrees of freedom
Tonne	A metric ton, is a unit of mass equaling 1000 kilograms
=TINV(5%, <i>df</i>)	MS Excel function for t-critical value at LOC 95% and <i>df</i>
<i>U</i>	Expanded uncertainty at a given k
<i>u</i>	Standard uncertainty at k=1
u_c	Combined standard uncertainty at $k=1$
μm	Micron, is an SI derived unit of length equaling 1×10^{-6} of a meter

5. Uncertified Concentration Values

Appendix 1 gives uncertified concentrations for other elements present in the CRM.

6. Units

All results for major oxides are reported as oxides in percentages. All results for major elements analyses reported in percentages or ppm. Results for Au and the platinum group elements are reported in g/t or ppm. Specific gravity (SG) is the ratio of the density of a substance to the density of a reference substance, *i.e.* equivalently; it is the ratio of the mass of a substance to the mass of a reference substance for the same given volume. Since specific gravity is a ratio of densities its units are therefore dimensionless.

7. Analytical and Physical Methods

A complete list of analytical and physical methods as generic method codes with a brief description of the methods is available on the AMIS web site www.amis.co.za

8. Origin of Material

The standard was made using material from the Ventersdorp Contact Reef of the Witwatersrand Basin.

9. Approximate Mineral and Chemical Composition

Throughout the Witwatersrand, gold ores occur in sheets or reefs originally deposited horizontally under water. The thickness of the reefs range between a lie of grit to several meters, with an average of 300mm. The sediments were also intersected by dykes and sills of dolomite, diabase and syenite intruding existing faults. The reefs can be in the form of either coarse conglomerates, rock pebbles and cemented in a silicate matrix. Pebbles usually derived from vein quartz, may also consist of quartzitic, chert jasper and quartz porphyry. The matrix consists of pure silica, but also contains minute flakes of muscovite and pyrophyllite as well as visible pyrite and other sulphides. The gold is usually confined to the matrix and tends to be concentrated along bedding planes (Vermeulen, 2001).

10. Quantitative Analysis by X-Ray Diffraction

Both natural and synthetic materials have a specific chemistry and atomic arrangement, known as phases. Phases can be identified and quantified using X-ray diffraction (XRD) which produces a plot of the intensity of X-rays scattered at different angles by crystalline phases in a material. Essentially, an X-ray diffraction pattern is the sum of the diffraction patterns produced by each phase. Simply put, an X-ray diffraction pattern is a fingerprint that allows the identification of what is in a target sample material. Knowledge of the mineral phase composition is useful in method development with techniques such as ICP-OES and XRF as potential matrix effects and spectral interferences can be recognised and accounted for. X-ray diffraction is effective in that it allows the identification of different phases of compounds that are identical in chemistry, but have a distinctly different the atoms, e.g. quartz, cristobalite, and glass are all different phases of SiO₂. Where quantitative XRD results do not correspond to results of other analytical techniques, it should be borne in mind that even though the data are quantitative they are meant to be used for indicative purposes in development of other analytical methods. Mineral names may not reflect the actual compositions of minerals identified, but rather the mineral group.

Quantification is determined from the chosen software package: this uses the full-profile Rietveld method of refining the profile of the calculated XRD pattern against the profile of the measured XRD pattern. The total calculated pattern is the sum of the calculated patterns of the individual phases. Results are given as weight % of the total crystalline phases and amorphous content. The amorphous content quantifies the amorphous material and unknown minerals or known minerals for which there is not a suitable crystal structure.

Corrections are incorporated into the process that allows for a more accurate description of the mineral's contribution to the measured pattern and to allow for variation due to atomic substitution, layer disordering, preferred orientation, and other factors that affect the acquisition of the XRD scan.

The limitations of qualitative XRD analysis are as follows:

- The detection of a phase may be dependent on its crystallinity.
- Where there exist multiple phases, overlap of diffracted reflections can occur, thus rendering some ambiguity into the interpretation.
- Overlapping reflections of a major phase can mask the presence of minor or trace phases.
- Some phases cannot be unambiguously identified as they are present in minor or trace amounts.

The limitations of quantitative XRD analysis by a full-profile Rietveld method are as follows:

- The limitations for qualitative XRD analysis apply.
- The method as described is standardless: it relies solely on the published crystallographic data available for each phase. Some data may not exactly describe the phases present.
- Particle size is important with respect to the absorption of the X-rays by the sample.
- Micronising reduces the particle size to that more suitable for quantitative analysis.

The accuracy of the analysis is dependent on sampling and sample preparation in addition to the calculated profiles being exactly representative of the chemistry of the component phases and their crystallinity. Some preferred orientation effects and reflection overlaps may occur which cannot be adequately resolved.

Table 5. Results of XRD analysis.

Phase	Formula	Composite %
Amorphous Content	—	11
Chlorite	(Fe,Al,Mg) ₆ (Si,Al) ₄ O ₁₀ (OH) ₈	5
Illite/Muscovite	K(Al,Mg,Fe) ₂ (Si,Al) ₄ O ₁₀ (OH) ₂	9
Pyrite	FeS ₂	1
Pyrophyllite	Al ₂ Si ₄ O ₁₀ (OH) ₂	2
Quartz	SiO ₂	71
Sodium plagioclase	NaAlSi ₃ O ₈	1
Total		100 %

For informational purposes only

11. Health and Safety

The material is a very fine powder coloured Pale brown (Corstor 5YR 6/2). Safety precautions for handling fine particulate matter are recommended, such as the use of safety glasses, breathing protection, gloves and a laboratory coat.

12. Method of Preparation

The particle size distribution for this material was shown to have a nominal top size of 54µm (95% passing 54µm). The procedure of preparation in brief is as follows: the material was crushed, dry-milled and air-classified to <54µm. It was then blended in a bi-conical mixer, systematically divided and sealed into 1kg Laboratory Packs. Explorer Packs are then subdivided from the Laboratory Packs as required. Final packaged units were then selected on a random basis and submitted for analysis to an independent laboratory accredited with the ISO17025 standard of general requirements for the competence of testing and calibration laboratories. The results obtained from this laboratory are then evaluated statistically by AMIS for homogeneity.

13. Particle Size Determination

The sample has been analysed using a Malvern Mastersizer 2000. Particles are passed through a focused laser beam that scatter light at an angle inversely proportional to their size. The intensity of light is measured and converted to a volume in particle size distribution. The results for this standard are presented in Table 8.

Table 6. Particle Size Determination by laser diffraction.

Size (µm)	Vol. Under %
<45um	96.1
<63um	97.8
<75um	98.3
<90um	98.8
<100um	99.1
<106um	99.2
<150um	99.9

For informational purposes only

14. Radioactivity

Shipments of this material require special labelling and placarding. AMIS0615 contains U (0.19 Bq/g) and Th 0.02 Bq/g) and is classified as EXEMPT MATERIAL in terms of "Safety Standards Series No. TS-R-1: Regulations for the Safe Transport of Radioactive Material, International Atomic Energy Agency, 2005, para 403, Table 1".

15. Handling

The material is packaged in Laboratory Packs and Explorer Packs that must be shaken or otherwise agitated before use. The analyte concentrations are quoted on a dry basis; therefore, the user needs to determine the moisture content to convert any obtained assay values to an air-dry basis (see Appendix 7 for an example calculation).

16. Storage information

The material should be stored in a cool dry place, in such a way that it does not compromise the integrity of the CRM. The material should be stored in conditions which will ensure it does not absorb moisture.

17. Methods of Analysis Requested

The following methods of analysis were requested:

- a) Au-Pb collection finished with either ICP-OES or ICP-MS or AAS or gravimetric
- b) Multi element scan to include all elements-4-acid total digestion including HF and/or peroxide fusion finished with either ICP-OES or ICP-MS or AAS
- c) LOI and all major oxides including U_3O_8 with XRF finish and/or Peroxide fusion and/or 4 acid digest including HF finished with either ICP-OES or ICP-MS or AAS (Please specify the temperature for LOI)
- d) SG – gas pycnometer
- e) S and C Combustion/LECO
- f) Moisture

18. Information Requested of Participating Laboratories

The following information was requested of the participating laboratories for the development of this CRM:

- a) State aliquots used for all determinations.
- b) All results for major elements to be reported as oxides in percentages.
- c) All results for multi-element scans and fusion to be reported in ppm.
- d) All results for Au and U to be reported in ppm.
- e) Report all QC data, to include replicates, blanks and certified reference materials used.
- f) All Round robin samples must be treated the same as routine test samples.
- g) All results must be reported to maximum decimal places i.e. dependent on laboratories capabilities
- h) Please ensure moisture content is determined and calculated. All results should be corrected by the moisture correction factor and this factor should be stated in the laboratory results. Please use the excel template provided by AMIS. If you require a copy, please email any of the email addresses below. Ensure all uncertainties are added to the results. Please send excel and PDF of all results. Ensure correct PPE is used i.e. gloves, dust masks and protective clothing. Analysis should be done under controlled environmental conditions. Ensure moisture content is determined and calculated. All results should be corrected by the moisture correction factor and this factor should be stated in the laboratory results. Use the excel template provided by AMIS. If you require a copy, please email any of the email addresses below. Ensure all uncertainties are added to the results. Send excel and PDF of all results.

19. Certification of Mean and Estimation of Measurement Uncertainty

The samples used in this certification process have been selected in such a way as to represent the entire batch of material and were taken from the final packaged units; therefore, all sources of uncertainty are included in the combined standard uncertainty determination. Initially the data submitted by all the laboratories are subjected to a z-score test, equation [1] to exclude outliers and the remaining data sets examined for their normality in distribution. This is followed by the exclusion of further outliers as defined by the IUPAC Harmonised Protocol of 1995 in which both Cochran and a Grubbs test are applied until all outliers are identified, equations [2] and [3]. A grand mean and standard deviation is recalculated using all remaining data (Thompson, 2008; Carr, 2011) (see Appendix 2)

20. Two Standard Deviations

AMIS reports two-standard deviations (2s) with all certified values. Two -standard deviations are calculated using the expression:

$$\text{Two standard deviations} = 2 (u_c)$$

Where u_c is the standard combined uncertainty (see Appendix, equation [14]).

21. Confidence Interval

AMIS reports a confidence interval (CI) with all certified values. Confidence interval as used by AMIS is:

$$\text{Confidence Interval (CI)} = \frac{(t_{critical})s}{\sqrt{N}}$$

Where, N is the number of laboratories (accepted laboratory data), $t_{critical}$ is a two-tailed value for $N - 1$ degrees of freedom (df) and s , is the standard deviation of the accepted laboratory means. A two-tailed critical value is found for $N - 1$ degrees of freedom from either a t -distribution table (Appendix 9) or MS Excel as =TINV (5%, df).

22. Expanded Uncertainty

ANOVA gives an estimate of the repeatability and the reproducibility of the data accepted for certification of the candidate reference material (see equations, [15] and [16], in the Appendix). Therefore, random variables (e.g., subsampling, instrument effects, interferences, operators and measurement conditions) that occur during the analysis of the candidate reference material by the various laboratories is considered. This approach does not necessarily quantify each individual source of uncertainty; however, the combined effect of random uncertainties is assessed (Ramsey & Ellison, 2007). A combined standard uncertainty is calculated from equation [14], which when multiplied by the t -critical value for $N - 1$ laboratories, gives an *expanded uncertainty* at a 95% level of confidence. The expanded uncertainty is a measure of the doubt around the certified value at a level of 95% confidence. The expanded uncertainty is used in the validation of accuracy (see equation [18]).

23. Confidence Interval and Expanded Uncertainty

A combined standard uncertainty will be greater than a combined CI . This is because ANOVA considers the within-lab repeatability (that is repeatability within each lab group) as well as the repeatability between each lab data set. This attends to random variables that contribute to the measurement of uncertainty, during the analysis of the test sample at the participating laboratories. The within-lab repeatability and the between lab repeatability is combined as the square root of the sum of squares of these two values giving a combined standard uncertainty, at a 68% confidence level. Multiplying the combined standard uncertainty by the t -critical value for $N - 1$, gives the expanded uncertainty at 95% level of confidence. It is recommended that the procedure described in Appendix 6, "Using the CRM in Quality Control" be used, in setting the limits of the CRM. Table 7 below shows mean gold values obtained by fire assay lead collection, for nine different laboratories, the confidence interval, two-standard deviations and expanded uncertainty.

Table 7. Example of replicate assay data in which the *CI*, *2s* and *U* are shown.

Lab No.	Mean Au (g/t)
1	0.268
2	0.273
3	0.270
4	0.288
5	0.274
6	0.256
7	0.263
8	0.258
9	0.288

<i>CI</i>	0.0088
<i>2s</i>	0.031
<i>U</i>	0.04

24. Participating Laboratories

The laboratories that are accredited with ISO17025 and provided timeous results are:

1. ALS Geochemistry Arabia
2. ALS Geochemistry Vancouver
3. Bureau Veritas Minerals Ultra Trace Pty Ltd
4. Chromatech Services
5. Intertek Mineral Tarkwa Ghana
6. Intertek Perth
7. SGS Geosol Laboratories Ltda (Brazil)
8. SGS South Africa
9. SGS Vancouver (Canada)
10. Shiva Analyticals India
11. Sibanye Gold SA

25. Accepted Assay Data

Data from the 11 laboratories used for certification are set out in Table 8.

Table 8. Data used to calculate the certified values after removal of outliers.

Pb Collection	Pb Collection	4A_MICP	4A_MICP	Combustion/LECO	Combustion/LECO	4A_MICP	LOI	SG
Au	Au	U	U ₃ O ₈	C	S	S	LOI	SG
g/t	g/t	ppm	%	ppm	%	%	%	Dimensionless
1.02	1.09	14.70	0.00	810	0.44	0.41	1.65	2.63
1.10	1.05	14.30	0.00	890	0.45	0.42	1.72	2.61
1.04	1.04	14.60	0.00	810	0.45	0.41	1.77	2.61
1.02	1.01	15.50	0.00	750	0.41	0.41	1.73	2.63
1.06	1.09	14.90	0.00	820	0.44	0.45	1.75	2.61
1.06	1.08	14.60	0.00	830	0.45	0.44	1.67	2.65
1.04	1.09	14.60	0.00	770	0.43	0.45	1.72	2.62
1.04	1.09	14.60	0.00	790	0.43	0.46	1.65	2.62
1.08	1.07	14.90	0.00	900	0.41	0.45	1.66	2.75
1.10	1.12	14.80	0.00	800	0.44	0.44	1.68	2.74
1.08	1.07	15.40	0.00	900	0.42	0.47	1.62	2.74
1.06	1.12	14.80	0.00	900	0.42	0.45	1.64	2.75
1.08	1.11	15.10	0.00	900	0.42	0.48	1.64	2.74
1.04	1.05	15.10	0.00	900	0.42	0.48	1.62	2.70
1.09	1.06	15.10	0.00	900	0.42	0.49	1.66	2.72
1.12	1.14	15.30	0.00	900	0.41	0.47	1.58	2.70
1.02	1.10	14.70	0.00	850	0.43	0.49	1.62	2.76
1.06	1.09	14.70	0.00	840	0.43	0.49	1.62	2.76
1.03	1.10	15.18	0.00	850	0.43	0.49	1.62	2.77
1.04	1.06	14.56	0.00	940	0.43	0.48	1.58	2.78
1.06	1.08	15.65	0.00	860	0.43	0.46	1.57	2.76
1.11	1.08	14.12	0.00	840	0.44	0.46	1.59	2.75
1.08	1.01	14.68	0.00	850	0.43	0.45	1.59	2.76
1.12	1.03	14.64	0.00	850	0.43	0.46	1.75	2.78
1.01	1.06	13.20	0.00	903	0.45	0.46	1.74	2.88
1.01	0.99	14.00	0.00	941	0.46	0.45	1.74	2.90
1.04	1.04	13.40	0.00	992	0.45	0.47	1.74	2.90
1.07	1.03	13.30	0.00	921	0.44	0.46	1.76	2.91
1.04	1.02	13.60	0.00	916	0.48	0.47	1.73	2.86
1.13	1.07	13.80	0.00	930	0.46	0.46	1.71	2.84
1.08	1.03	13.60	0.00	1020	0.45	0.49	1.77	2.89
1.11	1.04	13.90	0.00	926	0.46	0.47	1.76	2.88
1.09	1.09	15.75	0.00	800	0.47	0.47	1.62	
1.07	1.04	15.82	0.00	800	0.46	0.48	1.62	
1.06	1.09	15.52	0.00	800	0.46	0.47	1.71	
1.13		15.78	0.00	700	0.46	0.48	1.69	
1.14		15.68	0.00	800	0.45	0.43	1.67	
1.04		15.99	0.00	700	0.47	0.44	1.63	
1.10		15.64	0.00	700	0.45	0.46	1.63	
1.06		15.97	0.00	800	0.47	0.44	1.69	
					0.46	0.44	1.66	
					0.47	0.46	1.67	
					0.46	0.44	1.68	
					0.46	0.45	1.67	
					0.46	0.45	1.68	
					0.47	0.44	1.70	
					0.46	0.45	1.65	
						0.45		
						0.44		
						0.45		
						0.45		
						0.45		
						0.44		

Assay Data (Cont.)

XRF	XRF	XRF	XRF	XRF	XRF	XRF	XRF	XRF	XRF	XRF
Al ₂ O ₃	CaO	Cr ₂ O ₃	Fe ₂ O ₃	K ₂ O	MgO	Na ₂ O	P ₂ O ₅	SiO ₂	SO ₃	TiO ₂
%	%	%	%	%	%	%	%	%	%	%
7.31	0.14	0.04	2.48	1.10	0.83	0.10	0.02	86.14	1.04	0.19
7.35	0.15	0.04	2.47	1.12	0.84	0.10	0.02	85.49	1.05	0.19
7.30	0.14	0.04	2.49	1.12	0.83	0.10	0.02	86.19	1.03	0.19
7.31	0.15	0.04	2.48	1.12	0.83	0.09	0.02	86.16	1.04	0.19
7.33	0.15	0.04	2.49	1.12	0.83	0.09	0.02	85.52	1.06	0.19
7.31	0.15	0.04	2.46	1.12	0.82	0.10	0.02	85.93	1.05	0.19
7.36	0.15	0.04	2.49	1.12	0.83	0.10	0.02	85.70	1.03	0.19
7.36	0.15	0.04	2.43	1.10	0.83	0.09	0.02	85.42	1.03	0.19
7.32	0.15	0.04	2.44	1.13	0.84	0.11	0.02	85.12	1.11	0.19
7.39	0.15	0.04	2.43	1.13	0.84	0.10	0.02	85.71	1.11	0.19
7.30	0.15	0.04	2.43	1.13	0.84	0.10	0.02	85.08	1.09	0.19
7.37	0.15	0.04	2.43	1.11	0.84	0.09	0.02	85.46	1.09	0.19
7.33	0.14	0.04	2.42	1.14	0.84	0.10	0.02	85.34	1.11	0.18
7.36	0.15	0.04	2.45	1.12	0.85	0.10	0.02	85.23	1.11	0.18
7.33	0.14	0.04	2.44	1.13	0.84	0.10	0.02	85.97	1.11	0.19
7.29	0.14	0.04	2.49	1.14	0.84	0.11	0.02	85.81	1.11	0.19
7.33	0.15	0.04	2.48	1.12	0.92	0.10	0.02	85.99		0.19
7.33	0.15	0.04	2.48	1.12	0.91	0.09	0.02	85.85		0.20
7.34	0.14	0.04	2.49	1.12	0.92	0.09	0.02	85.64		0.20
7.33	0.14	0.04	2.49	1.13	0.92	0.10	0.02	85.97		0.19
7.32	0.14	0.04	2.50	1.13	0.91	0.09	0.02	86.05		0.19
7.34	0.15	0.04	2.48	1.13	0.91	0.10	0.02	86.22		0.20
7.32	0.14	0.04	2.48	1.12	0.82	0.09	0.02	85.50		0.20
7.33	0.14	0.04	2.42	1.12	0.81	0.10	0.02	85.61		0.20
7.38	0.14	0.04	2.41	1.12	0.81	0.11	0.02	84.90		0.21
7.39	0.14	0.04	2.41	1.12	0.81	0.10	0.02	85.00		0.20
7.34	0.14	0.04	2.41	1.12	0.82	0.10	0.02	85.00		0.20
7.36	0.14	0.04	2.41	1.12	0.81	0.10	0.02	85.00		0.20
7.36	0.14	0.04	2.40	1.12	0.81	0.10	0.02	84.90		0.20
7.35	0.14		2.42	1.12	0.81	0.11	0.02	85.00		0.21
7.33	0.14		2.47	1.12		0.11	0.02	84.80		0.20
7.35	0.14		2.43	1.12		0.10	0.02	84.90		0.20
7.35	0.14		2.45	1.07			0.02	86.50		0.19
7.37	0.14		2.46	1.07			0.02	86.60		0.19
7.36	0.14		2.46	1.07			0.02	86.60		0.19
7.37	0.14		2.46	1.07			0.02	86.60		0.19
7.39	0.14		2.46	1.13			0.02	86.70		0.20
7.37	0.14		2.46	1.12			0.02	86.70		0.20
7.37	0.14		2.45	1.13			0.02	86.70		0.19
7.38	0.15		2.44	1.12				86.80		0.19
	0.14		2.44	1.13				85.76		0.20
	0.14		2.44	1.13				85.80		0.20
	0.14		2.45	1.13				85.83		0.20
	0.14		2.44	1.13				85.79		0.20
	0.14		2.45					85.76		0.20
	0.14		2.45					85.89		0.20
	0.14							85.90		0.20
	0.14							85.83		0.20
	0.14									
	0.14									
	0.14									
	0.14									
	0.14									
	0.14									
	0.14									

Assay Data (Cont.)

4A MICP	4A MICP	4A MICP	4A MICP	4A MICP	4A MICP	4A MICP	4A MICP	4A MICP
Ag	Al	As	Ba	Be	Bi	Ca	Cd	Ce
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
0.5	36400	66	160	0.7	0.75	1000	0.19	37.56
0.49	34600	67	155	0.6	0.76	1000	0.18	36.24
0.46	34000	68	157	0.7	0.72	1000	0.2	36.08
0.48	36300	71	163	0.6	0.78	1000	0.2	37.65
0.51	36500	69	165	0.7	0.75	1000	0.21	36.9
0.41	34900	69	157	0.7	0.74	1000	0.19	36.68
0.45	35800	69	157	0.7	0.73	1100	0.21	36.48
0.55	35800	68	157	0.53	0.74	1000	0.17	37.16
0.79	38600	77	190	0.47	0.74	1100	0.21	38.8
0.66	38000	72	180	0.58	0.73	1100	0.2	36.5
0.71	38400	78.3	180	0.52	0.74	1100	0.22	38.6
0.67	39200	73.7	180	0.51	0.72	1100	0.22	37.9
0.67	38900	75.5	180	0.56	0.75	1100	0.21	38.3
0.64	37900	74.3	180	0.47	0.75	1100	0.2	37.2
0.67	40200	75	190	0.56	0.8	1200	0.21	38.6
0.78	38900	77.2	180	0.6	0.75	1100	0.22	39.8
0.72	37700	67	170	0.6	0.78	1000	0.21	33.92
0.68	37200	67	169	0.6	0.78	1000	0.2	33.68
0.76	38200	70	175	0.6	0.84	1000	0.19	33.76
0.7	36800	68	169	0.7	0.87	1000	0.21	35.37
0.68	37900	72	171	0.5	0.86	1000	0.19	33.42
0.79	37500	69	173	0.5	0.78	1000	0.21	34.07
0.75	37100	67	171	0.6	0.76	1000	0.2	33.06
0.8	37200	70	172	0.5	0.83	1000	0.2	38.1
0.5	36400	70	172	0.6	0.8	1200		38.1
0.5	37000	69	175	0.5	0.8	1200		38.4
0.5	37700	70	172	0.5	0.8	1200		37.6
0.5	36400	69	171	0.5	0.8	1200		37.5
0.5	36600	70	169	0.5	0.8	1200		37.4
0.5	37400	68	171	0.5	0.8	1200		37.2
0.5	36100	68	168	0.5	0.8	1100		37.6
0.5	36600	71	172	0.5	0.8	1100		35.2
0.8	37125	69	176	0.5	0.8	1100		36.4
0.7	37484	74	177	0.5	0.85	1100		35.7
0.7	36240	71	180	0.5	0.82	1100		35.6
0.8	37112	71	180	0.5	0.85	1100		34.3
0.7	36238	75	180	0.5	0.82	1100		36.3
0.8	37810	71	181	0.5	0.78	1100		36.3
0.7	37779	71	176	0.5	0.78	1100		34.5
0.7	37319	76	179	0.5	0.83	1100		
		65		0.5		1000		
		63		0.5		1100		
		65		0.5		1100		
		63		0.5		1100		
		65		0.5		1000		
		63		0.5		1100		
		72		0.5		1100		
		75		0.6		1100		
		74		0.6		1086		
		75		0.6		1085		
		74		0.6		1059		
		74		0.5		1087		
		73		0.6		1099		
		74		0.5		1097		
				0.6		1101		
						1082		

Assay Data (Cont.)

4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP
Co	Cs	Cu	Fe	Ga	Hf	K	La
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
22.6	1.89	41.5	16900	8.40	1.56	9500	19.60
22.6	1.77	41.5	16700	8.40	1.51	9200	19.10
22.5	1.83	41.4	16600	8.20	1.50	9400	18.90
23.5	1.82	42.8	17100	8.70	1.48	9600	19.90
22.9	1.91	42.1	17800	8.40	1.54	9600	19.50
23.4	1.87	41.8	16600	8.40	1.53	9300	19.50
22.9	1.85	42.8	17000	8.40	1.54	9500	19.20
22.9	1.89	41.1	16700	8.50	1.54	9500	19.50
25.5	2	42	17100	9.16	1.60	9500	19.80
23.6	2	39.5	16900	8.62	1.60	9600	18.50
24.8	2	41.4	17100	9.09	1.60	9300	19.60
24.1	2	41.3	17400	8.81	1.60	9500	19.40
24.8	2	42	17200	9.11	1.70	9500	19.60
24.8	2	40.2	16900	8.73	1.60	9500	20.00
24.6	2	41.5	17900	9.06	1.60	9400	20.10
24.8	2	41.8	17300	9.22	1.60	9300	17.60
24.2	1.9	39.7	17100	8.50	1.46	9300	17.90
23.4	1.8	39.5	16600	8.60	1.47	9300	18.00
24.1	1.8	39.8	17100	8.90	1.51	9300	17.70
25.3	1.8	39.2	16500	8.90	1.34	9400	18.20
24.9	1.8	39.6	16900	9.20	1.48	9400	18.10
23.8	1.8	38.7	16700	9.00	1.39	9400	18.30
23.8	1.8	38.9	16600	8.70	1.41	9300	17.30
25	1.8	39.6	16700	9.10	1.43	9400	18.50
24	1.74	44	17400	8.80	1.20	9400	18.70
24	1.78	43	17300	8.80	1.20	9700	18.80
25	1.8	45	18000	8.80	1.20	9300	18.40
22	1.69	44	17300	8.80	1.20	9400	18.20
24	1.66	45	17600	8.60	1.20	9600	18.30
25	1.76	44	17900	8.80	1.20	9300	18.00
24	1.77	43	17400	8.60	1.40	9400	18.40
24	1.73	44	17500	9.00	1.20	9200	17.70
24.5	1.8	40.7	17000	8.80	0.97	9400	18.70
25.4	1.9	41.7	17300	8.90	1.02	9500	18.00
25.4	1.7	42	17900	8.40	1.08	9300	18.10
24	1.8	41	17200	8.80	0.95	9200	17.00
23	1.8	41.2	17100	8.60	1.03	9500	18.00
24.4	1.8	40.7	17700	8.60	1.05	9500	18.70
25.1	1.8	39.5	17000	8.60	0.96	9300	17.40
24.3	1.9	41.1	17300	8.70	1.60	9443	
24.8		43	17600		1.50	9426	
24.9		43	17500		1.60	9547	
24.9		44	17200		1.60	9390	
24.9		43	17500		1.70	9450	
25.1		43	17700		1.50	9532	
25.3		44	18000		1.50	9550	
24.7		43	17800			9426	
25.8		44	17600				

Assay Data (Cont.)

4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP
Li	Lu	Mg	Mn	Mo	Na	Nb	Ni	P
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
24.20	0.09	4700	189	2.80	800	2.10	65.00	83.00
22.20	0.09	4700	184	2.71	800	2.30	64.70	83.00
23.30	0.09	4700	187	2.65	800	2.20	66.00	77.00
22.90	0.09	4900	190	2.65	800	2.10	68.00	93.00
23.60	0.08	5000	188	2.70	800	2.30	65.70	84.00
22.80	0.09	4700	183	2.77	800	2.20	64.20	70.00
23.00	0.09	4800	195	2.60	800	2.30	69.10	78.00
23.50	0.09	4700	189	2.64	800	2.20	66.00	86.00
23.00	0.12	4900	179	3.01	700	2.40	78.30	80.00
23.00	0.12	4800	176	2.94	700	2.30	73.10	90.00
23.00	0.12	4900	175	2.94	700	2.30	76.00	80.00
22.00	0.12	5000	176	3.02	700	2.30	75.90	90.00
23.00	0.12	4900	175	2.98	700	2.40	76.40	100.00
23.00	0.12	4800	170	3.36	700	2.30	74.60	90.00
23.00	0.12	5100	172	2.90	700	2.40	76.20	90.00
23.00	0.12	4900	174	2.97	700	2.30	76.70	90.00
22.00		5000	182	2.80	700	1.50	79.60	100.00
21.50		5000	172	2.75	700	1.50	78.30	100.00
21.50		5000	180	2.75	700	1.60	79.50	100.00
21.50		4900	172	2.69	700	1.50	78.00	100.00
22.00		4900	170	3.02	700	1.70	80.00	100.00
23.00		4900	172	2.70	700	1.40	79.60	90.00
23.00		4900	176	2.60	700	1.40	79.00	80.00
22.50		4900	178	2.87	700	1.50	78.70	90.00
23.00		5000	183	3.00	800	2.00	90.00	80.00
23.00		5000	179	3.00	800	1.90	90.00	90.00
23.00		4900	186	3.00	800	1.80	90.00	100.00
23.00		4900	184	3.00	800	2.00	90.00	80.00
22.00		5000	182	3.00	800	1.90	90.00	90.00
23.00		4900	187	3.00	800	2.10	73.00	93.00
22.00		4900	185	3.00	800	1.70	72.00	89.00
22.00		5000	183	3.00	800	1.80	78.00	88.00
22.70		5100	178	3.00	800		71.00	93.00
22.30		5100	178	3.00	800		75.00	87.00
22.40		5000	174	3.00	800		76.00	91.00
22.60		4800	178	3.00	800		72.00	93.00
23.00		5100	180	3.00	800		76.00	92.00
23.40		5000	180	3.00	800		77.00	
22.80		5000	181	3.04	800		76.00	
23.00		5111	180	2.99			76.00	
		5097		3.06			77.00	
		5005		2.85			75.00	
		5095		2.80			76.00	
		5165		3.02			78.00	
		5173		2.93			75.00	
		5187		2.93				
		5115		2.90				
				2.90				
				2.80				
				3.10				
				2.80				
				2.90				
				3.00				
				3.00				

Assay Data (Cont.)

4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP
Pb	Rb	Sb	Sc	Sn	Sr	Ta	Tb
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
34.00	37.60	2.74	5.80	1.10	18.40	0.36	0.23
34.70	36.60	3.13	5.60	1.00	18.00	0.37	0.20
33.60	36.60	3.31	5.60	1.10	18.20	0.34	0.21
35.00	38.10	2.97	6.10	1.00	18.50	0.32	0.22
33.60	37.10	3.07	5.80	1.00	18.00	0.35	0.21
33.40	36.70	3.23	5.80	1.00	18.30	0.33	0.21
33.00	36.80	3.96	5.80	1.00	18.30	0.34	0.21
34.30	37.00	3.13	5.90	1.00	17.30	0.31	0.21
36.00	36.30	2.97	6.30	1.00	16.00	0.34	0.19
34.50	38.90	3.18	5.90	1.00	16.80	0.34	0.19
35.90	37.80	3.35	6.10	0.90	16.20	0.32	0.19
34.60	38.00	3.15	6.00	1.00	17.00	0.34	0.18
36.00	37.10	3.93	6.20	1.00	16.60	0.33	0.19
36.40	38.10	3.53	6.00	0.90	16.60	0.34	0.19
36.30	38.40	3.66	6.10	0.90	16.90	0.34	0.19
36.60	36.70	3.53	6.20	0.90	16.90	0.22	0.18
37.20	36.60	2.76	5.90	1.00	16.70	0.22	0.26
36.40	37.10	2.80	5.90	1.00	16.90	0.25	0.26
37.00	36.10	2.82	5.50	1.00	16.70	0.23	0.26
35.20	37.80	3.03	5.80	1.00	16.80	0.25	0.28
37.60	36.60	3.11	6.10	1.00	16.80	0.20	0.28
36.00	36.70	2.71	5.60	1.00	16.70	0.20	0.26
36.10	37.60	3.08	5.70	1.00	16.70	0.22	0.26
36.20	38.00	2.85	6.00	1.20	16.00	0.26	0.26
37.00	37.80	3.00	6.00	1.10	16.00	0.28	
38.00	37.40	2.80	6.00	1.10	16.00	0.25	
38.00	37.40	3.00	6.00	1.10	15.50	0.29	
37.00	37.20	2.80	6.00	1.00	15.50	0.28	
36.00	37.00	2.80	6.00	1.00	15.50	0.28	
36.00	38.60	2.80	6.00	1.00	15.50	0.26	
36.00	36.80	2.80	6.00	1.00	16.00	0.30	
36.00	37.80	3.00	5.70	0.90	16.40		
38.00	37.20	3.28	5.50	1.10	17.00		
40.00	38.20	3.28	5.80	0.90	17.70		
37.00	38.20	3.48	5.50	1.00	17.30		
40.00	38.00	3.09	5.50	1.00	16.30		
39.00	36.80	2.90	5.70	1.10	17.50		
38.00	37.90	2.96	6.00	0.90	16.30		
40.00		2.70	6.00	1.10	16.80		
37.00		3.39	6.00		16.50		
39.00		3.00	6.00		16.90		
38.00		3.20	6.00		16.20		
39.00		3.20	6.00		16.90		
38.00		3.40	6.00		16.70		
39.00		3.30	6.00		16.70		
		3.20			16.40		
		3.40			17.10		
		3.30					

Assay Data (Cont.)

4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP
Th	Tl	V	W	Y	Yb	Zn	Zr
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
4.80	0.34	32.00	1.70	5.90	0.60	100.00	52.20
4.60	0.34	32.00	1.60	5.50	0.60	97.00	51.20
4.70	0.34	32.00	1.60	5.50	0.60	101.00	51.90
4.90	0.35	33.00	1.60	5.80	0.60	100.00	52.20
4.70	0.34	32.00	1.70	5.50	0.60	99.00	51.30
4.70	0.33	32.00	1.70	5.60	0.60	98.00	51.60
4.80	0.33	33.00	1.70	5.60	0.60	103.00	52.20
4.80	0.34	32.00	1.70	5.60	0.60	101.00	51.60
5.25	0.36	37.00	1.90	4.80	0.50	100.00	58.80
5.25	0.33	37.00	1.90	4.50	0.50	99.00	56.30
5.51	0.35	37.00	1.90	4.70	0.50	96.00	59.40
5.04	0.34	38.00	1.90	4.60	0.50	93.00	55.50
5.35	0.37	37.00	2.00	4.80	0.50	94.00	59.50
5.44	0.36	36.00	1.90	4.60	0.50	97.00	56.30
5.23	0.35	39.00	1.90	4.70	0.50	97.00	58.10
5.43	0.37	38.00	1.90	4.80	0.50	95.00	57.60
5.10	0.36	35.00	1.60	4.00	0.75	100.00	48.20
5.20	0.36	35.00	1.60	4.30	0.75	105.00	49.40
5.30	0.37	35.00	1.70	4.30	0.80	100.00	48.20
5.00	0.34	35.00	1.50	4.10	0.80	100.00	47.50
5.60	0.38	35.00	1.80	4.60	0.75	100.00	49.30
5.10	0.36	34.00	1.60	4.30	0.75	100.00	45.20
5.20	0.35	35.00	1.50	4.00	0.75	100.00	47.60
5.10	0.36	36.00	1.70	4.20	0.75	100.00	49.60
5.00	0.30	35.00	2.00	6.10		95.00	44.00
5.10	0.30	36.00	1.90	6.20		96.00	41.00
5.10	0.30	37.00	1.70	6.30		98.00	42.00
4.90	0.30	35.00	1.80	6.50		95.00	40.00
4.70	0.30	36.00	1.80	6.50		96.00	40.00
5.00	0.30	36.00	1.80	6.30		97.00	43.00
5.10	0.30	35.00	1.70	6.20		97.00	43.00
4.90	0.30	36.00	1.70	6.10		95.00	47.00
5.23	0.38	34.00		4.20		98.00	53.50
5.19	0.37	36.00		4.50		98.00	54.10
5.39	0.38	36.00		4.00		95.00	50.30
5.17	0.37	36.00		3.80		96.00	58.10
5.32	0.38	34.00		4.40		93.00	53.10
5.20	0.38	36.00		4.10		96.00	52.90
5.11	0.39	35.00		3.90		96.00	51.50
5.18	0.36	36.00		3.80		97.00	52.90
		38.00		4.90			
		37.00		4.50			
		37.00		4.50			
		37.00		4.50			
		36.00		4.60			
		38.00		4.90			
		37.00		4.60			
		36.00		4.60			

26. Reported Values

The certified values listed in this certificate fulfil the AMIS statistical criteria (see section 19) regarding agreement for certification and have been independently validated by Allan Fraser.

27. Validation of Accuracy (Trueness)

This CRM can be used to validate accuracy (trueness) as required in method validation as stated in the ISO17025 standard. See Appendix 3 for an example on the validation of accuracy using replicate data derived from the analysis of a CRM.

28. Limit of Detection and Limit of Quantitation in Gravimetric Fire Assay

In the determination of limit of detection (LOD) and limit of quantitation (LOQ) in gravimetric fire assay (*i.e.* lead collection and weighing of resulting gold prill), the minimum mass that an assay microbalance is capable of weighing and the original test sample mass determines the LOD and the LOQ in the assay (Fraser, 2015), (see Appendix 8 for an example of the calculation LOD and LOQ and Table 13 for a recommend reporting scheme for LOD and LOQ values).

29. Metrological Traceability

The values quoted herein are based on the consensus values derived from statistical analysis of the data from an inter-laboratory measurement program. Traceability to SI units is via the standards used by the individual laboratories the majority of which are accredited to the ISO17025 general requirements for the competence of testing and calibration laboratories and who have maintained measurement traceability during the analytical process.

30. Period of Validity

The certified values are valid for this product, while still sealed in its original packaging, until notification to the contrary. The stability of the material will be subject to continuous testing for the duration of the inventory. Should product stability become an issue, all customers will be notified and notification to that effect will be placed on the www.amis.co.za website.

31. Minimum Sample Size

Most of the laboratories reporting used a 0.5g sample size for the ICP-OES and a 30g sample size for the fire assay. These are the recommended minimum sample sizes for the use of this material.

32. Availability

This product is available in Laboratory Packs containing 1kg of material and Explorer Packs containing custom weights (from 50 to 250g) of material. The Laboratory Packs are sealed bottles delivered in sealed foil pouches. The Explorer Packs contain material in standard geochem envelopes, nitrogen flushed, and vacuum sealed in foil pouches.

33. Recommended use in Quality Control

Users should set their own limits *i.e.* 1, 2 and 3 standard deviations from an obtained mean value based on at least 10 replicate analyses using this CRM (see Appendix 6 for detail on the use of this CRM in quality control).

34. Legal Notice

This certificate and the reference material described in it have been prepared with due care and attention. However, AMIS, Makhosi Khoza, Daishnee Padayachee, Melesha Gopi Mungaroo and Allan Fraser; accept no liability for any decisions or actions taken following the use of the reference material.

Date of Version v1.00: 23 December 2020

Version: v1.00

Reason for Version v1.00: Amendment of Section 8 and 9- Origin of Material and Approximate Mineral and Chemical Composition

Version v1.00: replaces the original report of AMIS0615 Certification

Date of Version 000: 08 August 2019

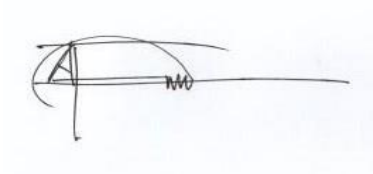
Version: 000

Approving Officer:

African Mineral Standards: _____

Melesha Gopi Mungaroo (Technical Manager)

Certifying Officer:

A handwritten signature in black ink, appearing to be 'AF', written over a horizontal line.

Geochemist: _____

Allan Fraser

M.Sc. (Geology), N.D. (Analytical Chem.),
Pr.Sci.Nat. Pr.Chem.SA

References

- Abzalov, M. (2011). Sampling Errors and Control of Assay Data Quality in Exploration and Mining Geology, Applications and Experiences of Quality Control. Ivanov, O., (Ed.), InTech, DOI: 10.5772/14965. Available from: <http://www.intechopen.com/books/applications-and-experiences-of-quality-control/sampling-errors-and-control-of-assay-data-quality-in-exploration-and-mining-geology>. Accessed 25th September, 2016.
- ASTM E122-09e1, Standard Practice for Calculating Sample Size to Estimate, With Specified Precision, the Average for a Characteristic of a Lot or Process, ASTM International, West Conshohocken, PA, 2011, www.astm.org
- Barwick V.J., Pritchard E. (Eds). (2011). Eurachem Guide; Terminology in Analytical Measurement – Introduction to VIM 3 (2011). ISBN 978-0-948926-29-7.
- Carr R.H. (2011). Estimating errors using graphs and taking good data. California State University Los Angeles. March 2011. 14. <http://web.calstatela.edu/faculty/kaniol/Err-Gph-Meas-IBooklet.pdf>. Retrieved: 29 August, 2016.
- Cochran, W.G. (1950). The Comparison of Percentages in Matched Samples. *Biometrika*, 37, 256-66.
- Ellison, S., Barwick, V., Duguid Farrant, T. (2009). *Practical statistics for the analytical scientist, a bench guide. 2nd Edition*. RSC Publishing. 25-172.
- ERM Application Note 1. (2005). Comparison of a measurement value with a certified value. European Commission – Joint Research Centre Institute for Reference Materials and Measurements (IRMM). 1-2.
- EURACHEM / CITAC Guide CG 4. (2012). Quantifying uncertainty in analytical measurement. 3rd Edition. Editors: Ellison, S. L. R., Williams, A. 4-121.
- Eurolab Technical Report No.1/2007. (2007). Measurement uncertainty revisited: Alternative approaches to uncertainty evaluation. European Federation of National Associations of Measurement, Testing and Laboratories. 43-46.
- Fraser, A.W. (2015). Minimising uncertainty in measurement and improving limit of detection in gold bearing ores from concentrations predicted by linear regression in atomic absorption spectrometry. M.Sc. thesis (unpublished), University of Johannesburg: 199 pages.
- Grubbs, F.E. (1969). Procedures for detecting outlying observations in samples. *Technometrics*, 11, 1969. 1-21.
- Horwitz, W., Albert, R. (2006). The Horwitz Ratio (HorRat): A useful index of method performance with respect to precision, *Journal of Association of Official Analytical Chemists International*, 89: 1095-1109.
- ISO/IEC 17025:2005(E). (2005). *General requirements for the competence of testing and calibration laboratories*. 2nd Edition: 14-59.
- ISO Guide 35 (2003). Certification of reference materials — General and Statistical principles, 3rd edition. ISO/REMCO WG 1.
- ISO 5725-2:1994. Accuracy (trueness and precision) of measurement methods and results -- Part 2: Basic method for the determination of repeatability and reproducibility of a standard measurement method.
- Long, J., Winefordner, J. (1983). Limit of detection - a closer look at the IUPAC definition. *Analytical Chemistry* 55: 712A - 724A.
- Miller, J., Miller, J. (2010). *Statistics for analytical chemistry*. 6th Edition. New York: Ellis Horwood. 36-126.
- Nelsen T.C., Wehling P. (2008). Collaborative studies for quantitative chemical analytical methods. AACC International Report. *Cereal Foods World*. September – October 2008, Vol. 53, No. 5. 285-288.
- Skoog, D., West, D. (1982). *Fundamentals of analytical chemistry*. 4th Edition. CBS College Publishing. Holt Saunders International Editions: 39-73.
- Thompson, M. (Ed.) (2008). Test for 'sufficient homogeneity' in a reference material. Analytical Methods Committee, AMCTB 17A, ISSN 1757-5958.
- Thompson, M. (Ed.) (2010). Internal quality control in routine analysis. AMC Technical Brief. Analytical Methods Committee. AMCTB No.46. 2010. Issn 1757-5958.
- Thompson, M., Lowthian, P. (2011). *Notes on statistics and data quality for analytical chemists*. Imperial College Press: 15-115.

Appendices

Appendix 1: Uncertified Element Statistics

Uncertified element statistics are shown in **Table 9**.

Table 9. Uncertified element concentrations statistics.

Element	Generic Method	n	Mean	SD	RSD %	Unit	Element	Generic Method	n	Mean	SD	RSD %	Unit
Al	FUS	8	38612.5	322.7	1	ppm	Nd	4A_MICP	15	14.3	0.4	3	ppm
Au	XRF	8	1.0	0.04	4	g/t	Nd	FUS	8	13.9	0.2	2	ppm
B	FUS	6	20.0	*	*	ppm	Pd	Pb Collection	15	0.004	*	24	g/t
BaO	XRF	8	0.01	0.01	38	%	Pr	4A_MICP	16	4.0	0.1	2	ppm
Cd	FUS	8	0.5	0.1	12	ppm	Pt	Pb Collection	10	0.01	*	34	g/t
Co	FUS	8	21.3	1.0	5	ppm	Rb	FUS	8	33.2	0.8	2	ppm
Cr	4A_MICP	32	219.2	65.3	30	ppm	Sb	FUS	8	2.3	0.1	4	ppm
Cr	FUS	8	312.5	23.1	7	ppm	Sc	FUS	7	7.0	0.00	*	ppm
Cs	FUS	8	1.8	0.1	5	ppm	Se	4A_MICP	15	1.9	1.0	50	ppm
Dy	4A_MICP	8	1.4	0.03	2	ppm	Si	FUS	8	387250.0	6227.8	2	ppm
Dy	FUS	8	1.8	0.1	4	ppm	Sm	4A_MICP	16	2.6	0.2	8	ppm
Er	4A_MICP	16	0.7	0.2	24	ppm	Sm	FUS	8	2.3	0.1	5	ppm
Er	FUS	8	0.9	0.1	7	ppm	SrO	XRF	1	0.01	*	*	%
Eu	4A_MICP	8	0.5	0.03	5	ppm	Tb	FUS	7	0.4	0.01	4	ppm
Fe	FUS	8	17437.5	244.6	1	ppm	Te	4A_MICP	13	0.1	0.01	23	ppm
Ga	FUS	8	7.0	*	*	ppm	Ti	4A_MICP	40	747.8	146.3	20	ppm
Gd	4A_MICP	7	2.0	*	*	ppm	Ti	FUS	8	1225.0	46.3	4	ppm
Ge	4A_MICP	16	0.5	0.4	87	ppm	Tm	4A_MICP	16	0.1	0.03	28	ppm
Ho	4A_MICP	8	0.3	0.01	3	ppm	Tm	FUS	8	0.1	0.01	10	ppm
In	4A_MICP	10	0.02	0.001	7	ppm	U	FUS	8	13.8	0.5	4	ppm
Mg	FUS	8	5037.5	106.1	2	ppm	U	XRF	16	16.9	5.2	31	ppm
MnO	XRF	40	0.02	0.01	23	%	U	XRF	16	16.9	5.2	31	ppm
Mo	FUS	8	4.0	*	*	ppm	U3O8	XRF	16	*	*	36	%
Moisture	Moisture	15	0.1	0.03	20	%	V2O5	XRF	2	0.02	0.01	47	%
Nb	FUS	8	4.0	*	*	ppm	Yb	FUS	8	0.8	0.1	9	ppm

* denotes that the results were too similar and SD and RSD% could not be calculated

Appendix 2 through 9, prepared by Allan Fraser.

Appendix 2. Certification of Reference Material and Estimation of Measurement Uncertainty

In the establishment of a consensus value for the CRM, outlier tests are carried out followed by performance statistics and the estimation of the measurement uncertainty. In practice, it is highly likely that data generated by multiple laboratories as an inter-laboratory comparison of material for certification, will contain erroneous as well as extreme measurements (outliers). The influence of outliers on summary statistics needs to be minimised by the application of procedures for outlier identification on raw data. The use of z-scoring, Cochran's test for suspect repeatability variances, along with Grubbs test for suspect measurement values allows for the detection of outliers (IUPAC, 1995). Method performance in terms of precision as relative standard deviation is judged by the application of the Horwitz ratio, which gives an indication of whether the observed relative standard deviation at the concentration levels of analyte determined are acceptable (Horwitz & Albert, 2006).

In the absence of an extensive uncertainty budget, measurement uncertainty is estimated from the reproducibility standard deviation from inter-laboratory data and reported as an expanded uncertainty at a level of confidence of 95% (Miller & Miller, 2010).

The steps below give detail on the establishment of a consensus value through the elimination of outliers, method performance and estimation of measurement uncertainty using standard uncertainties and the analysis of variance.

Z-Score

A z-score is calculated using equation [1]:

$$z = \frac{x - x_a}{s_p} \quad [1]$$

Where, x is the result of a submitted sample, x_a is the mean and s_p is the standard deviation of the submitted results from all the participating laboratories. Z-Scores are interpreted as follows:

$|z| \leq 2$ satisfactory performance
 $2 < |z| \leq 3$ questionable performance
 $|z| > 3$ unsatisfactory performance

(Thompson & Lowthian, 2011)

Data with z-scores exceeding two are discarded and are not included for further assessment.

Cochran's Test

The test of Cochran (1950) as shown in equation [2] is applied to any suspect repeatability variances:

$$C_{calc} = \frac{s_{max}^2}{\sum_{i=1}^l s_i^2} \quad [2]$$

Where, C_{calc} , s_{max}^2 and $\sum_{i=1}^l s_i^2$, are the calculated values for Cochran's test, data set with the maximum variance and the sum of the variances of all of the participating, l laboratory datasets. The C_{calc} value is compared with a critical value, C_{crit} at a level of confidence of 95% and an alpha of 0.05% (see Ellison, *et al.*, 2009, Appendix A, Table A.3a, page 209 for a table of critical values for the test of Cochran at LOC 95%).

According to ISO 5725-2 (1999), results from a laboratory with a suspect repeatability variance can be excluded if it is shown by the Cochran test to be an outlier. Therefore, if $C_{calc} > C_{crit}$, the laboratory with the maximum variance is removed. The data found to be excluded should not be $>2/9$, or 22% of the total data.

Grubbs Test

The test of Grubbs (1969) calculates a test statistic, G_{calc} and in the detection of a single outlier, G_1 is found by using

$$G_{1\,calc} = \frac{|Suspect\ value - \bar{x}|}{s} \quad [3]$$

Where, the sample mean and standard deviation, \bar{x} and s , are calculated with the suspect value included. The $G_{1\,calc}$ statistic is compared to a critical value for N measurements. See Ellison, *et al.*, 2009, Appendix A, Table A.2, page 208 for a table of critical values for the test of Cochran at LOC 95%.

Method Performance

The Horwitz function is used to assess the performance of the data under consideration, with respect to precision (Horwitz & Albert, 2006). A calculated %RSD is found using the Horwitz expression

$$\%RSD = \pm 2^{(1-0.5\log C)} \quad [4]$$

where, C is the analyte concentration in percent divided by 100 and \log is the natural logarithm. The observed %RSD is calculated as

$$Observed\ \%RSD = \frac{s}{Mean} \times 100 \quad [5]$$

where s is the standard deviation of n replicates.

The ratio of the observed %RSD and the calculated %RSD gives the Horwitz ratio (HorRat):

$$HorRat = \frac{\%RSD\ Observed}{\%RSD\ Calculated} \quad [6]$$

A HorRat <2 indicates that the method is of adequate precision. Should the HorRat be >2 the overall data are discarded, and the candidate material considered not suitable for certification as the precision is excessive for the concentration of the analyte being determined (Nelsen & Wehling, 2008).

Grand Mean

The grand mean ($\bar{\bar{x}}$) *i.e.* the certified value of a dataset is the total of all the data values divided by the total sample size (n):

$$\bar{\bar{x}} = \sum \frac{x}{n} \quad [7]$$

Certified Value

From ANOVA as per the description in section 19, an 'appropriate precision' as shown in [8] is calculated for sufficient homogeneity (Thompson, 2008):

$$s_r \leq 0.3u_c \quad [8]$$

Where, s_r is the within laboratory repeatability, as determined from [14]. Once [8] is satisfied, a grand mean [7] is calculated and this is taken to be the certified value.

Total Variation (SST)

The total variation (not the variance) comprises the sum of the squares of the differences of each mean with the grand mean.

$$SST = \sum (x - \bar{x})^2 \quad [9]$$

Between Group Variation (SSB)

The *variation* due to the between the laboratories is denoted SSB or Sum of Squares Between laboratories and given by [10]. If the laboratory means are close to each other (and therefore the Grand Mean) SSB will be a small value. There are P samples involved with one datum value for each sample (the sample mean), so there are P-1 degrees of freedom.

$$SSB = \sum n(\bar{x} - \bar{\bar{x}})^2 \quad [10]$$

The *variance* due to the interaction between the laboratories is denoted MSB for Mean Square Between groups and is the SSB divided by its degrees of freedom.

$$MS = \frac{SSB}{n - 1} \quad [11]$$

Within Group Variation (SSW)

The variation due to differences within individual samples is denoted SSW for Sum of Squares Within laboratories. The degrees of freedom are equal to the sum of the individual degrees of freedom for each sample. Since each sample has degrees of freedom (*df*) equal to one less than their sample sizes, and there are *k* samples, the total degrees of freedom is P less than the total sample size: $df = n - P$.

$$SSW = \sum df \cdot s^2 \quad [12]$$

The variance due to the differences within individual samples is denoted MSW for Mean Square Within groups. This is the within group variation divided by its degrees of freedom:

$$MSW = \frac{SSW}{P - n} \quad [13]$$

From equations [9] through [13], the ANOVA table as shown in Table 10 is developed.

Table 10. A single-factor ANOVA table showing key elements. Where P is the total number of groups, or laboratories. P-1 is 1 less than number of laboratories, P (n-1) is the number of data values minus number of groups (equals degrees of freedom for each group added together), and P-1 + P(n-1) is 1 less than the number of data points. MS is the mean squares of between laboratories and within laboratories. After Ellison *et al.*, (2009), Table 6.2, page 61.

Source	Sum of Squares	df	Mean Sum of Squares	F	p	F _{crit}
Between Laboratories	SSB	P-1	MSB=SSB/df	MSB/MSW	=FDIST(x,df,df)	F-table
Within Laboratories	SSW	P(n-1)	MSW=SSW/df	–	–	–
Total	SSB+SSW	P-1 + P(n-1)	–	–	–	–

Combined Standard Uncertainty

The combined standard uncertainty (u_c) represents the effects of random events such as days, instruments, and analysts on the precision of the analytical procedures of all accepted data of the participating laboratories. Using the output from ANOVA, the combined standard uncertainty (u_c) is determined from the square root of the sum of squares of the variances of the within laboratory repeatability, s_r , and the between laboratory precision, s_s :

$$u_c = \sqrt{s_r^2 + s_s^2} \quad [14]$$

Within laboratory repeatability is determined as

$$s_r = \sqrt{MSB} \quad [15]$$

and, the between laboratory precision as

$$s_s = \sqrt{\frac{(MSW - MSB)}{n}} \quad [16]$$

where MSW is the mean squares of the within laboratory variance, MSB is the mean squares for the between laboratories and n in this case, is the number of replicates in a group of the accepted data (Thompson & Lowthian, 2011).

Expanded Uncertainty

The expanded uncertainty (U) at a confidence level of 95% is determined by multiplication of the combined uncertainty (u_c) by a coverage factor (k) found from $N-1$ degrees of freedom (df), where N is the number of laboratory means accepted in the establishment of the certified value. The t-critical value for 5% significance can be found in a t-critical table (see Appendix 9, or from MS Excel as =TINV (5%, df)).

Uncertainty Statement

Typically, an uncertainty statement is presented as follows: Au =0.77±0.04 g/t, where the number following the symbol ± is the numerical value of an expanded uncertainty, $U = ku_c$, with U determined from a combined standard uncertainty multiplied by a coverage factor $k = 2$ or, a t-critical value for $N-1$ accepted laboratories. Since it can be assumed that the possible estimated values of the standard are approximately normally distributed with standard uncertainty, u_c , the certified value of the CRM is believed to lie in the interval defined by U with a level of confidence of approximately 95 %, e.g. a mean value of 0.77±0.04g/t will have intervals of: 0.73≤0.77≤0.81 g/t.

Appendix 3. Example: Comparison of Mean and Certified Value for Validation of Accuracy

According to ERM (2005); Eurolab (2007); Abzalov (2011) and Carr (2011), the validation of accuracy for a given mean and certified value requires the inclusion of the measurement uncertainty of the CRM in a t-test for statistical significance. The classical Student's t-test as shown in [17], does not consider the measurement uncertainty of the CRM. To compensate for this, Eurolab Technical Report No.1/2007 recommends equation [18] for the validation of CRMs with stated measurement uncertainties.

$$t_{calc} = \frac{|\bar{x} - \mu|}{\frac{s}{\sqrt{n}}} \quad [17]$$

$$t_{calc} = \frac{|\bar{x} - \mu|}{\sqrt{(u_{\mu})^2 + \frac{s^2}{n}}} \quad [18]$$

Where, t_{calc} is the calculated t-statistic, \bar{x} the mean of n replicates with a standard deviation of s for a CRM of μ certified value. The standard uncertainty u is the stated expanded uncertainty (U) of the CRM divided by the coverage factor (k) as stated on the certificate of analysis. Note that the $| \quad |$ bars indicate that the absolute value between the mean and the certified value is to be used, *i.e.* ignore the sign.

An example in which [18] is used for validation of accuracy is given below.

Example

A CRM is independently replicated nine times for Al_2O_3 concentration by XRF analysis, *i.e.* 9 individual fused glass beads were prepared. The observed mean and standard deviation of the replicate data are shown with the certified value and expanded uncertainty in Table 11. In validation of accuracy, the hypothesis question is: Is the difference between the observed mean and the certified value statistically significant at a level of confidence of 95%? Alternatively put, is there sufficient evidence to conclude that the data *i.e.* replicates generated, are inaccurate?

The relevant hypotheses are:

Null hypothesis: H_0 : Mean = Certified value of CRM with stated measurement uncertainty. The acceptance of H_0 means that accuracy is demonstrated; *i.e.* insufficient evidence to reject H_0 ;

Alternate hypothesis: H_1 : Mean \neq Certified value of CRM with stated measurement uncertainty. The acceptance of H_1 means that accuracy is not demonstrated, *i.e.* there is sufficient evidence to accept H_1 ;

Table 11. CRM certified value quoted expanded uncertainty U , the coverage factor for the CRM, $k=2.25$ and mean for $n=9$ replicates and corresponding standard deviation for the replicate data.

CRM Certified Value	Expanded % (U)	Coverage Factor (k)	Mean ($n=9$)	n	Standard Deviation (s)
4.62%	0.08%	2.25	4.59%	9	0.01015

The standard uncertainty (u) is found by dividing the expanded uncertainty by the coverage factor:

$$u = \frac{0.08}{2.25} = 0.0356 \%$$

Using the observed mean for the replicate data ($n=9$) obtained for the CRM and substituting into [18]:

$$t_{calc} = \frac{|\bar{x} - \mu|}{\sqrt{0.0356^2 + \frac{0.01015^2}{9}}} = \frac{|4.59 - 4.62|}{\sqrt{0.00126 + 0.00001145}} = 0.84$$

Therefore, $t_{calc} = 0.84$ and $t_{crit}(5\%, 8) = 2.31$ (df is 8, therefore, $t_{crit}=2.31$, see Appendix 9, page 35) which is >0.84 . Similarly, the p -value= 0.43 which is >0.05 . This is strong evidence in favour of accepting the null hypothesis that there is no significant statistical difference between the certified value and the observed mean. Therefore, under the conditions that the uncertainty associated with the certified value is known, the accuracy is validated for the CRM tested. If the null hypothesis is accepted that the mean obtained is not statistically different from the certified value, then the principle of traceability has been conformed to.

Appendix 4. Two-standard Deviations

Two-standard deviations are calculated using the expression:

$$\text{Two standard deviations} = 2 (u_c) \quad [19]$$

Where, u_c is the standard combined uncertainty (equation [14]).

Appendix 5. Confidence Interval

Confidence interval is calculated as:

$$\text{Confidence Interval (CI)} = \frac{(t_{critical})s}{\sqrt{N}} \quad [20]$$

Where, N is the number of laboratories (accepted laboratory data), $t_{critical}$ is a two-tailed value for $N - 1$ degrees of freedom (df) and s , is the standard deviation of the accepted laboratory means. A two-tailed critical value is found for $N - 1$ degrees of freedom from either a t -distribution table (Appendix 9) or MS Excel as =TINV (5%, df).

Appendix 6. Using the CRM in Quality Control

QC chart control limits should not be determined by the certified value and stated measurement uncertainty of the certified reference material used. These parameters although “certified” will never be known; it is only the corresponding statistical estimates, *i.e.* standard deviation and the mean calculated from replicated results that are known and these should be used in quality control charts. However, should the laboratory choose to use the certified value as the mean then the quoted $2s$, or CI value for the CRM can be used in the quality control chart.

It is recommended that a Shewhart chart be developed for the use if this CRM is to be used as a control sample in laboratory quality control. A Shewhart chart is a plot of sequential assay results obtained from quality control material such as an AMIS CRM. The warning and control limits are based on the standard deviation obtained from the mean of the replicates of a CRM (Ellison, *et al.*, 2009; Thompson, 2010). The procedure in preparing a Shewhart chart is as follows:

1. Analyse 10 to 15 replicates or more of the AMIS CRM.
2. Apply the Grubbs test for outliers.
3. Determine the mean of the replicates after application of the Grubbs test.
4. Determine the standard deviation, using equation [21], of the replicates following Grubbs test.
5. Calculate the standard deviation, s from:

$$s = \sqrt{\frac{\sum(x_i - \bar{x})^2}{n - 1}} \quad [21]$$

where, x_i is an individual measurement in the data set, \bar{x} is the mean of the data set at $n-1$ degrees of freedom (df) and n is the number of replicates. The sample standard deviation can be found using the MS Excel formula “=stdev.s (number1;)”.

6. Verify accuracy of the mean value using equation [18].
7. Once accuracy is verified, calculate $\pm 2s$ and $\pm 3s$, where s is the standard deviation calculated from [21].
8. Construct the Shewhart control chart around the mean of n replicates.
9. Use $\pm 2s$ as the warning limits.
10. Use $\pm 3s$ as the control limits.
11. It is recommended that if 2 to 3 points are outside the warning limits analyse another sample and if it is then within warning limits, continue. If it is outside the warning limits, stop and troubleshoot.
12. It is recommended that if any point is outside control limits, analyse another portion (sample) of the CRM. If it is within control limits, continue. If it is outside control limits, stop and troubleshoot.
13. For reference purposes, the CRM certified value can be plotted on the Shewhart chart alongside the mean value.

On a regular basis the accuracy of the replicates of the CRM should be assessed in terms of the certified value of the CRM using equation [18].

Appendix 7. Conversion to Air-dry Basis (Prepared by Allan Fraser)

Since AMIS certified analyte values are reported on a dry-basis, the user laboratory is required to dry a portion (accurately weigh out 1.0 grams in duplicate) of the CRM material in air at 105°C in a drying oven to constant mass to determine the moisture content. Use a crucible with a flat inner surface with a surface area not smaller than 10 cm² with the CRM material spread evenly over same; this represents a 0.1 gram spread per cm². In correcting the certified value for moisture content, a moisture correction factor is calculated:

$$\text{Moisture correction factor (MCF)} = \frac{100 - \% \text{Moisture at } 105^{\circ}C}{100} \quad [22]$$

$$\text{Air dry basis concentration} = \text{MCF} \times \text{certified value on a dry basis} \quad [23]$$

Example

The moisture content determined at 105°C on a CRM is 0.500%. The certified analyte concentration for the CRM is 12.62±0.52% (dry basis). Calculating the moisture correction factor using [22] gives:

$$\text{Moisture correction factor} = \frac{100 - 0.500}{100} = 0.995$$

Multiplying the factor of 0.995 by the certified value as stated on the certificate of analysis on a dry basis (as in [23]) gives the analyte concentration on an air-dry basis:

$$0.995 \times 12.62\% = 12.56\%$$

The stated measurement uncertainty also needs to be corrected using [22] and [23], e.g. $0.995 \times 0.52 = 0.51_{(7)}$, rounded to 0.52%. The air-dry basis concentration *i.e.* $12.56 \pm 0.52\%$ is to be used as the certified value with its corresponding measurement of uncertainty.

Appendix 8. Example of Determination of LOD and LOQ in Fire Assay

The limit of detection (LOD) is the minimum detectable quantity of the analyte of interest (Skoog & West, 1985). To determine the LOD in fire assay by lead collection, the minimum mass that an assay microbalance is capable of weighing (m in micrograms, and the original test sample mass, $Mass_{assay}$ in grams) determines the LOD. The smallest prill mass most assay microbalances can measure is $1\mu\text{g}$ or 0.001mg . Even with a microscope it may be difficult to locate and pick up a prill weighing ten times that amount (*i.e.* 0.01mg or $10\mu\text{g}$) and weigh it. If an analyst can weigh a prill of $1\mu\text{g}$ then the LOD becomes $1\mu\text{g}$. However, the concentration factor would be 50 times for a 50-gram assay sample and therefore the LOD in g/t becomes $1\mu\text{g}$ divided by the original mass of the sample in grams taken for fire assay [24]. Therefore, the LOD in fire assay is computed as:

$$LOD = \frac{m (\mu\text{g})}{Mass_{assay} (g)} (\text{g/t}) \quad [24]$$

The limit of quantitation (LOQ), is simply the LOD multiplied by 10 (Long & Winefordner, 1983):

$$LOQ = 10 \cdot \frac{m (\mu\text{g})}{Mass_{assay} (g)} (\text{g/t}) \quad [25]$$

Therefore, with a sample mass of 50g taken for fire assay, the limit of detection would be 0.02g/t . *i.e.* $1\mu\text{g} = 1\text{g/t}$, therefore $1\mu\text{g}/50\text{g} = 0.02\text{g/t}$. If no prill was found, then the LOD result would be $<0.02\text{g/t}$ or "not detected". Using a larger assay sample mass improves the LOD and LOQ (Table 11). Table 13 gives a recommended reporting scheme for LOD and LOQ.

Table 12. Mass of assay sample and corresponding limit of detection and limit of quantitation for an assay microbalance capability of smallest prill mass of $1\mu\text{g}$ or 0.001mg .

Mass Assay Sample (g)	LOD (g/t)	LOQ (g/t)
30	0.03	0.3
50	0.02	0.2
100	0.01	0.1

Table 13. Recommended reporting scheme for LOD and LOQ in fire assay.

Data	Report as
<LOD	Not detected
<LOQ	Detected
\geq LOQ	Report assay result

Appendix 9. T-distribution table

Table 14. T-distribution table for t-critical values (t crit.) for a two-tailed t-test at a 95% level of confidence.

<i>df</i>	Two-tailed	<i>df</i>	Two-tailed
1	12.71	23	2.06
2	4.30	24	2.06
3	3.18	25	2.06
4	2.78	26	2.05
5	2.57	27	2.05
6	2.44	28	2.04
7	2.36	29	2.04
8	2.30	30	2.04
9	2.26	35	2.03
10	2.22	40	2.02
11	2.20	45	2.01
12	2.17	50	2.00
13	2.16	55	2.00
14	2.14	60	2.00
15	2.13	70	1.99
16	2.12	80	1.98
17	2.11	90	1.98
18	2.10	100	1.98
19	2.09	120	1.98
20	2.08	Infinity	1.96
21	2.08		
22	2.07		

End of certificate