



AMIS_Documents

Doc: ADOC_074

Originator:
Quality Specialist

Approver:
Managing
Director

Revision No: 006

Revision Date:26.05.2019

Issued By: Management Rep

Certificate

AMIS0644

Certified Reference Material

Gold and Uranium Wits Supergroup Tailings, South
Africa

Certificate of Analysis

AMIS

A: 11 Avalon Road, West Lake View, Ext 11, Modderfontein, 1609, South Africa

P: PO Box 856, Isando, 1600, Gauteng, South Africa

T: +27 (0) 11 923 0800

W: www.amis.co.za



Table of Contents

Table of Contents	2
List of Tables	3
List of Appendices	3
Summary Statistics	4
1. Certified Concentrations and Uncertainties	6
2. Statistical Comparison of Means	8
3. Intended Use.....	8
4. Abbreviations and Symbols.....	9
5. Uncertified Concentration Values	10
6. Units	10
7. Analytical and Physical Methods	10
8. Origin of Material.....	10
9. Approximate Mineral and Chemical Composition	10
10. Quantitative Analysis by X-Ray Diffraction	11
11. Health and Safety	12
12. Method of Preparation.....	12
13. Particle Size Determination	12
14. Radioactivity.....	13
15. Handling.....	13
16. Storage information.....	13
17. Methods of Analysis Requested.....	13
18. Information Requested of Participating Laboratories	14
19. Certification of Mean and Estimation of Measurement Uncertainty.....	14
20. Two Standard Deviations	14
21. Confidence Interval	14
22. Expanded Uncertainty	15
23. Confidence Interval and Expanded Uncertainty.....	15
24. Participating Laboratories.....	15
25. Accepted Assay Data	16
26. Reported Values	22
27. Validation of Accuracy (Trueness).....	22
28. Limit of Detection and Limit of Quantitation in Gravimetric Fire Assay	22
29. Metrological Traceability.....	22
30. Period of Validity	22
31. Minimum Sample Size	22
32. Availability	22
33. Recommended use in Quality Control	23
34. Legal Notice	23
References	24
Appendices	25
End of certificate	34

List of Tables

Table 1. Certified concentrations, two standard deviations, combined and expanded uncertainty.	6
Table 2. Certified major oxides concentrations, two standard deviations, combined and expanded uncertainty.....	7
Table 3. The results of a two-sample equal or unequal variance t-test (two-tailed) data sets in which different analytical methods /instrumentation were used.....	8
Table 4. Abbreviations, symbols and descriptions.	9
Table 5. Results of XRD analysis.	12
Table 6. Particle Size Determination by laser diffraction.....	13
Table 7. Example of replicate assay data in which the <i>CI</i> , <i>2s</i> and <i>U</i> are shown.....	15
Table 8. Data used to calculate the certified values after removal of outliers.	16
Table 9. Uncertified element concentrations statistics.	25
Table 10. A single-factor ANOVA table showing key elements. Where <i>P</i> is the total number of groups, or laboratories. <i>P</i> -1 is 1 less than number of laboratories, <i>P</i> (<i>n</i> -1) is the number of data values minus number of groups (equals degrees of freedom for each group added together), and <i>P</i> -1 + <i>P</i> (<i>n</i> -1) is 1 less than the number of data points. <i>MS</i> is the mean squares of between laboratories and within laboratories. After Ellison <i>et al.</i> , (2009), Table 6.2, page 61.....	29
Table 11. CRM certified value quoted expanded uncertainty <i>U</i> , the coverage factor for the CRM, <i>k</i> =2.25 and mean for <i>n</i> =9 replicates and corresponding standard deviation for the replicate data.	30
Table 12. Mass of assay sample and corresponding limit of detection and limit of quantitation for an assay microbalance capability of smallest prill mass of 1µg or 0.001mg.	33
Table 13. Recommended reporting scheme for LOD and LOQ in fire assay.....	33
Table 14. T-distribution table for t-critical values (t crit.) for a two-tailed t-test at a 95% level of confidence.....	34

List of Appendices

Appendix 1: Uncertified Element Statistics.....	25
Appendix 2. Certification of Reference Material and Estimation of Measurement Uncertainty.....	26
Appendix 3. Example: Comparison of Mean and Certified Value for Validation of Accuracy	30
Appendix 4. Two-standard Deviations	31
Appendix 5. Confidence Interval.....	31
Appendix 6. Using the CRM in Quality Control.....	31
Appendix 7. Conversion to Air-dry Basis	32
Appendix 8. Example of Determination of LOD and LOQ in Fire Assay	33
Appendix 9. T-distribution table	34

Summary Statistics

Recommended Concentrations and Limits (at two Standard Deviations) Certified Concentrations

Analyte	Method	Certified (μ) ⁸	(2s) ¹⁰ \pm	Unit
Au	Pb Collection ¹	0.151	0.015	g/t
U	4A_MICP ²	53	6	ppm
U	XRF ³	58	14	ppm
Moisture	Moisture ⁴	0.28	0.02	%
SG	SG ⁵	2.78	0.16	Dimensionless
S	Combustion/LECO ⁶	0.53	0.04	%
S	4A_MICP	0.57	0.03	%
Al	4A_MICP	5.85	0.38	%
As	4A_MICP	110	9	ppm
Ba	4A_MICP	310	36	ppm
Be	4A_MICP	1	0.1	ppm
Bi	4A_MICP	0.9	0.05	ppm
Ca	4A_MICP	5693	454	ppm
Co	4A_MICP	36	3	ppm
Cr	4A_MICP	235	18	ppm
Cs	4A_MICP	2	0.2	ppm
Cu	4A_MICP	63	5	ppm
Fe	4A_MICP	2.05	0.18	%
Fe	FUS ⁷	2.07	0.27	%
Ga	4A_MICP	13	0.8	ppm
Hf	4A_MICP	3	0.4	ppm
K	4A_MICP	1.34	0.056	%
K	FUS	1.25	0.14	%
Li	4A_MICP	9	0.8	ppm
Mg	4A_MICP	2331	188	ppm
Mg	FUS	2359	305	ppm
Mn	4A_MICP	561	39	ppm
Mn	FUS	537	59	ppm
Mo	4A_MICP	3	1	ppm
Na	4A_MICP	1536	150	ppm
Nb	4A_MICP	4	0.7	ppm
Ni	4A_MICP	92	9	ppm
P	4A_MICP	158	20	ppm
Pb	4A_MICP	88	4	ppm
Rb	4A_MICP	57	2	ppm
Sb	4A_MICP	3	0.5	ppm
Sr	4A_MICP	53	8	ppm
Ta	4A_MICP	0.6	0.2	ppm
Th	4A_MICP	23	4	ppm
Ti	4A_MICP	1038	127	ppm
Ti	FUS	1639	95	ppm
Tl	4A_MICP	0.5	0.03	ppm
V	4A_MICP	42	5	ppm
W	4A_MICP	0.6	0.2	ppm
Y	4A_MICP	9	2	ppm
Zn	4A_MICP	144	10	ppm
Zr	4A_MICP	100	31	ppm

Major Oxides
Certified Concentrations (at two Standard Deviations)

Analyte	Method	Certified (μ) ⁸	(2s) ¹⁰ \pm	Unit
U ₃ O ₈	4A_MICP	0.006	0.0007	%
U ₃ O ₈	XRF	0.007	0.00	%
Al ₂ O ₃	XRF	11.08	0.22	%
CaO	XRF	0.75	0.02	%
Cr ₂ O ₃	XRF	0.037	0.01	%
Fe ₂ O ₃	FUS	2.96	0.39	%
Fe ₂ O ₃	XRF	2.87	0.076	%
K ₂ O	FUS	1.51	0.17	%
K ₂ O	XRF	1.58	0.021	%
MgO	FUS	0.39	0.05	%
MgO	XRF	0.39	0.02	%
MnO	FUS	0.069	0.01	%
Na ₂ O	XRF	0.20	0.03	%
P ₂ O ₅	XRF	0.037	0.01	%
SiO ₂	XRF	78.97	0.89	%
SO ₃	XRF	1.36	0.018	%
TiO ₂	FUS	0.27	0.02	%
TiO ₂	XRF	0.28	0.01	%

1. Certified Concentrations and Uncertainties

AMIS0644 is a new standard material, developed and certified in June 2019. Table 1 gives the certified concentrations, confidence interval, combined and expanded uncertainty for the certified reference material. Table 2 shows the certified major oxides concentrations, two standard deviations, confidence interval, combined and expanded uncertainty.

Table 1. Certified concentrations, two standard deviations, combined and expanded uncertainty.

Analyte	Method	Certified (μ) ⁸	N	n	k	% RSD	(u_c) ⁹	(2s) ¹⁰ \pm	(CI) ¹¹ 95%	(U) ¹² \pm	Unit
Au	Pb Collection ¹	0.151	8	62	2.365	5	0.0076	0.015	0.0050	0.02	g/t
U	4A_MICP ²	53	3	24	4.303	5	3	6	7	12	ppm
U	XRF ³	58	2	16	12.706	12	7	14	41	86	ppm
Moisture	Moisture ⁴	0.28	2	16	12.706	4	0.01	0.02	0.06	0.1	%
SG	SG ⁵	2.78	4	31	3.182	3	0.079	0.16	0.13	0.3	Dimensionless
S	Combustion/LECO ⁶	0.53	5	40	2.776	4	0.02	0.04	0.02	0.05	%
S	4A_MICP	0.57	4	30	3.182	2	0.01	0.03	0.02	0.04	%
Al	4A_MICP	5.85	4	31	3.182	3	0.19	0.38	0.27	0.6	%
As	4A_MICP	110	4	31	3.182	4	5	9	5	14	ppm
Ba	4A_MICP	310	5	40	2.776	6	18	36	21	50	ppm
Be	4A_MICP	1	4	30	3.182	5	0.06	0.1	0.07	0.2	ppm
Bi	4A_MICP	0.9	2	15	12.706	3	0.02	0.05	0.2	0.3	ppm
Ca	4A_MICP	5693	4	31	3.182	4	227	454	323	722	ppm
Co	4A_MICP	36	4	30	3.182	4	1	3	2	4	ppm
Cr	4A_MICP	235	3	24	4.303	4	9	18	18	38	ppm
Cs	4A_MICP	2	3	24	4.303	4	0.1	0.2	0.1	0.4	ppm
Cu	4A_MICP	63	5	39	2.776	4	3	5	2	7	ppm
Fe	4A_MICP	2.05	4	31	3.182	4	0.088	0.18	0.13	0.3	%
Fe	FUS ⁷	2.07	3	23	4.303	7	0.14	0.27	0.34	0.6	%
Ga	4A_MICP	13	3	23	4.303	3	0.4	0.8	0.4	2	ppm
Hf	4A_MICP	3	3	23	4.303	8	0.2	0.4	0.5	1	ppm
K	4A_MICP	1.34	3	24	4.303	2	0.028	0.056	0.060	0.1	%
K	FUS	1.25	2	16	12.706	6	0.071	0.14	0.64	0.9	%
Li	4A_MICP	9	3	23	4.303	4	0.4	0.8	0.4	2	ppm
Mg	4A_MICP	2331	4	31	3.182	4	94	188	133	299	ppm
Mg	FUS	2359	3	23	4.303	6	153	305	352	656	ppm
Mn	4A_MICP	561	4	32	3.182	3	20	39	28	62	ppm
Mn	FUS	537	3	23	4.303	6	30	59	38	128	ppm
Mo	4A_MICP	3	4	32	3.182	19	0.5	1	0.7	2	ppm
Na	4A_MICP	1536	5	38	2.776	5	75	150	86	208	ppm
Nb	4A_MICP	4	3	22	4.303	8	0.3	0.7	0.7	1	ppm
Ni	4A_MICP	92	4	32	3.182	5	4	9	6	14	ppm
P	4A_MICP	158	4	32	3.182	6	10	20	14	33	ppm
Pb	4A_MICP	88	5	37	2.776	2	2	4	1	6	ppm
Rb	4A_MICP	57	2	15	12.706	2	0.9	2	4	12	ppm
Sb	4A_MICP	3	3	24	4.303	9	0.2	0.5	0.5	1	ppm
Sr	4A_MICP	53	4	32	3.182	7	4	8	6	13	ppm
Ta	4A_MICP	0.6	3	24	4.303	18	0.1	0.2	0.2	0.4	ppm
Th	4A_MICP	23	3	24	4.303	9	2	4	5	9	ppm
Ti	4A_MICP	1038	4	31	3.182	6	63	127	76	202	ppm
Ti	FUS	1639	2	16	12.706	3	48	95	334	604	ppm
Tl	4A_MICP	0.5	2	16	12.706	3	0.02	0.03	0.1	0.2	ppm
V	4A_MICP	42	5	37	2.776	7	3	5	3	8	ppm
W	4A_MICP	0.6	3	22	4.303	15	0.09	0.2	0.2	0.4	ppm
Y	4A_MICP	9	3	23	4.303	10	0.9	2	2	4	ppm
Zn	4A_MICP	144	5	37	2.776	3	5	10	6	13	ppm
Zr	4A_MICP	100	3	23	4.303	16	16	31	40	68	ppm

Table 2. Certified major oxides concentrations, two standard deviations, combined and expanded uncertainty.

Analyte	Method	Certified (μ) ⁸	N	n	k	% RSD	(u_c) ⁹	(2s) ¹⁰ ±	(CI) ¹¹ 95%	(U) ¹² ±	Unit
U ₃ O ₈	4A_MICP	0.006	3	24	4.303	5	0.0003	0.0007	0.001	0.001	%
U ₃ O ₈	XRF	0.007	2	15	12.706	11	0.001	0.00	0.004	0.01	%
Al ₂ O ₃	XRF	11.08	5	40	2.776	1	0.11	0.22	0.13	0.3	%
CaO	XRF	0.75	6	46	2.571	2	0.01	0.02	0.01	0.03	%
Cr ₂ O ₃	XRF	0.037	5	37	2.776	12	0.004	0.01	0.005	0.01	%
Fe ₂ O ₃	FUS	2.96	3	23	4.303	7	0.19	0.39	0.48	0.8	%
Fe ₂ O ₃	XRF	2.87	6	48	2.571	1	0.038	0.076	0.038	0.1	%
K ₂ O	FUS	1.51	2	16	12.706	6	0.085	0.17	0.77	1	%
K ₂ O	XRF	1.58	6	46	2.571	1	0.011	0.021	0.010	0.03	%
MgO	FUS	0.39	3	23	4.303	7	0.03	0.05	0.06	0.1	%
MgO	XRF	0.39	5	39	2.776	3	0.01	0.02	0.01	0.03	%
MnO	FUS	0.069	3	23	4.303	6	0.004	0.01	0.005	0.02	%
Na ₂ O	XRF	0.20	6	46	2.571	7	0.01	0.03	0.01	0.04	%
P ₂ O ₅	XRF	0.037	5	40	2.776	8	0.003	0.01	0.003	0.01	%
SiO ₂	XRF	78.97	5	40	2.776	1	0.44	0.89	0.51	1	%
SO ₃	XRF	1.36	2	16	12.706	1	0.0092	0.018	0.052	0.1	%
TiO ₂	FUS	0.27	2	16	12.706	3	0.008	0.02	0.05	0.1	%
TiO ₂	XRF	0.28	5	37	2.776	2	0.005	0.01	0.006	0.01	%

1. Pb Collection
2. 4A_MICP is a Multi-acid digestion with either ICPOES/ICPMS/AAS finish
3. XRF is X-ray Fluorescence
4. Moisture
5. SG is Specific Gravity
6. Combustion/LECO
7. FUS is Fusion digestion with ICP finish
8. The certified value μ , is an unweighted grand mean of the means of N accepted sets of data from different laboratories and n number of test sample replicates. The certified value is traceable to SI units and is reported on a dry basis.
9. The combined uncertainty of the certified value is the within-laboratory reproducibility standard deviation derived from the analysis of variance of results from N number of laboratories and n number of sample replicates. (u_c)
10. Two standard deviations (2s)
11. Confidence interval at 95% level of confidence.
12. Expanded uncertainty (U) at a confidence level of 95% is determined by multiplication of the combined uncertainty (u_c) with a coverage factor (k) found from N-1 degrees of freedom (see Appendix 7 for t-distribution table). Example: $U = 2.36 \times 0.23 = 0.5\%$

2. Statistical Comparison of Means

A comparison of means for replicate data for the same element concentration determined by different analytical methods is done equating the variances between the two data sets; if the variances are found to be equal (F-test, p -value >0.05), then an equal variance t-test is applied. Should the variances be statistically significant, i.e. $p<0.05$, then an unequal variance t-test is performed. For either t-test, if the obtained p -value ≥ 0.05 , the null hypothesis that the means (certified values) are equal is accepted (Table 3). This gives the analyst confidence in the certified values reported by different analytical methods on the same analyte.

Table 3. The results of a two-sample equal or unequal variance t-test (two-tailed) data sets in which different analytical methods /instrumentation were used.

Method	Certified value	Method	Certified value	F-Test Outcome	p -value (t-test)	t-test Outcome
U 4A_MICP	0.005 %	U XRF	0.006 %	Equal Variance ($p=0.251$)	0.207	Accept H_0 ; certified values are equal
U ₃ O ₈ 4A_MICP	0.006 %	U ₃ O ₈ XRF	0.007 %	Equal Variance ($p=0.3$)	0.157	Accept H_0 ; certified values are equal
S Combustion/LECO	0.528 %	S 4A_MICP	0.566 %	Equal Variance ($p=0.191$)	0.008	Reject H_0 ; certified values are <i>not</i> equal
Fe ₂ O ₃ FUS	2.96 %	Fe ₂ O ₃ XRF	2.87 %	Unequal Variance ($p=0.002$)	0.526	Accept H_0 ; certified values are equal
MgO FUS	0.393 %	MgO XRF	0.392 %	Unequal Variance ($p=0.037$)	0.977	Accept H_0 ; certified values are equal
TiO ₂ FUS	0.273 %	TiO ₂ XRF	0.277 %	Equal Variance ($p=0.264$)	0.435	Accept H_0 ; certified values are equal
Fe FUS	2.07 %	Fe 4A_MICP	2.05 %	Equal Variance ($p=0.211$)	0.874	Accept H_0 ; certified values are equal
K FUS	1.25 %	K 4A_MICP	1.34 %	Equal Variance ($p=0.1$)	0.126	Accept H_0 ; certified values are equal
Mg FUS	0.236 %	Mg 4A_MICP	0.233 %	Equal Variance ($p=0.201$)	0.748	Accept H_0 ; certified values are equal
Mn FUS	0.054 %	Mn 4A_MICP	0.056 %	Equal Variance ($p=0.455$)	0.126	Accept H_0 ; certified values are equal
Ti FUS	0.164 %	Ti 4A_MICP	0.104 %	Equal Variance ($p=0.507$)	0.001	Reject H_0 ; certified values are <i>not</i> equal

3. Intended Use

AMIS0644 is a matrix matched Certified Reference Material, fit for use as a control sample in routine assay laboratory quality control when inserted within runs of test samples and measured in parallel to test samples. This material can also be used for method development, use as independent calibration verification check standard (i.e. if not used as a calibration standard in an instrument calibration), or for validation of accuracy in a method validation exercise (see Appendix 3). The recommend procedure for the use of this CRM as a control standard in laboratory quality control is to develop a Shewhart chart, where a mean value and corresponding 1, 2 and 3 standard deviations are derived from replicate measurements of the CRM (see Appendix 6). This CRM can also be used to assess inter-laboratory or

instrument bias and establish within-laboratory precision and within-laboratory reproducibility. The certified concentrations and expanded uncertainty for this material are property values based on an inter-laboratory measurement campaign and reflect consensus results from the laboratories that took part in the exercise.

4. Abbreviations and Symbols

Abbreviations and symbols used in this document are shown in Table 4.

Table 4. Abbreviations, symbols and descriptions.

Abbreviation/Symbol	Description
Alpha (α)	Significance level (denoted by alpha, ' α ') of 0.05 or 5%
ANOVA	Analysis of variance by statistical means
Bq	The becquerel is the SI derived unit of radioactivity.
BIF	Banded iron formation
CRM	Certified reference material
df	Degrees of freedom, typically, $n-1$, or $N-1$
F_{calc}	Calculated F statistic from ANOVA or Fisher's test
F-critical or F_{crit}	F-critical value from F-distribution table
GOI	Gain on ignition
H_0	Null hypothesis
H_1	Alternate hypothesis
g/t	Grams per tonne
k	Coverage factor, e.g. $k=2$ for 95% level of confidence
LOC	Level of confidence or confidence level
LOD	Limit of detection
LOQ	Limit of quantitation
LOI	Loss on ignition
MS	Mean squares (ANOVA)
MSb	Mean squares between (ANOVA)
MSw	Mean squares within (ANOVA)
N	Number of labs
n	Number of replicates
μ	Property or certified value of a CRM
p	' p -value' a measure of the strength of evidence against H_0
P	Total number of data points in ANOVA
ppm	Parts per million. Equivalent to g/t
RSD	Relative standard deviation usually expressed as % at a 68% LOC
Replicates	Replication is the repetition of an experimental condition so that the variability associated with an analysis can be estimated (ASTM E1847)
s	Standard deviation
s_r	Within laboratory repeatability as derived from ANOVA
s_s	Between laboratory standard deviation as derived from ANOVA
SS	Sum of squares in ANOVA

Abbreviation/Symbol	Description
SST	Total variation in ANOVA
SSB	Between group (laboratory) variance
SSW	Within group (laboratory) variance
2s	Two times standard deviation
SI	Standard International system of units
t_{calc}	Calculated t statistic from a one-sample, two-tailed t-test
t-critical or t_{crit}	t-critical value at given alpha and degrees of freedom
Tonne	A metric ton, is a unit of mass equaling 1000 kilograms
=TINV(5%, <i>df</i>)	MS Excel function for t-critical value at LOC 95% and <i>df</i>
U	Expanded uncertainty at a given k
u	Standard uncertainty at k=1
u_c	Combined standard uncertainty at k=1
μm	Micron, is an SI derived unit of length equaling 1×10^{-6} of a meter

5. Uncertified Concentration Values

Appendix 1 gives uncertified concentrations for other elements present in the CRM.

6. Units

All results for major oxides are reported as oxides in percentages. All results for major elements analyses reported in percentages or ppm. Results for Au and the platinum group elements are reported in g/t or ppm. Specific gravity (SG) is the ratio of the density of a substance to the density of a reference substance, *i.e.* equivalently; it is the ratio of the mass of a substance to the mass of a reference substance for the same given volume. Since specific gravity is a ratio of densities its units are therefore dimensionless.

7. Analytical and Physical Methods

A complete list of analytical and physical methods as generic method codes with a brief description of the methods is available on the AMIS web site www.amis.co.za

8. Origin of Material

This standard is made from uranium bearing gold ore originating from the Ventersdorp, Contact Reef, Carbon Leader Reef and Vaal Reef. The uranium bearing ore was blended with low grade material of similar matrix to obtain the required grade.

9. Approximate Mineral and Chemical Composition

Throughout the Witwatersrand, gold ores occur in sheets or reefs originally deposited horizontally under water. The thickness of the reefs range between a lie of grit to several meters, with an average of 300mm. The sediments were also intersected by dykes and sills of dolomite, diabase and syenite intruding existing faults. The reefs can be in the form of either coarse conglomerates, rock pebbles and cemented in a silicate matrix. Pebbles usually derived from vein quartz, may also consist of quartzitic, chert jasper and quartz porphyry. The matrix consists of pure silica, but also contains minute flakes of muscovite and pyrophyllite as well as visible pyrite and other sulphides. The gold is usually confined to the matrix and tends to be concentrated along bedding planes (Vermeulen, 2001).

10. Quantitative Analysis by X-Ray Diffraction

Both natural and synthetic materials have a specific chemistry and atomic arrangement, known as phases. Phases can be identified and quantified using X-ray diffraction (XRD) which produces a plot of the intensity of X-rays scattered at different angles by crystalline phases in a material. Essentially, an X-ray diffraction pattern is the sum of the diffraction patterns produced by each phase. Simply put, an X-ray diffraction pattern is a fingerprint that allows the identification of what is in a target sample material. Knowledge of the mineral phase composition is useful in method development with techniques such as ICP-OES and XRF as potential matrix effects and spectral interferences can be recognised and accounted for. X-ray diffraction is effective in that it allows the identification of different phases of compounds that are identical in chemistry, but have a distinctly different the atoms, e.g. quartz, cristobalite, and glass are all different phases of SiO₂. Where quantitative XRD results do not correspond to results of other analytical techniques, it should be borne in mind that even though the data are quantitative they are meant to be used for indicative purposes in development of other analytical methods. Mineral names may not reflect the actual compositions of minerals identified, but rather the mineral group.

Quantification is determined from the chosen software package: this uses the full-profile Rietveld method of refining the profile of the calculated XRD pattern against the profile of the measured XRD pattern. The total calculated pattern is the sum of the calculated patterns of the individual phases. Results are given as weight % of the total crystalline phases and amorphous content. The amorphous content quantifies the amorphous material and unknown minerals or known minerals for which there is not a suitable crystal structure.

Corrections are incorporated into the process that allows for a more accurate description of the mineral's contribution to the measured pattern and to allow for variation due to atomic substitution, layer disordering, preferred orientation, and other factors that affect the acquisition of the XRD scan.

The limitations of qualitative XRD analysis are as follows:

- The detection of a phase may be dependent on its crystallinity.
- Where there exist multiple phases, overlap of diffracted reflections can occur, thus rendering some ambiguity into the interpretation.
- Overlapping reflections of a major phase can mask the presence of minor or trace phases.
- Some phases cannot be unambiguously identified as they are present in minor or trace amounts.

The limitations of quantitative XRD analysis by a full-profile Rietveld method are as follows:

- The limitations for qualitative XRD analysis apply.
- The method as described is standardless: it relies solely on the published crystallographic data available for each phase. Some data may not exactly describe the phases present.
- Particle size is important with respect to the absorption of the X-rays by the sample.
- Micronising reduces the particle size to that more suitable for quantitative analysis.

The accuracy of the analysis is dependent on sampling and sample preparation in addition to the calculated profiles being exactly representative of the chemistry of the component phases and their crystallinity. Some preferred orientation effects and reflection overlaps may occur which cannot be adequately resolved.

Table 5. Results of XRD analysis.

Phase	Formula	Composite%
Amorphous Content		13
Chlorite	$(\text{Fe,Al,Mg})_6(\text{Si,Al})_4\text{O}_{10}(\text{OH})_8$	2
Gypsum	$\text{CaSO}_{4.2}\text{H}_2\text{O}$	1
Illite/Muscovite	$\text{K}(\text{Al,Mg,Fe})_2(\text{Si,Al})_4\text{O}_{10}(\text{OH})_2$	17
Pyrite	FeS_2	1
Pyrophyllite	$\text{Al}_2\text{Si}_4\text{O}_{10}(\text{OH})_2$	13
Quartz	SiO_2	54
Rutile	TiO_2	<1
Total		101%

For informational purposes only

11. Health and Safety

The material is a very fine powder coloured pale red (Corstor 5YR 6/2). Safety precautions for handling fine particulate matter are recommended, such as the use of safety glasses, breathing protection, gloves and a laboratory coat.

12. Method of Preparation

The particle size distribution for this material was shown to have a nominal top size of 54µm (95% passing 54µm). The procedure of preparation in brief is as follows: the material was crushed, dry-milled and air-classified to <54µm. It was then blended in a bi-conical mixer, systematically divided and sealed into 1kg Laboratory Packs. Explorer Packs are then subdivided from the Laboratory Packs as required. Final packaged units were then selected on a random basis and submitted for analysis to an independent laboratory accredited with the ISO17025 standard of general requirements for the competence of testing and calibration laboratories. The results obtained from this laboratory are then evaluated statistically by AMIS for homogeneity.

13. Particle Size Determination

The sample has been analysed using a Malvern Mastersizer 2000. Particles are passed through a focused laser beam that scatter light at an angle inversely proportional to their size. The intensity of light is measured and converted to a volume in particle size distribution. The results for this standard are presented in Table 6.

Table 6. Particle Size Determination by laser diffraction.

Size (µm)	Vol. Under %
<45um	93.0
<63um	95.4
<75um	96.2
<90um	97.1
<100um	97.6
<106um	97.8
<150um	99.3

For informational purposes only

14. Radioactivity

Shipments of this material require special labelling and placarding. AMIS0644 contains U (0.66 Bq/g) and Th 0.09 Bq/g) and is classified as EXEMPT MATERIAL in terms of "Safety Standards Series No. TS-R-1: Regulations for the Safe Transport of Radioactive Material, International Atomic Energy Agency, 2005, para 403, Table 1".

15. Handling

The material is packaged in Laboratory Packs and Explorer Packs that must be shaken or otherwise agitated before use. The analyte concentrations are quoted on a dry basis; therefore, the user needs to determine the moisture content to convert any obtained assay values to an air-dry basis (see Appendix 7 for an example calculation).

16. Storage information

The material should be stored in a cool dry place, in such a way that it does not compromise the integrity of the CRM. The material should be stored in conditions which will ensure it does not absorb moisture.

17. Methods of Analysis Requested

The following methods of analysis were requested:

- a) Au-Pb collection finished with either ICP-OES or ICP-MS or AAS or gravimetric
- b) Multi element scan to include all elements especially U-4-acid total digestion including HF and/or peroxide fusion finished with either ICP-OES or ICP-MS or AAS
- c) LOI and all major oxides including U3O8 with XRF finish and/or Peroxide fusion and/or 4 acid digest including HF finished with either ICP-OES or ICP-MS or AAS (Please specify the temperature for LOI)
- d) SG – gas pycnometer
- e) S and C Combustion/LECO
- f) Moisture

18. Information Requested of Participating Laboratories

The following information was requested of the participating laboratories for the development of this CRM:

- a) State aliquots used for all determinations.
- b) All results for major elements to be reported as oxides in percentages.
- c) All results for multi-element scans and fusion to be reported in ppm.
- d) All results for Au to be reported in ppm.
- e) Report all QC data, to include replicates, blanks and certified reference materials used.
- f) All Round robin samples must be treated the same as routine test samples.
- g) All results must be reported to maximum decimal places i.e. dependent on laboratories capabilities
- h) Please ensure moisture content is determined and calculated. All results should be corrected by the moisture correction factor and this factor should be stated in the laboratory results.
- i) Please use the excel template provided by AMIS. If you require a copy, please email any of the email addresses below. Ensure all uncertainties are added to the results.
- j) Please send excel and PDF of all results.
- k) Ensure correct PPE is used i.e. gloves, dust masks and protective clothing.
- l) Analysis should be done under controlled environmental conditions.

19. Certification of Mean and Estimation of Measurement Uncertainty

The samples used in this certification process have been selected in such a way as to represent the entire batch of material and were taken from the final packaged units; therefore, all sources of uncertainty are included in the combined standard uncertainty determination. Initially the data submitted by all the laboratories are subjected to a z-score test, equation [1] to exclude outliers and the remaining data sets examined for their normality in distribution. This is followed by the exclusion of further outliers as defined by the IUPAC Harmonized Protocol of 1995 in which both Cochran and a Grubbs test are applied until all outliers are identified, equations [2] and [3]. A grand mean and standard deviation is re-calculated using all remaining data (Thompson, 2008; Carr, 2011) (see Appendix 2)

20. Two Standard Deviations

AMIS reports two-standard deviations (2s) with all certified values. Two -standard deviations are calculated using the expression:

$$\text{Two standard deviations} = 2 (u_c)$$

Where u_c is the standard combined uncertainty (see Appendix, equation [14]).

21. Confidence Interval

AMIS reports a confidence interval (CI) with all certified values. Confidence interval as used by AMIS is:

$$\text{Confidence Interval (CI)} = \frac{(t_{critical})S}{\sqrt{N}}$$

Where, N is the number of laboratories (accepted laboratory data), $t_{critical}$ is a two-tailed value for $N - 1$ degrees of freedom (df) and s , is the standard deviation of the accepted laboratory means. A two-tailed critical value is found for $N - 1$ degrees of freedom from either a t -distribution table (Appendix 9) or MS Excel as =TINV (5%, df).

22. Expanded Uncertainty

ANOVA gives an estimate of the repeatability and the reproducibility of the data accepted for certification of the candidate reference material (see equations, [15] and [16], in the Appendix). Therefore, random variables (e.g., subsampling, instrument effects, interferences, operators and measurement conditions) that occur during the analysis of the candidate reference material by the various laboratories is considered. This approach does not necessarily quantify each individual source of uncertainty; however, the combined effect of random uncertainties is assessed (Ramsey & Ellison, 2007). A combined standard uncertainty is calculated from equation [14], which when multiplied by the t -critical value for $N-1$ laboratories, gives an *expanded uncertainty* at a 95% level of confidence. The expanded uncertainty is a measure of the doubt around the certified value at a level of 95% confidence. The expanded uncertainty is used in the validation of accuracy (see equation [18]).

23. Confidence Interval and Expanded Uncertainty

A combined standard uncertainty will be greater than a combined CI . This is because ANOVA considers the within-lab repeatability (that is repeatability within each lab group) as well as the repeatability between each lab data set. This attends to random variables that contribute to the measurement of uncertainty, during the analysis of the test sample at the participating laboratories. The within-lab repeatability and the between lab repeatability is combined as the square root of the sum of squares of these two values giving a combined standard uncertainty, at a 68% confidence level. Multiplying the combined standard uncertainty by the t -critical value for $N-1$, gives the expanded uncertainty at 95% level of confidence. It is recommended that the procedure described in Appendix 6, “Using the CRM in Quality Control” be used, in setting the limits of the CRM. Table 7 below shows mean gold values obtained by fire assay lead collection, for nine different laboratories, the confidence interval, two-standard deviations and expanded uncertainty.

Table 7. Example of replicate assay data in which the CI , $2s$ and U are shown.

Lab No.	Mean Au (g/t)	CI	0.0088
1	0.268	$2s$	0.031
2	0.273	U	0.04
3	0.270		
4	0.288		
5	0.274		
6	0.256		
7	0.263		
8	0.258		
9	0.288		

24. Participating Laboratories

The laboratories that are accredited with ISO17025 and provided timeous results are:

1. ALS Geochemistry Arabia
2. ALS Geochemistry Vancouver
3. Bureau Veritas Minerals Ultra Trace Pty Ltd
4. Intertek Mineral Tarkwa Ghana
5. Intertek Perth

6. SGS Mineral Services Lakefield (Canada)
7. SGS Vancouver (Canada)
8. Shiva Analyticals India
9. Super Laboratory Services (Pty) Ltd.

25. Accepted Assay Data

Data from the 9 laboratories used for certification are set out in Table 8.

Table 8. Data used to calculate the certified values after removal of outliers.

4A_MICP	XRF	XRF	XRF	XRF	FUS	XRF	FUS	XRF
U ₂ O ₅	U ₃ O ₈	Al ₂ O ₃	CaO	Cr ₂ O ₃	Fe ₂ O ₃	Fe ₂ O ₃	K ₂ O	K ₂ O
%	%	%	%	%	%	%	%	%
0.01	0.01	11.15	0.74	0.04	3.17	2.81	1.45	1.58
0.01	0.01	11.10	0.73	0.04	3.21	2.83	1.45	1.58
0.01	0.01	11.10	0.73	0.04	3.19	2.82	1.45	1.58
0.01	0.01	11.10	0.73	0.04	3.27	2.81	1.45	1.58
0.01	0.01	11.25	0.73	0.04	3.14	2.82	1.45	1.58
0.01	0.01	11.15	0.73	0.04	3.20	2.82	1.45	1.58
0.01	0.01	11.15	0.73	0.04	3.09	2.83	1.45	1.58
0.01	0.01	11.15	0.74	0.04	2.87	2.82	1.45	1.58
0.01	0.01	11.00	0.75	0.04	2.86	2.87	1.57	1.56
0.01	0.01	11.00	0.74	0.04	2.85	2.90	1.57	1.56
0.01	0.01	11.00	0.75	0.04	2.80	2.89	1.57	1.57
0.01	0.01	11.00	0.76	0.04	2.83	2.90	1.57	1.58
0.01	0.01	11.00	0.75	0.04	2.85	2.87	1.57	1.57
0.01	0.01	11.00	0.74	0.04	2.92	2.88	1.57	1.57
0.01	0.01	11.00	0.76	0.04	2.92	2.88	1.57	1.56
0.01		11.00	0.75	0.04	2.76	2.89	1.57	1.57
0.01		11.09	0.74	0.04	2.83	2.86		1.59
0.01		11.09	0.75	0.04	2.85	2.87		1.59
0.01		11.10	0.75	0.04	2.87	2.86		1.59
0.01		11.08	0.74	0.04	2.85	2.86		1.59
0.01		11.09	0.75	0.04	2.85	2.86		1.59
0.01		11.08	0.75	0.03	2.79	2.86		1.59
0.01		11.09	0.75	0.03	2.85	2.87		1.59
0.01		11.10	0.75	0.03		2.87		1.59
		10.90	0.75	0.03		2.84		1.56
		10.98	0.76	0.03		2.84		1.57
		11.04	0.76	0.03		2.87		1.57
		10.97	0.76	0.03		2.83		1.56
		10.94	0.75	0.03		2.84		1.58
		10.97	0.75	0.04		2.87		1.56
		10.91	0.75	0.04		2.82		1.60
		10.90	0.75	0.04		2.83		1.58
		11.23	0.74	0.04		2.89		1.59
		11.21	0.73	0.04		2.91		1.58
		11.16	0.73	0.04		2.92		1.58
		11.21	0.73	0.04		2.90		1.59
		11.23	0.74	0.04		2.89		1.58
		11.17	0.74			2.93		1.57
		11.21	0.73			2.93		1.58
		11.22	0.76			2.89		1.58
			0.76			2.91		1.58
			0.76			2.91		1.58
			0.75			2.90		1.57
			0.76			2.92		1.58
			0.76			2.91		1.58
			0.76			2.92		1.57
						2.91		
						2.90		

Assay Data (Cont.)

FUS	XRF	FUS	XRF	XRF	XRF	XRF	FUS	XRF
MgO	MgO	MnO	Na ₂ O	P ₂ O ₅	SiO ₂	SO ₃	TiO ₂	TiO ₂
%	%	%	%	%	%	%	%	%
0.44	0.39	0.07	0.19	0.03	78.50	1.36	0.28	0.28
0.41	0.39	0.06	0.19	0.04	78.30	1.36	0.28	0.28
0.44	0.39	0.07	0.20	0.03	78.50	1.37	0.28	0.28
0.42	0.39	0.07	0.19	0.03	78.60	1.36	0.27	0.28
0.40	0.39	0.07	0.19	0.03	78.50	1.36	0.27	0.28
0.41	0.40	0.07	0.19	0.04	78.50	1.36	0.28	0.28
0.44	0.39	0.06	0.19	0.04	78.30	1.36	0.28	0.28
0.39	0.38	0.07	0.20	0.03	78.50	1.36	0.28	0.28
0.39	0.38	0.06	0.19	0.04	79.00	1.38	0.27	0.28
0.38	0.39	0.07	0.18	0.04	78.80	1.37	0.27	0.28
0.37	0.38	0.07	0.18	0.04	79.20	1.36	0.27	0.28
0.37	0.39	0.07	0.17	0.04	79.00	1.35	0.28	0.28
0.37	0.40	0.07	0.17	0.04	79.00	1.36	0.27	0.28
0.39	0.39	0.07	0.18	0.04	79.20	1.37	0.27	0.28
0.38	0.39	0.07	0.18	0.03	78.90	1.38	0.27	0.28
0.38	0.39	0.07	0.18	0.04	79.40	1.37	0.27	0.28
0.38	0.39	0.08	0.22	0.04	79.26			0.28
0.38	0.38	0.07	0.21	0.03	79.42			0.28
0.38	0.39	0.07	0.22	0.03	79.34			0.28
0.36	0.38	0.07	0.22	0.03	79.22			0.28
0.38	0.38	0.07	0.22	0.04	79.19			0.28
0.36	0.38	0.07	0.22	0.04	79.27			0.28
0.38	0.38	0.07	0.22	0.03	79.31			0.28
	0.38		0.22	0.03	79.40			0.27
	0.40		0.20	0.04	78.62			0.28
	0.38		0.20	0.04	78.88			0.27
	0.39		0.20	0.04	79.13			0.27
	0.39		0.20	0.04	78.73			0.27
	0.40		0.20	0.04	78.50			0.28
	0.39		0.20	0.04	78.73			0.27
	0.41		0.20	0.04	78.15			0.28
	0.41		0.20	0.04	78.34			0.27
	0.41		0.20	0.04	79.50			0.27
	0.41		0.21	0.04	79.40			0.27
	0.41		0.20	0.04	79.46			0.27
	0.40		0.21	0.04	79.58			0.27
	0.40		0.21	0.04	79.29			0.27
	0.40		0.20	0.04	79.42			
	0.40		0.19	0.04	79.32			
			0.18	0.04	79.29			
			0.19					
			0.19					
			0.19					
			0.19					
			0.20					
			0.19					

Assay Data (Cont.)

Pb Collection	4A_MICP	XRF	Moisture	SG	Combustion/LECO	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP
Au	U	U	Moisture	SG	S	S	Al	As	Ba	Be	Bi
g/t	ppm	ppm	%	Dimensionless	%	%	ppm	ppm	ppm	ppm	ppm
0.14	50.00	59.00	0.29	2.70	0.50	0.57	58100	114	320	1.10	0.91
0.14	50.00	59.00	0.30	2.71	0.50	0.59	61200	116	340	1.20	0.87
0.15	50.00	51.00	0.29	2.71	0.50	0.59	60300	119	340	1.10	0.85
0.14	50.00	59.00	0.27	2.70	0.51	0.58	59000	114	330	1.10	0.87
0.14	50.00	51.00	0.29	2.71	0.50	0.57	59400	113	330	1.10	0.87
0.15	50.00	58.00	0.29	2.69	0.50	0.58	58500	111	330	1.10	0.87
0.15	50.00	51.00	0.28	2.71	0.52	0.57	57700	108	320	1.10	0.88
0.15	50.00	51.00	0.28	2.72	0.50	0.56	58700	114	330	1.10	0.86
0.15	54.96	60.00	0.28	2.85	0.55	0.56	59942	105	283	1.10	0.90
0.15	55.21	60.00	0.27	2.87	0.55	0.57	60793	105	280	1.20	0.90
0.16	55.84	60.00	0.28	2.90	0.55	0.57	59159	105	282	1.10	0.90
0.15	55.69	60.00	0.28	2.91	0.55	0.57	59537	106	285	1.20	0.90
0.16	55.48	70.00	0.28	2.84	0.55	0.56	60155	105	283	1.20	0.90
0.16	55.95	50.00	0.28	2.89	0.55	0.58	61359	105	284	1.20	0.90
0.15	55.92	60.00	0.27	2.91	0.55	0.55	61206	105	284	1.20	0.90
0.15	55.74	70.00	0.27	2.77	0.56	0.55	60666	107	280	1.20	
0.15	53.10			2.77	0.52	0.56	56000	103	310	1.18	
0.15	53.00			2.73	0.52	0.55	55400	108	307	1.07	
0.15	54.20			2.73	0.52	0.55	55900	110	306	1.04	
0.15	53.60			2.71	0.52	0.57	55200	111	314	1.11	
0.15	54.60			2.71	0.51	0.57	56800	109	308	1.03	
0.15	53.50			2.72	0.52	0.58	56500	115	310	1.11	
0.15	53.40			2.68	0.51	0.55	58000	114	312	1.11	
0.16	51.20			2.78	0.51	0.55	57500	111	311	1.23	
0.16				2.78	0.54	0.56	58300	112	310	1.20	
0.16				2.78	0.53	0.56	58900	110	310	1.20	
0.16				2.80	0.53	0.55	58800	112	310	1.20	
0.16				2.78	0.53	0.56	56700	107	310	1.20	
0.16				2.80	0.53	0.57	57400	111	310	1.20	
0.17				2.79	0.53	0.57	59600	107	320	1.20	
0.16				2.78	0.53		58800	118	320		
0.15					0.53				330		
0.15					0.52				309		
0.16					0.54				310		
0.15					0.54				320		
0.15					0.54				312		
0.15					0.55				307		
0.15					0.53				321		
0.15					0.53				308		
0.15					0.54				309		
0.14											
0.14											
0.15											
0.15											
0.15											
0.14											
0.16											
0.16											
0.16											
0.17											
0.16											
0.16											
0.16											
0.16											

Assay Data (Cont.)

4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	FUS	4A_MICP	4A_MICP	4A_MICP	FUS	4A_MICP
Ca	Co	Cr	Cs	Cu	Fe	Fe	Ga	Hf	K	K	Li
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
5900	37.00	233	2.60	62.00	21100	22181	13.50	2.80	13500	12000	9.20
6100	38.00	239	2.50	66.00	21900	22482	13.60	2.80	13900	12000	9.20
6100	39.00	245	2.50	64.00	21500	22281	13.40	2.70	13900	12000	9.20
6000	37.00	235	2.50	63.00	21700	22875	13.00	2.80	13600	12000	9.40
6000	37.00	235	2.50	63.00	21400	21996	13.90	2.80	13600	12000	9.20
6000	36.00	238	2.60	62.00	20900	22383	13.30	2.70	13600	12000	9.30
5800	36.00	229	2.50	61.00	21500	21585	13.50	2.90	13500	12000	9.40
5900	37.00	232	2.60	62.00	19800	20100	13.60	2.90	13500	12000	9.30
5675	35.15	243	2.29	65.30	19700	20000	13.85	2.70	13219	13000	9.00
5660	33.95	235	2.32	62.63	19700	19900	12.80	2.70	13214	13000	8.50
5577	33.55	239	2.41	62.47	19900	19600	13.10	2.70	13075	13000	8.80
5646	35.21	237	2.42	65.53	19600	19800	13.20	2.70	13250	13000	8.50
5685	35.15	241	2.37	67.14	20300	19900	12.80	2.80	13278	13000	9.30
5779	35.60	248	2.52	65.61	19800	20400	13.60	2.90	13555	13000	9.40
5734	34.33	250	2.55	64.01	20000	20400	13.95	2.90	13469	13000	9.90
5713	36.40	245	2.65	67.18	19600	19300	12.80	2.90	13457	13000	9.50
5600	36.40	222	2.50	66.00	19500	19800	13.60	3.20	12900		10.00
5500	36.70	225	2.60	62.00	19800	19900	12.80	3.20	13200		9.00
5600	37.00	230	2.40	63.00	19400	20100	13.40	3.20	13300		9.50
5500	36.20	221	2.50	64.00	19300	19900	13.20	3.20	13200		9.00
5500	36.00	223	2.50	63.00	20300	19900	13.00	3.20	13100		9.00
5700	36.80	230	2.40	64.00	20200	19500	12.80	3.20	13100		9.00
5700	37.10	231	2.30	63.00	20600	19900	13.40	2.80	13200		10.00
5800	35.00	238	2.60	65.00	21000				13300		
5500	35.00			59.30	21300						
5700	35.00			59.30	21300						
5600	35.00			58.30	20800						
5600	35.00			59.50	20600						
5300	35.00			59.30	21000						
5300	35.00			60.60	21200						
5500				62.70	21400						
				64.10							
				65.00							
				60.00							
				60.00							
				60.00							
				60.00							
				60.00							
				60.00							
				60.00							

Assay Data (Cont.)

4A_MICP	FUS	4A_MICP	FUS	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP
Mg	Mg	Mn	Mn	Mo	Na	Nb	Ni	P	Pb	Rb	Sb
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
2400	2611	565	539	2.00	1600	4.40	92.00	170	85.00	57.00	2.70
2500	2448	570	502	3.00	1600	4.80	98.00	150	91.00	56.50	2.60
2500	2631	575	556	2.00	1600	4.50	97.00	160	92.00	57.60	2.70
2400	2527	581	575	2.00	1600	4.50	94.00	170	91.00	57.30	2.70
2400	2390	564	517	2.00	1600	4.50	94.00	160	89.00	56.60	2.70
2400	2438	568	516	2.00	1600	4.90	93.00	160	90.00	56.90	2.50
2400	2616	557	468	2.00	1447	4.80	93.00	170	84.00	57.60	2.70
2400	2340	585	543	2.00	1407	4.30	94.00	170	89.00	57.30	2.90
2307	2340	575	465	3.50	1404	4.10	86.93	168	85.24	57.20	2.50
2307	2280	576	543	3.00	1441	4.20	85.01	170	84.53	56.00	2.44
2279	2220	568	543	2.90	1436	4.30	86.76	164	91.05	55.80	2.39
2306	2220	577	543	3.00	1411	4.30	86.44	169	87.85	58.00	2.73
2305	2220	580	543	3.10	1425	4.20	86.55	168	86.82	55.80	2.36
2347	2340	590	543	3.00	1439	4.40	85.24	166	86.71	58.20	2.73
2339	2280	582	543	3.00	1570	4.50	86.24	167	87.00	54.80	2.69
2334	2300	584	526	3.10	1511	4.90	86.79	165.6	87.00		2.72
2200	2300	528	585	2.86	1537	4.00	95.00	155	88.00		2.20
2200	2300	536	548	2.92	1561	4.00	95.00	147	88.00		2.20
2200	2300	538	554	2.84	1596	4.00	94.00	148	89.00		2.20
2200	2200	530	558	2.86	1601	4.00	95.00	148	88.00		2.40
2200	2300	528	554	3.42	1585	4.00	95.00	153	88.00		2.40
2300	2200	550	544	3.00	1586	4.00	95.00	156	88.00		2.60
2200	2300	544	564	3.43	1500		96.00	150	90.10		2.40
2300		561		3.47	1500		97.00	152	90.00		2.20
2400		542		2.50	1500		95.00	150	88.80		
2300		550		3.00	1500		90.00	140	86.30		
2400		552		3.00	1500		90.00	150	87.80		
2400		562		3.00	1500		90.00	150	92.20		
2300		546		3.00	1500		90.00	140	91.20		
2300		556		3.00	1600		90.00	160	88.00		
2400		560		3.00	1500		95.00	160	87.00		
		560		2.50	1600		90.00	150	89.00		
					1600				87.00		
					1600				86.00		
					1500				86.00		
					1600				88.00		
					1600				87.00		
					1600						

Assay Data (Cont.)

4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	FUS	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP	4A_MICP
Sc	Sr	Ta	Th	Ti	Ti	Ti	V	W	Y	Zn	Zr
ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm	ppm
6.00	56.00	0.61	23.02	1100	1680	0.51	43.00	0.70	8.60	153	89
6.00	59.00	0.70	23.12	1100	1680	0.52	45.00	0.70	8.90	153	90
6.00	60.00	0.70	23.32	1000	1680	0.52	45.00	0.60	8.90	153	86
6.00	57.00	0.63	23.30	1100	1620	0.52	44.00	0.60	8.30	147	89
6.00	59.00	0.64	23.22	1100	1620	0.52	44.00	0.70	8.30	152	88
6.00	58.00	0.71	23.79	1100	1680	0.52	44.00	0.70	8.70	141	88
6.00	57.00	0.75	23.33	1000	1680	0.54	43.00	0.60	8.90	144	86
6.00	58.00	0.62	23.45	1100	1680	0.52	44.00	0.60	9.00	141	87
6.00	47.92	0.56	26.00	1051	1600	0.50	38.64	0.60	8.40	143	93
6.00	48.38	0.55	25.00	1122	1600	0.50	37.57	0.60	8.30	146	97
6.00	48.54	0.54	24.40	1064	1600	0.50	36.95	0.60	8.50	142	95
6.00	49.58	0.51	24.60	986	1700	0.50	40.11	0.60	8.60	142	92
6.00	48.51	0.56	25.00	1092	1600	0.50	38.82	0.60	8.30	144	93
6.00	49.03	0.55	26.40	1147	1600	0.50	38.36	0.70	9.00	141	94
6.00	48.63	0.57	25.90	1072	1600	0.50	37.09	0.50	8.90	143	95
6.00	49.10	0.60	26.40	1010	1600	0.50	39.10	0.50	9.40	141	96
6.00	53.70	0.40	21.40	970			43.00	0.50	10.60	144	117
6.00	53.40	0.50	21.30	970			44.00	0.50	10.50	144	117
6.00	53.40	0.40	21.60	940			44.00	0.50	10.40	149	120
6.00	54.00	0.50	21.90	930			43.00	0.50	10.60	146	119
6.00	52.80	0.50	21.60	980			45.00	0.50	9.80	146	119
6.00	53.00	0.60	22.00	980			46.00	0.50	9.80	140	120
6.00	53.90	0.60	22.00	1030			45.00		10.00	140	115
	53.20	0.50	20.70	1000			44.00			141	
	52.50			1000			41.00			138	
	53.00			1100			41.00			139	
	52.00			1000			42.00			143	
	54.00			1000			41.00			145	
	51.00			1000			41.00			147	
	52.00			1000			43.00			140	
	51.00			1100			42.00			145	
	56.00						43.00			140	
							40.00			145	
							40.00			135	
							40.00			135	
							40.00			140	
							40.00			140	

26. Reported Values

The certified values listed in this certificate fulfil the AMIS statistical criteria (see section 19) regarding agreement for certification and have been independently validated by Allan Fraser.

27. Validation of Accuracy (Trueness)

This CRM can be used to validate accuracy (trueness) as required in method validation as stated in the ISO17025 standard. See Appendix 3 for an example on the validation of accuracy using replicate data derived from the analysis of a CRM.

28. Limit of Detection and Limit of Quantitation in Gravimetric Fire Assay

In the determination of limit of detection (LOD) and limit of quantitation (LOQ) in gravimetric fire assay (*i.e.* lead collection and weighing of resulting gold prill), the minimum mass that an assay microbalance is capable of weighing and the original test sample mass determines the LOD and the LOQ in the assay (Fraser, 2015), (see Appendix 8 for an example of the calculation LOD and LOQ and Table 13 for a recommend reporting scheme for LOD and LOQ values).

29. Metrological Traceability

The values quoted herein are based on the consensus values derived from statistical analysis of the data from an inter-laboratory measurement program. Traceability to SI units is via the standards used by the individual laboratories the majority of which are accredited to the ISO17025 general requirements for the competence of testing and calibration laboratories and who have maintained measurement traceability during the analytical process.

30. Period of Validity

The certified values are valid for this product, while still sealed in its original packaging, until notification to the contrary. The stability of the material will be subject to continuous testing for the duration of the inventory. Should product stability become an issue, all customers will be notified and notification to that effect will be placed on the www.amis.co.za website.

31. Minimum Sample Size

Most of the laboratories reporting used a 0.5g sample size for the ICP-OES and a 30g sample size for the fire assay. These are the recommended minimum sample sizes for the use of this material.

32. Availability

This product is available in Laboratory Packs containing 1kg of material and Explorer Packs containing custom weights (from 50 to 250g) of material. The Laboratory Packs are sealed bottles delivered in sealed foil pouches. The Explorer Packs contain material in standard geochem envelopes, nitrogen flushed, and vacuum sealed in foil pouches.

33. Recommended use in Quality Control

Users should set their own limits *i.e.* 1, 2 and 3 standard deviations from an obtained mean value based on at least 10 replicate analyses using this CRM (see Appendix 6 for detail on the use of this CRM in quality control).

34. Legal Notice

This certificate and the reference material described in it have been prepared with due care and attention. However, AMIS, Makhosi Khoza, and Allan Fraser; accept no liability for any decisions or actions taken following the use of the reference material.

Date of Version 000: 18 June 2019

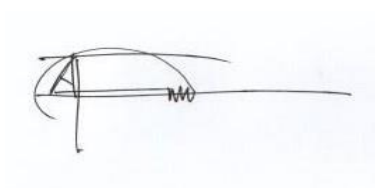
Version: 000

Approving Officer:

African Mineral Standards: _____

Makhosi Khoza (Quality Specialist)

Certifying Officer:

A handwritten signature in black ink, appearing to be 'A. Fraser', written over a horizontal line.

Geochemist: _____

Allan Fraser

M.Sc. (Geology), N.D. (Analytical Chem.),
Pr.Sci.Nat. Pr.Chem.SA

References

- Abzalov, M. (2011). Sampling Errors and Control of Assay Data Quality in Exploration and Mining Geology, Applications and Experiences of Quality Control. Ivanov, O., (Ed.), InTech, DOI: 10.5772/14965. Available from: <http://www.intechopen.com/books/applications-and-experiences-of-quality-control/sampling-errors-and-control-of-assay-data-quality-in-exploration-and-mining-geology>. Accessed 25th September, 2016.
- ASTM E122-09e1, Standard Practice for Calculating Sample Size to Estimate, With Specified Precision, the Average for a Characteristic of a Lot or Process, ASTM International, West Conshohocken, PA, 2011, www.astm.org
- Barwick V.J., Pritchard E. (Eds). (2011). Eurachem Guide; Terminology in Analytical Measurement – Introduction to VIM 3 (2011). ISBN 978-0-948926-29-7.
- Carr R.H. (2011). Estimating errors using graphs and taking good data. California State University Los Angeles. March 2011. 14. <http://web.calstatela.edu/faculty/kaniol/Err-Gph-Meas-IBooklet.pdf>. Retrieved: 29 August, 2016.
- Cochran, W.G. (1950). The Comparison of Percentages in Matched Samples. *Biometrika*, 37, 256-66.
- Ellison, S., Barwick, V., Duguid Farrant, T. (2009). *Practical statistics for the analytical scientist, a bench guide. 2nd Edition*. RSC Publishing. 25-172.
- ERM Application Note 1. (2005). Comparison of a measurement value with a certified value. European Commission – Joint Research Centre Institute for Reference Materials and Measurements (IRMM). 1-2.
- EURACHEM / CITAC Guide CG 4. (2012). Quantifying uncertainty in analytical measurement. 3rd Edition. Editors: Ellison, S. L. R., Williams, A. 4-121.
- Eurolab Technical Report No.1/2007. (2007). Measurement uncertainty revisited: Alternative approaches to uncertainty evaluation. European Federation of National Associations of Measurement, Testing and Laboratories. 43-46.
- Fraser, A.W. (2015). Minimising uncertainty in measurement and improving limit of detection in gold bearing ores from concentrations predicted by linear regression in atomic absorption spectrometry. M.Sc. thesis (unpublished), University of Johannesburg: 199 pages.
- Grubbs, F.E. (1969). Procedures for detecting outlying observations in samples. *Technometrics*, 11, 1969. 1-21.
- Horwitz, W., Albert, R. (2006). The Horwitz Ratio (HorRat): A useful index of method performance with respect to precision, *Journal of Association of Official Analytical Chemists International*, 89: 1095-1109.
- ISO/IEC 17025:2005(E). (2005). *General requirements for the competence of testing and calibration laboratories*. 2nd Edition: 14-59.
- ISO Guide 35 (2003). Certification of reference materials — General and Statistical principles, 3rd edition. ISO/REMCO WG 1.
- ISO 5725-2:1994. Accuracy (trueness and precision) of measurement methods and results -- Part 2: Basic method for the determination of repeatability and reproducibility of a standard measurement method.
- Long, J., Winefordner, J. (1983). Limit of detection - a closer look at the IUPAC definition. *Analytical Chemistry* 55: 712A - 724A.
- Miller, J., Miller, J. (2010). *Statistics for analytical chemistry*. 6th Edition. New York: Ellis Horwood. 36-126.
- Nelsen T.C., Wehling P. (2008). Collaborative studies for quantitative chemical analytical methods. AACC International Report. *Cereal Foods World*. September – October 2008, Vol. 53, No. 5. 285-288.
- Skoog, D., West, D. (1982). *Fundamentals of analytical chemistry*. 4th Edition. CBS College Publishing. Holt Saunders International Editions: 39-73.
- Thompson, M. (Ed.) (2008). Test for 'sufficient homogeneity' in a reference material. Analytical Methods Committee, AMCTB 17A, ISSN 1757-5958.
- Thompson, M. (Ed.) (2010). Internal quality control in routine analysis. AMC Technical Brief. Analytical Methods Committee. AMCTB No.46. 2010. Issn 1757-5958.
- Thompson, M., Lowthian, P. (2011). *Notes on statistics and data quality for analytical chemists*. Imperial College Press: 15-115.

Appendices

Appendix 1: Uncertified Element Statistics

Uncertified element statistics are shown in **Table 9**.

Table 9. Uncertified element concentrations statistics.

Element	Generic Method	n	Mean	SD	RSD %	Unit
Ag	4A_MICP	15	0.1	0.02	14	ppm
Al	FUS	24	58866.5	3615.4	6	ppm
Al ₂ O ₃	FUS	24	11.1	0.7	6	%
As	FUS	16	126.8	43.8	35	ppm
Ba	FUS	8	307.9	2.5	1	ppm
BaO	XRF	8	0.03	*	*	%
Bi	FUS	7	0.8	*	*	ppm
C	Combustion/LECO	38	1355.0	151.5	11	ppm
Ca	FUS	8	5000.0	*	*	ppm
CaO	FUS	16	0.7	0.1	11	%
Cd	4A_MICP	15	0.3	0.02	6	ppm
Cd	FUS	8	0.3	0.1	18	ppm
Ce	4A_MICP	16	86.6	4.6	5	ppm
Ce	FUS	7	79.4	0.9	1	ppm
Co	FUS	15	31.3	7.6	24	ppm
Cr	FUS	16	284.1	16.8	6	ppm
Cs	FUS	8	2.7	0.2	6	ppm
Cu	FUS	16	71.0	12.3	17	ppm
Dy	4A_MICP	8	2.4	0.1	4	ppm
Dy	FUS	7	3.4	0.1	2	ppm
Er	4A_MICP	8	1.3	0.1	7	ppm
Er	FUS	8	1.8	0.03	2	ppm
Eu	4A_MICP	7	0.9	*	*	ppm
Eu	FUS	8	1.0	0.03	3	ppm
Ga	FUS	8	12.9	0.8	6	ppm
Gd	4A_MICP	8	3.9	0.1	3	ppm
Gd	FUS	8	4.1	0.1	3	ppm
Ge	4A_MICP	16	0.4	0.3	72	ppm
Hf	FUS	8	3.6	0.5	14	ppm
Ho	4A_MICP	8	0.4	0.02	5	ppm
Ho	FUS	8	0.6	0.01	2	ppm
In	4A_MICP	16	0.03	0.01	26	ppm
La	4A_MICP	16	47.0	3.6	8	ppm
La	FUS	8	40.9	1.5	4	ppm
Li	FUS	3	10.7	0.6	5	ppm
LOI	LOI	47	2.9	0.2	7	%
Lu	4A_MICP	7	0.1	*	*	ppm
Lu	FUS	8	0.2	0.02	7	ppm
MnO	XRF	0	*	*	*	%
Mo	FUS	8	3.0	*	*	ppm
Nb	FUS	8	6.8	0.5	7	ppm
Nd	4A_MICP	8	33.0	1.0	3	ppm

Element	Generic Method	n	Mean	SD	RSD %	Unit
Nd	FUS	8	30.3	1.2	4	ppm
Ni	FUS	16	86.3	14.9	17	ppm
P	FUS	2	100.0	*	*	ppm
Pb	FUS	8	95.5	1.9	2	ppm
Pd	Pb Collection	14	0.01	*	*	g/t
Pr	4A_MICP	8	9.3	0.1	1	ppm
Pr	FUS	8	8.8	0.2	2	ppm
Pt	Pb Collection	15	0.01	0.002	22	g/t
Rb	FUS	8	60.1	4.2	7	ppm
Re	4A_MICP	2	0.003	0.001	28	ppm
S	FUS	8	0.6	0.02	3	%
Sb	FUS	8	2.5	0.3	14	ppm
Sc	4A_MICP	23	6.0	*	*	ppm
Sc	FUS	8	6.8	0.5	7	ppm
Se	4A_MICP	7	1.0	*	*	ppm
SiO ₂	CALC	8	78.0	0.1	0.1	%
SiO ₂	FUS	8	79.4	0.9	1	%
Sm	4A_MICP	8	5.7	0.1	3	ppm
Sm	FUS	8	5.0	0.1	2	ppm
Sn	4A_MICP	16	1.3	0.4	34	ppm
Sn	FUS	8	2.0	*	*	ppm
Sr	FUS	8	53.0	1.6	3	ppm
SrO	XRF	8	0.01	*	*	%
Ta	FUS	8	0.8	0.1	8	ppm
Tb	4A_MICP	8	0.5	0.02	4	ppm
Tb	FUS	8	0.6	0.03	5.1	ppm
Te	4A_MICP	12	0.1	0.02	20	ppm
Th	FUS	8	21.8	0.9	4	ppm
Tl	FUS	4	0.5	*	*	ppm
Tm	4A_MICP	8	0.2	0.01	6	ppm
Tm	FUS	8	0.3	0.01	5	ppm
U	FUS	8	53.6	1.5	3	ppm
U	XRF	16	58.1	6.2	11	ppm
U ₃ O ₈	XRF	15	0.01	0.001	10	%
V	FUS	8	44.3	0.7	2	ppm
V ₂ O ₅	XRF	3	0.01	*	*	%
Y	FUS	8	16.6	1.1	6	ppm
Yb	4A_MICP	8	1.2	0.1	5	ppm
Yb	FUS	8	1.6	0.1	3	ppm
Zn	FUS	15	131.1	5.6	4	ppm
Zr	FUS	8	138.4	6.7	5	ppm

* denotes that the results were too similar and SD and RSD% could not be calculated

Appendix 2 through 9, prepared by Allan Fraser.

Appendix 2. Certification of Reference Material and Estimation of Measurement Uncertainty

In the establishment of a consensus value for the CRM, outlier tests are carried out followed by performance statistics and the estimation of the measurement uncertainty. In practice, it is highly likely that data generated by multiple laboratories as an inter-laboratory comparison of material for certification, will contain erroneous as well as extreme measurements (outliers). The influence of outliers on summary statistics needs to be minimised by the application of procedures for outlier identification on raw data. The use of z-scoring, Cochran's test for suspect repeatability variances, along with Grubbs test for suspect measurement values allows for the detection of outliers (IUPAC, 1995). Method performance in terms of precision as relative standard deviation is judged by the application of the Horwitz ratio, which gives an indication of whether the observed relative standard deviation at the concentration levels of analyte determined are acceptable (Horwitz & Albert, 2006).

In the absence of an extensive uncertainty budget, measurement uncertainty is estimated from the reproducibility standard deviation from inter-laboratory data and reported as an expanded uncertainty at a level of confidence of 95% (Miller & Miller, 2010).

The steps below give detail on the establishment of a consensus value through the elimination of outliers, method performance and estimation of measurement uncertainty using standard uncertainties and the analysis of variance.

Z-Score

A z-score is calculated using equation [1]:

$$z = \frac{x - x_a}{s_p} \quad [1]$$

Where, x is the result of a submitted sample, x_a is the mean and s_p is the standard deviation of the submitted results from all the participating laboratories. Z-Scores are interpreted as follows:

$|z| \leq 2$ satisfactory performance
 $2 < |z| \leq 3$ questionable performance
 $|z| > 3$ unsatisfactory performance

(Thompson & Lowthian, 2011)

Data with z-scores exceeding two are discarded and are not included for further assessment.

Cochran's Test

The test of Cochran (1950) as shown in equation [2] is applied to any suspect repeatability variances:

$$C_{calc} = \frac{s_{max}^2}{\sum_{i=1}^l s_i^2} \quad [2]$$

Where, C_{calc} , s_{max}^2 and $\sum_{i=1}^l s_i^2$, are the calculated values for Cochran's test, data set with the maximum variance and the sum of the variances of all of the participating, l laboratory datasets. The C_{calc} value is compared with a critical value, C_{crit} at a level of confidence of 95% and an alpha of 0.05% (see Ellison, *et al.*, 2009, Appendix A, Table A.3a, page 209 for a table of critical values for the test of Cochran at LOC 95%).

According to ISO 5725-2 (1999), results from a laboratory with a suspect repeatability variance can be excluded if it is shown by the Cochran test to be an outlier. Therefore, if $C_{calc} > C_{crit}$, the laboratory with the maximum variance is removed. The data found to be excluded should not be $>2/9$, or 22% of the total data.

Grubbs Test

The test of Grubbs (1969) calculates a test statistic, G_{calc} and in the detection of a single outlier, G_1 is found by using

$$G_{1\,calc} = \frac{|Suspect\ value - \bar{x}|}{s} \quad [3]$$

Where, the sample mean and standard deviation, \bar{x} and s , are calculated with the suspect value included. The $G_{1\,calc}$ statistic is compared to a critical value for N measurements. See Ellison, *et al.*, 2009, Appendix A, Table A.2, page 208 for a table of critical values for the test of Cochran at LOC 95%.

Method Performance

The Horwitz function is used to assess the performance of the data under consideration, with respect to precision (Horwitz & Albert, 2006). A calculated %RSD is found using the Horwitz expression

$$\%RSD = \pm 2^{(1-0.5\log C)} \quad [4]$$

where, C is the analyte concentration in percent divided by 100 and \log is the natural logarithm. The observed %RSD is calculated as

$$Observed\ \%RSD = \frac{s}{Mean} \times 100 \quad [5]$$

where s is the standard deviation of n replicates.

The ratio of the observed %RSD and the calculated %RSD gives the Horwitz ratio (HorRat):

$$HorRat = \frac{\%RSD\ Observed}{\%RSD\ Calculated} \quad [6]$$

A HorRat <2 indicates that the method is of adequate precision. Should the HorRat be >2 the overall data are discarded, and the candidate material considered not suitable for certification as the precision is excessive for the concentration of the analyte being determined (Nelsen & Wehling, 2008).

Grand Mean

The grand mean ($\bar{\bar{x}}$) *i.e.* the certified value of a dataset is the total of all the data values divided by the total sample size (n):

$$\bar{\bar{x}} = \sum \frac{x}{n} \quad [7]$$

Certified Value

From ANOVA as per the description in section 19, an 'appropriate precision' as shown in [8] is calculated for sufficient homogeneity (Thompson, 2008):

$$s_r \leq 0.3u_c \quad [8]$$

Where, s_r is the within laboratory repeatability, as determined from [14]. Once [8] is satisfied, a grand mean [7] is calculated and this is taken to be the certified value.

Total Variation (SST)

The total variation (not the variance) comprises the sum of the squares of the differences of each mean with the grand mean.

$$SST = \sum (x - \bar{x})^2 \quad [9]$$

Between Group Variation (SSB)

The *variation* due to the between the laboratories is denoted SSB or Sum of Squares Between laboratories and given by [10]. If the laboratory means are close to each other (and therefore the Grand Mean) SSB will be a small value. There are P samples involved with one datum value for each sample (the sample mean), so there are P-1 degrees of freedom.

$$SSB = \sum n(\bar{x} - \bar{\bar{x}})^2 \quad [10]$$

The *variance* due to the interaction between the laboratories is denoted MSB for Mean Square Between groups and is the SSB divided by its degrees of freedom.

$$MS = \frac{SSB}{n - 1} \quad [11]$$

Within Group Variation (SSW)

The variation due to differences within individual samples is denoted SSW for Sum of Squares Within laboratories. The degrees of freedom are equal to the sum of the individual degrees of freedom for each sample. Since each sample has degrees of freedom (*df*) equal to one less than their sample sizes, and there are *k* samples, the total degrees of freedom is P less than the total sample size: $df = n - P$.

$$SSW = \sum df \cdot s^2 \quad [12]$$

The variance due to the differences within individual samples is denoted MSW for Mean Square Within groups. This is the within group variation divided by its degrees of freedom:

$$MSW = \frac{SSW}{P - n} \quad [13]$$

From equations [9] through [13], the ANOVA table as shown in Table 10 is developed.

Table 10. A single-factor ANOVA table showing key elements. Where P is the total number of groups, or laboratories. P-1 is 1 less than number of laboratories, P (n-1) is the number of data values minus number of groups (equals degrees of freedom for each group added together), and P-1 + P(n-1) is 1 less than the number of data points. MS is the mean squares of between laboratories and within laboratories. After Ellison *et al.*, (2009), Table 6.2, page 61.

Source	Sum of Squares	df	Mean Sum of Squares	F	p	F _{crit}
Between Laboratories	SSB	P-1	MSB=SSB/df	MSB/MSW	=FDIST(x,df,df)	F-table
Within Laboratories	SSW	P(n-1)	MSW=SSW/df	–	–	–
Total	SSB+SSW	P-1 + P(n-1)	–	–	–	–

Combined Standard Uncertainty

The combined standard uncertainty (u_c) represents the effects of random events such as days, instruments, and analysts on the precision of the analytical procedures of all accepted data of the participating laboratories. Using the output from ANOVA, the combined standard uncertainty (u_c) is determined from the square root of the sum of squares of the variances of the within laboratory repeatability, s_r , and the between laboratory precision, s_s :

$$u_c = \sqrt{s_r^2 + s_s^2} \quad [14]$$

Within laboratory repeatability is determined as

$$s_r = \sqrt{MSB} \quad [15]$$

and, the between laboratory precision as

$$s_s = \sqrt{\frac{(MSW - MSB)}{n}} \quad [16]$$

where MSW is the mean squares of the within laboratory variance, MSB is the mean squares for the between laboratories and n in this case, is the number of replicates in a group of the accepted data (Thompson & Lowthian, 2011).

Expanded Uncertainty

The expanded uncertainty (U) at a confidence level of 95% is determined by multiplication of the combined uncertainty (u_c) by a coverage factor (k) found from $N-1$ degrees of freedom (df), where N is the number of laboratory means accepted in the establishment of the certified value. The t-critical value for 5% significance can be found in a t-critical table (see Appendix 9, or from MS Excel as =TINV (5%, df)).

Uncertainty Statement

Typically, an uncertainty statement is presented as follows: Au =0.77±0.04 g/t, where the number following the symbol ± is the numerical value of an expanded uncertainty, $U = ku_c$, with U determined from a combined standard uncertainty multiplied by a coverage factor $k = 2$ or, a t-critical value for $N-1$ accepted laboratories. Since it can be assumed that the possible estimated values of the standard are approximately normally distributed with standard uncertainty, u_c , the certified value of the CRM is believed to lie in the interval defined by U with a level of confidence of approximately 95 %, e.g. a mean value of 0.77±0.04g/t will have intervals of: 0.73≤0.77≤0.81 g/t.

Appendix 3. Example: Comparison of Mean and Certified Value for Validation of Accuracy

According to ERM (2005); Eurolab (2007); Abzalov (2011) and Carr (2011), the validation of accuracy for a given mean and certified value requires the inclusion of the measurement uncertainty of the CRM in a t-test for statistical significance. The classical Student's t-test as shown in [17], does not consider the measurement uncertainty of the CRM. To compensate for this, Eurolab Technical Report No.1/2007 recommends equation [18] for the validation of CRMs with stated measurement uncertainties.

$$t_{calc} = \frac{|\bar{x} - \mu|}{\frac{s}{\sqrt{n}}} \quad [17]$$

$$t_{calc} = \frac{|\bar{x} - \mu|}{\sqrt{(u_{\mu})^2 + \frac{s^2}{n}}} \quad [18]$$

Where, t_{calc} is the calculated t-statistic, \bar{x} the mean of n replicates with a standard deviation of s for a CRM of μ certified value. The standard uncertainty u is the stated expanded uncertainty (U) of the CRM divided by the coverage factor (k) as stated on the certificate of analysis. Note that the $| \quad |$ bars indicate that the absolute value between the mean and the certified value is to be used, *i.e.* ignore the sign.

An example in which [18] is used for validation of accuracy is given below.

Example

A CRM is independently replicated nine times for Al_2O_3 concentration by XRF analysis, *i.e.* 9 individual fused glass beads were prepared. The observed mean and standard deviation of the replicate data are shown with the certified value and expanded uncertainty in Table 11. In validation of accuracy, the hypothesis question is: Is the difference between the observed mean and the certified value statistically significant at a level of confidence of 95%? Alternatively put, is there sufficient evidence to conclude that the data *i.e.* replicates generated, are inaccurate?

The relevant hypotheses are:

Null hypothesis: H_0 : Mean = Certified value of CRM with stated measurement uncertainty. The acceptance of H_0 means that accuracy is demonstrated; *i.e.* insufficient evidence to reject H_0 ;

Alternate hypothesis: H_1 : Mean \neq Certified value of CRM with stated measurement uncertainty. The acceptance of H_1 means that accuracy is not demonstrated, *i.e.* there is sufficient evidence to accept H_1 ;

Table 11. CRM certified value quoted expanded uncertainty U , the coverage factor for the CRM, $k=2.25$ and mean for $n=9$ replicates and corresponding standard deviation for the replicate data.

CRM Certified Value	Expanded % (U)	Coverage Factor (k)	Mean ($n=9$)	n	Standard Deviation (s)
4.62%	0.08%	2.25	4.59%	9	0.01015

The standard uncertainty (u) is found by dividing the expanded uncertainty by the coverage factor:

$$u = \frac{0.08}{2.25} = 0.0356 \%$$

Using the observed mean for the replicate data ($n=9$) obtained for the CRM and substituting into [18]:

$$t_{calc} = \frac{|\bar{x} - \mu|}{\sqrt{0.0356^2 + \frac{0.01015^2}{9}}} = \frac{|4.59 - 4.62|}{\sqrt{0.00126 + 0.00001145}} = 0.84$$

Therefore, $t_{calc} = 0.84$ and $t_{crit}(5\%, 8) = 2.31$ (df is 8, therefore, $t_{crit}=2.31$, see Appendix 9, page 34) which is >0.84 . Similarly, the p -value= 0.43 which is >0.05 . This is strong evidence in favour of accepting the null hypothesis that there is no significant statistical difference between the certified value and the observed mean. Therefore, under the conditions that the uncertainty associated with the certified value is known, the accuracy is validated for the CRM tested. If the null hypothesis is accepted that the mean obtained is not statistically different from the certified value, then the principle of traceability has been conformed to.

Appendix 4. Two-standard Deviations

Two-standard deviations are calculated using the expression:

$$\text{Two standard deviations} = 2 (u_c) \quad [19]$$

Where, u_c is the standard combined uncertainty (equation [14]).

Appendix 5. Confidence Interval

Confidence interval is calculated as:

$$\text{Confidence Interval (CI)} = \frac{(t_{critical})s}{\sqrt{N}} \quad [20]$$

Where, N is the number of laboratories (accepted laboratory data), $t_{critical}$ is a two-tailed value for $N - 1$ degrees of freedom (df) and s , is the standard deviation of the accepted laboratory means. A two-tailed critical value is found for $N - 1$ degrees of freedom from either a t -distribution table (Appendix 9) or MS Excel as =TINV (5%, df).

Appendix 6. Using the CRM in Quality Control

QC chart control limits should not be determined by the certified value and stated measurement uncertainty of the certified reference material used. These parameters although “certified” will never be known; it is only the corresponding statistical estimates, *i.e.* standard deviation and the mean calculated from replicated results that are known and these should be used in quality control charts. However, should the laboratory choose to use the certified value as the mean then the quoted $2s$, or CI value for the CRM can be used in the quality control chart.

It is recommended that a Shewhart chart be developed for the use if this CRM is to be used as a control sample in laboratory quality control. A Shewhart chart is a plot of sequential assay results obtained from quality control material such as an AMIS CRM. The warning and control limits are based on the standard deviation obtained from the mean of the replicates of a CRM (Ellison, *et al.*, 2009; Thompson, 2010). The procedure in preparing a Shewhart chart is as follows:

1. Analyse 10 to 15 replicates or more of the AMIS CRM.
2. Apply the Grubbs test for outliers.
3. Determine the mean of the replicates after application of the Grubbs test.
4. Determine the standard deviation, using equation [21], of the replicates following Grubbs test.
5. Calculate the standard deviation, s from:

$$s = \sqrt{\frac{\sum(x_i - \bar{x})^2}{n - 1}} \quad [21]$$

where, x_i is an individual measurement in the data set, \bar{x} is the mean of the data set at $n-1$ degrees of freedom (df) and n is the number of replicates. The sample standard deviation can be found using the MS Excel formula “=stdev.s (number1;)”.

6. Verify accuracy of the mean value using equation [18].
7. Once accuracy is verified, calculate $\pm 2s$ and $\pm 3s$, where s is the standard deviation calculated from [21].
8. Construct the Shewhart control chart around the mean of n replicates.
9. Use $\pm 2s$ as the warning limits.
10. Use $\pm 3s$ as the control limits.
11. It is recommended that if 2 to 3 points are outside the warning limits analyse another sample and if it is then within warning limits, continue. If it is outside the warning limits, stop and troubleshoot.
12. It is recommended that if any point is outside control limits, analyse another portion (sample) of the CRM. If it is within control limits, continue. If it is outside control limits, stop and troubleshoot.
13. For reference purposes, the CRM certified value can be plotted on the Shewhart chart alongside the mean value.

On a regular basis the accuracy of the replicates of the CRM should be assessed in terms of the certified value of the CRM using equation [18].

Appendix 7. Conversion to Air-dry Basis (Prepared by Allan Fraser)

Since AMIS certified analyte values are reported on a dry-basis, the user laboratory is required to dry a portion (accurately weigh out 1.0 grams in duplicate) of the CRM material in air at 105°C in a drying oven to constant mass to determine the moisture content. Use a crucible with a flat inner surface with a surface area not smaller than 10 cm² with the CRM material spread evenly over same; this represents a 0.1 gram spread per cm². In correcting the certified value for moisture content, a moisture correction factor is calculated:

$$\text{Moisture correction factor (MCF)} = \frac{100 - \% \text{Moisture at } 105^{\circ}C}{100} \quad [22]$$

$$\text{Air dry basis concentration} = \text{MCF} \times \text{certified value on a dry basis} \quad [23]$$

Example

The moisture content determined at 105°C on a CRM is 0.500%. The certified analyte concentration for the CRM is 12.62±0.52% (dry basis). Calculating the moisture correction factor using [22] gives:

$$\text{Moisture correction factor} = \frac{100 - 0.500}{100} = 0.995$$

Multiplying the factor of 0.995 by the certified value as stated on the certificate of analysis on a dry basis (as in [23]) gives the analyte concentration on an air-dry basis:

$$0.995 \times 12.62\% = 12.56\%$$

The stated measurement uncertainty also needs to be corrected using [22] and [23], e.g. $0.995 \times 0.52 = 0.51_{(7)}$, rounded to 0.52%. The air-dry basis concentration *i.e.* $12.56 \pm 0.52\%$ is to be used as the certified value with its corresponding measurement of uncertainty.

Appendix 8. Example of Determination of LOD and LOQ in Fire Assay

The limit of detection (LOD) is the minimum detectable quantity of the analyte of interest (Skoog & West, 1985). To determine the LOD in fire assay by lead collection, the minimum mass that an assay microbalance is capable of weighing (m in micrograms, and the original test sample mass, $Mass_{assay}$ in grams) determines the LOD. The smallest prill mass most assay microbalances can measure is $1\mu\text{g}$ or 0.001mg . Even with a microscope it may be difficult to locate and pick up a prill weighing ten times that amount (*i.e.* 0.01mg or $10\mu\text{g}$) and weigh it. If an analyst can weigh a prill of $1\mu\text{g}$ then the LOD becomes $1\mu\text{g}$. However, the concentration factor would be 50 times for a 50-gram assay sample and therefore the LOD in g/t becomes $1\mu\text{g}$ divided by the original mass of the sample in grams taken for fire assay [24]. Therefore, the LOD in fire assay is computed as:

$$LOD = \frac{m (\mu\text{g})}{Mass_{assay} (g)} (\text{g/t}) \quad [24]$$

The limit of quantitation (LOQ), is simply the LOD multiplied by 10 (Long & Winefordner, 1983):

$$LOQ = 10 \cdot \frac{m (\mu\text{g})}{Mass_{assay} (g)} (\text{g/t}) \quad [25]$$

Therefore, with a sample mass of 50g taken for fire assay, the limit of detection would be 0.02g/t . *i.e.* $1\mu\text{g} = 1\text{g/t}$, therefore $1\mu\text{g}/50\text{g} = 0.02\text{g/t}$. If no prill was found, then the LOD result would be $<0.02\text{g/t}$ or "not detected". Using a larger assay sample mass improves the LOD and LOQ (Table 11). Table 13 gives a recommended reporting scheme for LOD and LOQ.

Table 12. Mass of assay sample and corresponding limit of detection and limit of quantitation for an assay microbalance capability of smallest prill mass of $1\mu\text{g}$ or 0.001mg .

Mass Assay Sample (g)	LOD (g/t)	LOQ (g/t)
30	0.03	0.3
50	0.02	0.2
100	0.01	0.1

Table 13. Recommended reporting scheme for LOD and LOQ in fire assay.

Data	Report as
<LOD	Not detected
<LOQ	Detected
\geq LOQ	Report assay result

Appendix 9. T-distribution table

Table 14. T-distribution table for t-critical values (t crit.) for a two-tailed t-test at a 95% level of confidence.

df	Two-tailed	df	Two-tailed
1	12.71	23	2.06
2	4.30	24	2.06
3	3.18	25	2.06
4	2.78	26	2.05
5	2.57	27	2.05
6	2.44	28	2.04
7	2.36	29	2.04
8	2.30	30	2.04
9	2.26	35	2.03
10	2.22	40	2.02
11	2.20	45	2.01
12	2.17	50	2.00
13	2.16	55	2.00
14	2.14	60	2.00
15	2.13	70	1.99
16	2.12	80	1.98
17	2.11	90	1.98
18	2.10	100	1.98
19	2.09	120	1.98
20	2.08	Infinity	1.96
21	2.08		
22	2.07		

End of certificate