



AMIS_Documents

Doc: ADOC_074

Originator:

Quality
Specialist

Approver:
Management Rep

Revision No: 011

Revision Date:22.01.2021

Issued By: Management Rep

Certificate

AMIS0770

Certified Reference Material

PGM UG2, South Africa

Certificate of Analysis

AMIS

A: 11 Avalon Road, West Lake View, Ext 11, Modderfontein, 1609, South Africa

P: PO Box 856, Isando, 1600, Gauteng, South Africa

T: +27 (0) 11 923 0800

W: www.amis.co.za



Table of Contents

| | |
|--|----|
| Table of Contents | 2 |
| List of Tables | 3 |
| List of Appendices | 3 |
| Summary Statistics | 4 |
| 1. Certified Concentrations and Uncertainties | 6 |
| 2. Statistical Comparison of Means | 8 |
| 3. Intended Use | 9 |
| 4. Abbreviations and Symbols | 10 |
| 5. Uncertified Concentration Values | 11 |
| 6. Units | 11 |
| 7. Analytical and Physical Methods | 11 |
| 8. Origin of Material | 11 |
| 9. Approximate Mineral and Chemical Composition | 11 |
| 10. Quantitative Analysis by X-Ray Diffraction | 12 |
| 11. Health and Safety | 13 |
| 12. Method of Preparation | 13 |
| 13. Particle Size Determination | 13 |
| 14. Handling | 14 |
| 15. Storage information | 14 |
| 16. Methods of Analysis Requested | 14 |
| 17. Information Requested of Participating Laboratories | 14 |
| 18. Certification of Mean and Estimation of Measurement Uncertainty | 15 |
| 19. Two Standard Deviations | 15 |
| 20. Confidence Interval | 15 |
| 21. Expanded Uncertainty | 15 |
| 22. Confidence Interval and Expanded Uncertainty | 15 |
| 23. Participating Laboratories | 16 |
| 24. Accepted Assay Data | 17 |
| 25. Reported Values | 23 |
| 26. Validation of Accuracy (Trueness) | 23 |
| 27. Limit of Detection and Limit of Quantitation in Gravimetric Fire Assay | 23 |
| 28. Metrological Traceability | 23 |
| 29. Period of Validity | 23 |
| 30. Minimum Sample Size | 23 |
| 31. Availability | 23 |
| 32. Recommended use in Quality Control | 23 |
| 33. Legal Notice | 24 |
| References | 25 |
| Appendices | 26 |
| End of certificate | 36 |

List of Tables

| | |
|--|----|
| Table 1. Certified concentrations, two standard deviations, combined and expanded uncertainty. | 6 |
| Table 2. Certified major oxides concentrations, two standard deviations, combined and expanded uncertainty..... | 7 |
| Table 3. The results of a two-sample equal or unequal variance t-test (two-tailed) data sets in which different analytical methods /instrumentation were used..... | 8 |
| Table 4. Abbreviations, symbols and descriptions. | 10 |
| Table 5. Results of XRD analysis. | 12 |
| Table 6. Particle Size Determination by laser diffraction..... | 13 |
| Table 7. Example of replicate assay data in which the <i>CI</i> , <i>2s</i> and <i>U</i> are shown..... | 16 |
| Table 8. Data used to calculate the certified values after removal of outliers. | 17 |
| Table 9. Uncertified element concentrations statistics. | 26 |
| Table 10. A single-factor ANOVA table showing key elements. Where <i>P</i> is the total number of groups, or laboratories. <i>P</i> -1 is 1 less than number of laboratories, <i>P</i> (<i>n</i> -1) is the number of data values minus number of groups (equals degrees of freedom for each group added together), and <i>P</i> -1 + <i>P</i> (<i>n</i> -1) is 1 less than the number of data points. <i>MS</i> is the mean squares of between laboratories and within laboratories. After Ellison <i>et al.</i> , (2009), Table 6.2, page 61..... | 31 |
| Table 11. CRM certified value, quoted expanded uncertainty <i>U</i> , the coverage factor for the CRM, <i>k</i> =2.25 and mean for <i>n</i> =9 replicates and corresponding standard deviation for the replicate data. | 32 |
| Table 12. Mass of assay sample and corresponding limit of detection and limit of quantitation for an assay microbalance capability of smallest prill mass of 1µg or 0.001mg. | 35 |
| Table 13. Recommended reporting scheme for LOD and LOQ in fire assay..... | 35 |
| Table 14. T-distribution table for t-critical values (t crit.) for a two-tailed t-test at a 95% level of confidence..... | 36 |

List of Appendices

| | |
|---|----|
| Appendix 1. Uncertified Element Statistics..... | 26 |
| Appendix 2. Certification of Reference Material and Estimation of Measurement Uncertainty..... | 28 |
| Appendix 3. Example: Comparison of Mean and Certified Value for Validation of Accuracy | 32 |
| Appendix 4. Two-standard Deviations..... | 33 |
| Appendix 5. Confidence Interval..... | 33 |
| Appendix 6. Using the CRM in Quality Control..... | 33 |
| Appendix 7. Conversion to Air-dry Basis..... | 34 |
| Appendix 8. Example of Determination of LOD and LOQ in Fire Assay | 35 |
| Appendix 9. T-distribution table | 36 |

Summary Statistics

Recommended Concentrations and Limits (at two Standard Deviations) Certified Concentrations

| Analyte | Method | ⁷ (μ) | ⁹ (2s) ± | Unit |
|---------|----------------------------|------------------|---------------------|------|
| Au | ¹ NiS | 0.038 | 0.011 | g/t |
| Au | ² Pb Collection | 0.041 | 0.014 | g/t |
| Ir | NiS | 0.116 | 0.042 | g/t |
| Pd | NiS | 0.774 | 0.14 | g/t |
| Pd | Pb Collection | 0.770 | 0.13 | g/t |
| Pt | NiS | 1.49 | 0.32 | g/t |
| Pt | Pb Collection | 1.46 | 0.15 | g/t |
| Rh | NiS | 0.260 | 0.043 | g/t |
| Ru | NiS | 0.479 | 0.093 | g/t |
| Co | ³ 2A_MICP | 15 | 2 | ppm |
| Co | ⁴ 4A_MICP | 150 | 10 | ppm |
| Co | ⁵ FUS | 154 | 17 | ppm |
| Cu | 2A_MICP | 117 | 15 | ppm |
| Cu | 4A_MICP | 112 | 10 | ppm |
| Ni | 2A_MICP | 279 | 50 | ppm |
| Ni | 4A_MICP | 865 | 58 | ppm |
| Ni | FUS | 873 | 92 | ppm |
| S | 4A_MICP | 0.050 | 0.0009 | % |
| S | XRF | 0.052 | 0.006 | % |
| Al | 4A_MICP | 7.91 | 0.28 | % |
| Al | FUS | 7.59 | 0.32 | % |
| Al | ⁶ XRF | 7.70 | 0.046 | % |
| As | 4A_MICP | 2 | 0.4 | ppm |
| Ba | 4A_MICP | 57 | 6 | ppm |
| Be | 4A_MICP | 0.5 | 0.05 | ppm |
| Bi | 4A_MICP | 0.1 | 0.03 | ppm |
| Ca | 4A_MICP | 4.05 | 0.12 | % |
| Ca | FUS | 4.16 | 0.22 | % |
| Ca | XRF | 4.17 | 0.018 | % |
| Ce | 4A_MICP | 6 | 0.2 | ppm |
| Cr | XRF | 10.20 | 0.16 | % |
| Cs | 4A_MICP | 0.3 | 0.06 | ppm |
| Fe | 4A_MICP | 11.52 | 0.59 | % |
| Fe | FUS | 11.33 | 1.2 | % |
| Fe | XRF | 11.61 | 0.22 | % |
| Ga | 4A_MICP | 28 | 4 | ppm |
| Hf | 4A_MICP | 0.6 | 0.04 | ppm |
| K | 4A_MICP | 1403 | 235 | ppm |
| Li | 4A_MICP | 7 | 2 | ppm |
| Mg | 4A_MICP | 7.67 | 0.66 | % |
| Mg | FUS | 7.57 | 0.47 | % |
| Mg | XRF | 7.68 | 0.059 | % |
| Mn | 4A_MICP | 1469 | 49 | ppm |
| Mn | FUS | 1585 | 136 | ppm |
| Na | 4A_MICP | 5946 | 696 | ppm |
| Na | XRF | 5884 | 347 | ppm |
| Nb | 4A_MICP | 141 | 5 | ppm |
| P | 4A_MICP | 112 | 7 | ppm |
| P | XRF | 98 | 30 | ppm |
| Pb | 4A_MICP | 7 | 0.8 | ppm |
| Rb | 4A_MICP | 5 | 0.9 | ppm |
| Sb | 4A_MICP | 2 | 0.3 | ppm |
| Sc | 4A_MICP | 18 | 3 | ppm |
| Si | XRF | 15.65 | 0.16 | % |
| Sr | 4A_MICP | 130 | 29 | ppm |
| Ta | 4A_MICP | 150 | 4 | ppm |
| Th | 4A_MICP | 0.5 | 0.06 | ppm |
| Ti | 4A_MICP | 3281 | 98 | ppm |
| Ti | FUS | 3561 | 335 | ppm |
| U | 4A_MICP | 0.5 | 0.1 | ppm |
| V | 4A_MICP | 857 | 69 | ppm |

| Analyte | Method | ⁷ (μ) | ⁹ (2s) ± | Unit |
|---------|---------|------------------|---------------------|------|
| W | 4A_MICP | 2 | 0.2 | ppm |
| Y | 4A_MICP | 5 | 0.8 | ppm |
| Zn | 4A_MICP | 305 | 5 | ppm |
| Zn | FUS | 330 | 14 | ppm |
| Zr | 4A_MICP | 15 | 3 | ppm |

Calculation 3E = Platinum (Pb Collection) + Palladium (Pb Collection) + Gold (Pb Collection) = 2.271g/t

Calculation 6E = Platinum (NiS) + Palladium (NiS) + Gold (NiS) + Iridium (NiS)+ Ruthenium (NiS) + Rhodium (NiS) = 3.157g/t

Major Oxides Certified Concentrations (at two Standard Deviations)

| Analyte | Method | ⁷ (μ) | ⁹ (2s) ± | Unit |
|--------------------------------|---------|------------------|---------------------|------|
| Al ₂ O ₃ | XRF | 14.54 | 0.073 | % |
| CaO | XRF | 5.84 | 0.029 | % |
| Cr ₂ O ₃ | XRF | 14.91 | 0.23 | % |
| Fe ₂ O ₃ | FUS | 16.64 | 0.80 | % |
| Fe ₂ O ₃ | XRF | 16.60 | 0.32 | % |
| MgO | XRF | 12.77 | 0.056 | % |
| Na ₂ O | 4A_MICP | 0.80 | 0.08 | % |
| Na ₂ O | XRF | 0.79 | 0.05 | % |
| SiO ₂ | XRF | 33.39 | 0.14 | % |
| SO ₃ | XRF | 0.13 | 0.01 | % |

1. Certified Concentrations and Uncertainties

AMIS0770 is a new standard material, developed and certified in July 2021. Table 1 gives the certified concentrations, confidence interval, combined and expanded uncertainty for the certified reference material. Table 2 shows the certified major oxides concentrations, two standard deviations, confidence interval, combined and expanded uncertainty.

Table 1. Certified concentrations, two standard deviations, combined and expanded uncertainty.

| Analyte | Method | ⁷ (μ) | N | n | k | % RSD | ⁸ (u_c) | ⁹ (2s) \pm | ¹⁰ (CI) 95% | ¹¹ (U) \pm | Unit |
|---------|----------------------------|------------------------|---|----|--------|-------|------------------------|-------------------------|------------------------|-------------------------|------|
| Au | ¹ NiS | 0.038 | 4 | 31 | 3.182 | 14 | 0.0053 | 0.011 | 0.0080 | 0.02 | g/t |
| Au | ² Pb Collection | 0.041 | 4 | 30 | 3.182 | 17 | 0.0071 | 0.014 | 0.010 | 0.02 | g/t |
| Ir | NiS | 0.116 | 4 | 29 | 3.182 | 18 | 0.021 | 0.042 | 0.037 | 0.07 | g/t |
| Pd | NiS | 0.774 | 5 | 38 | 2.776 | 9 | 0.068 | 0.14 | 0.085 | 0.2 | g/t |
| Pd | Pb Collection | 0.770 | 3 | 23 | 4.303 | 9 | 0.067 | 0.13 | 0.15 | 0.3 | g/t |
| Pt | NiS | 1.49 | 5 | 38 | 2.776 | 11 | 0.16 | 0.32 | 0.19 | 0.4 | g/t |
| Pt | Pb Collection | 1.46 | 3 | 23 | 4.303 | 5 | 0.077 | 0.15 | 0.097 | 0.3 | g/t |
| Rh | NiS | 0.260 | 4 | 31 | 3.182 | 8 | 0.021 | 0.043 | 0.031 | 0.07 | g/t |
| Ru | NiS | 0.479 | 5 | 36 | 2.776 | 10 | 0.046 | 0.093 | 0.055 | 0.1 | g/t |
| Co | ³ 2A_MICP | 15 | 3 | 24 | 4.303 | 8 | 1 | 2 | 3 | 5 | ppm |
| Co | ⁴ 4A_MICP | 150 | 2 | 16 | 12.706 | 3 | 5 | 10 | 43 | 63 | ppm |
| Co | ⁵ FUS | 154 | 6 | 46 | 2.571 | 5 | 8 | 17 | 7 | 22 | ppm |
| Cu | 2A_MICP | 117 | 2 | 16 | 12.706 | 7 | 8 | 15 | 68 | 98 | ppm |
| Cu | 4A_MICP | 112 | 3 | 23 | 4.303 | 5 | 5 | 10 | 12 | 22 | ppm |
| Ni | 2A_MICP | 279 | 2 | 16 | 12.706 | 9 | 25 | 50 | 222 | 316 | ppm |
| Ni | 4A_MICP | 865 | 3 | 23 | 4.303 | 3 | 29 | 58 | 68 | 125 | ppm |
| Ni | FUS | 873 | 5 | 40 | 2.776 | 5 | 46 | 92 | 56 | 128 | ppm |
| S | 4A_MICP | 0.050 | 2 | 16 | 12.706 | 0.9 | 0.0004 | 0.0009 | 0.003 | 0.006 | % |
| S | XRF | 0.052 | 3 | 22 | 4.303 | 6 | 0.003 | 0.006 | 0.007 | 0.01 | % |
| Al | 4A_MICP | 7.91 | 3 | 22 | 4.303 | 2 | 0.14 | 0.28 | 0.32 | 0.6 | % |
| Al | FUS | 7.59 | 3 | 24 | 4.303 | 2 | 0.16 | 0.32 | 0.34 | 0.7 | % |
| Al | ⁶ XRF | 7.70 | 3 | 22 | 4.303 | 0.3 | 0.023 | 0.046 | 0.037 | 0.1 | % |
| As | 4A_MICP | 2 | 2 | 16 | 12.706 | 10 | 0.2 | 0.4 | 2 | 3 | ppm |
| Ba | 4A_MICP | 57 | 3 | 24 | 4.303 | 5 | 3 | 6 | 8 | 13 | ppm |
| Be | 4A_MICP | 0.5 | 3 | 24 | 4.303 | 5 | 0.02 | 0.05 | 0.02 | 0.1 | ppm |
| Bi | 4A_MICP | 0.1 | 3 | 23 | 4.303 | 14 | 0.02 | 0.03 | 0.03 | 0.07 | ppm |
| Ca | 4A_MICP | 4.05 | 3 | 24 | 4.303 | 1 | 0.058 | 0.12 | 0.086 | 0.3 | % |
| Ca | FUS | 4.16 | 3 | 22 | 4.303 | 3 | 0.11 | 0.22 | 0.16 | 0.5 | % |
| Ca | XRF | 4.17 | 2 | 16 | 12.706 | 0.2 | 0.0091 | 0.018 | 0.068 | 0.1 | % |
| Ce | 4A_MICP | 6 | 2 | 16 | 12.706 | 1 | 0.08 | 0.2 | 0.3 | 1 | ppm |
| Cr | XRF | 10.20 | 3 | 23 | 4.303 | 0.8 | 0.078 | 0.16 | 0.12 | 0.3 | % |
| Cs | 4A_MICP | 0.3 | 3 | 24 | 4.303 | 9 | 0.03 | 0.06 | 0.07 | 0.1 | ppm |
| Fe | 4A_MICP | 11.52 | 3 | 24 | 4.303 | 3 | 0.30 | 0.59 | 0.71 | 1 | % |
| Fe | FUS | 11.33 | 5 | 36 | 2.776 | 5 | 0.60 | 1.2 | 0.89 | 2 | % |
| Fe | XRF | 11.61 | 2 | 16 | 12.706 | 1 | 0.11 | 0.22 | 0.98 | 1 | % |
| Ga | 4A_MICP | 28 | 3 | 24 | 4.303 | 7 | 2 | 4 | 4 | 8 | ppm |
| Hf | 4A_MICP | 0.6 | 3 | 22 | 4.303 | 3 | 0.02 | 0.04 | 0.04 | 0.09 | ppm |
| K | 4A_MICP | 1403 | 4 | 31 | 3.182 | 8 | 118 | 235 | 186 | 374 | ppm |
| Li | 4A_MICP | 7 | 3 | 24 | 4.303 | 10 | 0.8 | 2 | 2 | 3 | ppm |
| Mg | 4A_MICP | 7.67 | 4 | 30 | 3.182 | 4 | 0.33 | 0.66 | 0.55 | 1 | % |
| Mg | FUS | 7.57 | 3 | 23 | 4.303 | 3 | 0.24 | 0.47 | 0.55 | 1 | % |
| Mg | XRF | 7.68 | 2 | 16 | 12.706 | 0.4 | 0.029 | 0.059 | 0.23 | 0.4 | % |
| Mn | 4A_MICP | 1469 | 3 | 24 | 4.303 | 2 | 24 | 49 | 55 | 105 | ppm |
| Mn | FUS | 1585 | 4 | 32 | 3.182 | 4 | 68 | 136 | 95 | 217 | ppm |
| Na | 4A_MICP | 5946 | 5 | 38 | 2.776 | 6 | 348 | 696 | 433 | 967 | ppm |
| Na | XRF | 5884 | 3 | 23 | 4.303 | 3 | 174 | 347 | 379 | 747 | ppm |
| Nb | 4A_MICP | 141 | 2 | 16 | 12.706 | 2 | 2 | 5 | 8 | 32 | ppm |
| P | 4A_MICP | 112 | 3 | 22 | 4.303 | 3 | 3 | 7 | 6 | 14 | ppm |
| P | XRF | 98 | 2 | 15 | 12.706 | 15 | 15 | 30 | 135 | 192 | ppm |
| Pb | 4A_MICP | 7 | 2 | 16 | 12.706 | 6 | 0.4 | 0.8 | 3 | 5 | ppm |
| Rb | 4A_MICP | 5 | 2 | 16 | 12.706 | 8 | 0.4 | 0.9 | 4 | 5 | ppm |
| Sb | 4A_MICP | 2 | 3 | 23 | 4.303 | 9 | 0.2 | 0.3 | 0.2 | 0.6 | ppm |
| Sc | 4A_MICP | 18 | 3 | 24 | 4.303 | 9 | 2 | 3 | 4 | 7 | ppm |
| Si | XRF | 15.65 | 2 | 16 | 12.706 | 0.5 | 0.082 | 0.16 | 0.69 | 1 | % |
| Sr | 4A_MICP | 130 | 3 | 24 | 4.303 | 11 | 14 | 29 | 35 | 61 | ppm |
| Ta | 4A_MICP | 150 | 2 | 16 | 12.706 | 1 | 2 | 4 | 9 | 27 | ppm |
| Th | 4A_MICP | 0.5 | 3 | 23 | 4.303 | 6 | 0.03 | 0.06 | 0.06 | 0.1 | ppm |
| Ti | 4A_MICP | 3281 | 2 | 16 | 12.706 | 1 | 49 | 98 | 375 | 622 | ppm |
| Ti | FUS | 3561 | 3 | 22 | 4.303 | 5 | 168 | 335 | 354 | 721 | ppm |

| Analyte | Method | ⁷ (μ) | N | n | k | % RSD | ⁸ (u_c) | ⁹ (2s) \pm | ¹⁰ (CI) 95% | ¹¹ (U) \pm | Unit |
|---------|---------|------------------------|---|----|--------|-------|------------------------|-------------------------|------------------------|-------------------------|------|
| U | 4A_MICP | 0.5 | 2 | 16 | 12.706 | 9 | 0.05 | 0.1 | 0.4 | 0.6 | ppm |
| V | 4A_MICP | 857 | 3 | 24 | 4.303 | 4 | 35 | 69 | 85 | 149 | ppm |
| W | 4A_MICP | 2 | 3 | 22 | 4.303 | 4 | 0.08 | 0.2 | 0.2 | 0.3 | ppm |
| Y | 4A_MICP | 5 | 3 | 24 | 4.303 | 8 | 0.4 | 0.8 | 1 | 2 | ppm |
| Zn | 4A_MICP | 305 | 2 | 16 | 12.706 | 0.7 | 2 | 5 | 8 | 29 | ppm |
| Zn | FUS | 330 | 2 | 16 | 12.706 | 2 | 7 | 14 | 32 | 87 | ppm |
| Zr | 4A_MICP | 15 | 3 | 24 | 4.303 | 9 | 1 | 3 | 3 | 5 | ppm |

Table 2. Certified major oxides concentrations, two standard deviations, combined and expanded uncertainty.

| Analyte | Method | ⁷ (μ) | N | n | k | % RSD | ⁸ (u_c) | ⁹ (2s) \pm | ¹⁰ (CI) 95% | ¹¹ (U) \pm | Unit |
|--------------------------------|---------|------------------------|---|----|--------|-------|------------------------|-------------------------|------------------------|-------------------------|------|
| Al ₂ O ₃ | XRF | 14.54 | 3 | 22 | 4.303 | 0.3 | 0.036 | 0.073 | 0.040 | 0.2 | % |
| CaO | XRF | 5.84 | 2 | 16 | 12.706 | 0.2 | 0.014 | 0.029 | 0.11 | 0.2 | % |
| Cr ₂ O ₃ | XRF | 14.91 | 3 | 23 | 4.303 | 0.8 | 0.11 | 0.23 | 0.18 | 0.5 | % |
| Fe ₂ O ₃ | FUS | 16.64 | 4 | 32 | 3.182 | 2 | 0.40 | 0.80 | 0.58 | 1 | % |
| Fe ₂ O ₃ | XRF | 16.60 | 2 | 16 | 12.706 | 1 | 0.16 | 0.32 | 1.0 | 2 | % |
| MgO | XRF | 12.77 | 2 | 16 | 12.706 | 0.2 | 0.028 | 0.056 | 0.032 | 0.4 | % |
| Na ₂ O | 4A_MICP | 0.80 | 5 | 36 | 2.776 | 5 | 0.04 | 0.08 | 0.05 | 0.1 | % |
| Na ₂ O | XRF | 0.79 | 3 | 23 | 4.303 | 3 | 0.02 | 0.05 | 0.05 | 0.1 | % |
| SiO ₂ | XRF | 33.39 | 2 | 16 | 12.706 | 0.2 | 0.070 | 0.14 | 0.31 | 0.9 | % |
| SO ₃ | XRF | 0.13 | 3 | 22 | 4.303 | 6 | 0.007 | 0.01 | 0.02 | 0.03 | % |

1. NiS is fire assay with either ICP/AA/Grav finish
2. Pb Collection is fire assay lead collection with either ICP/AA/Grav finish
3. 2A_MICP is aqua regia digestion with either ICPOES/ICPMS/AAS finish
4. 4A_MICP is a multi-acid digestion with either ICPOES/ICPMS/AAS finish
5. FUS is Fusion digestion with either ICPOES/ICPMS/AAS finish
6. XRF is X-ray Fluorescence
7. The certified value μ , is an unweighted grand mean of the means of N accepted sets of data from different laboratories and n number of test sample replicates. The certified value is traceable to SI units and is reported on a dry basis.
8. The combined uncertainty of the certified value is the within-laboratory reproducibility standard deviation derived from the analysis of variance of results from N number of laboratories and n number of sample replicates. (u_c)
9. Two standard deviations (2s)
10. Confidence interval at 95% level of confidence.
11. Expanded uncertainty (U) at a confidence level of 95% is determined by multiplication of the combined uncertainty (u_c) with a coverage factor (k) found from N-1 degrees of freedom (see Appendix 7 for t-distribution table). Example: $U = 2.36 \times 0.23 = 0.5\%$

2. Statistical Comparison of Means

A comparison of means for replicate data for the same element concentration determined by different analytical methods is done equating the variances between the two data sets; if the variances are found to be equal (F-test, p -value >0.05), then an equal variance t-test is applied. Should the variances be statistically significant, i.e. $p<0.05$, then an unequal variance t-test is performed. For either t-test, if the obtained p -value ≥ 0.05 , the null hypothesis that the means (certified values) are equal is accepted (Table 3). This gives the analyst confidence in the certified values reported by different analytical methods on the same analyte.

Table 3. The results of a two-sample equal or unequal variance t-test (two-tailed) data sets in which different analytical methods /instrumentation were used.

| Method | Certified value | Method | Certified value | F-Test Outcome | p -value (t-test) | t-test Outcome |
|------------------------------------|-----------------|------------------------------------|-----------------|--------------------------------|---------------------|--|
| Au NiS | 0.038 g/t | Au Pb Collection | 0.041 g/t | Equal Variance ($p=0.339$) | 0.496 | Accept H_0 ; certified values are equal |
| Pd NiS | 0.774 g/t | Pd Pb Collection | 0.770 g/t | Equal Variance ($p=0.490$) | 0.928 | Accept H_0 ; certified values are equal |
| Pt NiS | 1.49 g/t | Pt Pb Collection | 1.46 g/t | Equal Variance ($p=0.064$) | 0,75 | Accept H_0 ; certified values are equal |
| Cu 2A_MICP | 0.012 % | Cu 4A_MICP | 0.011 % | Equal Variance ($p=0.248$) | 0.382 | Accept H_0 ; certified values are equal |
| S 4A_MICP | 0.050 % | S XRF | 0.052 % | Equal Variance ($p=0.097$) | 0.387 | Accept H_0 ; certified values are equal |
| Na 4A_MICP | 0.595 % | Na XRF | 0.588 % | Equal Variance ($p=0.167$) | 0.784 | Accept H_0 ; certified values are equal |
| P 4A_MICP | 0.011 % | P XRF | 0,01 % | Unequal Variance ($p=0.022$) | 0.404 | Accept H_0 ; certified values are equal |
| Mn 4A_MICP | 0.147 % | Mn FUS | 0.159 % | Equal Variance ($p=0.122$) | 0.025 | Reject H_0 ; certified values are <i>not</i> equal |
| Zn 4A_MICP | 0.030 % | Zn FUS | 0.033 % | Equal Variance ($p=0.155$) | 0,011 | Reject H_0 ; certified values are <i>not</i> equal |
| Fe ₂ O ₃ FUS | 16.64 % | Fe ₂ O ₃ XRF | 16.60 % | Equal Variance ($p=0.304$) | 0.910 | Accept H_0 ; certified values are equal |
| Na ₂ O 4A_MICP | 0.804 % | Na ₂ O XRF | 0.794 % | Equal Variance ($p=0.206$) | 0.733 | Accept H_0 ; certified values are equal |

Table 4. Comparison of the mean values for which three different analytical methods are used to determine their concentration.

| Method | Certified value | Method | Certified value | Method | Certified value | F-ANOVA p -value | F ANOVA Outcome |
|------------|-----------------|------------|-----------------|--------|-----------------|--------------------|------------------------------|
| Co 2A_MICP | 0.002 % | Co 4A_MICP | 0.015 % | Co FUS | 0.015 % | <0.001 | At least one mean is unequal |
| Ni 2A_MICP | 0.028 % | Ni 4A_MICP | 0.087 % | Ni FUS | 0.087% | <0.001 | At least one mean is unequal |
| Al 4A_MICP | 7.91 % | Al FUS | 7.59 % | Al XRF | 7.70 % | 0.030 | At least one mean is unequal |
| Ca 4A_MICP | 4.05 % | Ca FUS | 4.16 % | Ca XRF | 4.17 % | 0.052 | Means are Equal |
| Fe 4A_MICP | 11.5 % | Fe FUS | 11.3 % | Fe XRF | 11.6 % | 0.810 | Means are Equal |
| Mg 4A_MICP | 7.67 % | Mg FUS | 7.57 % | Mg XRF | 7.68 % | 0.866 | Means are Equal |
| Ti 4A_MICP | 0.328 % | Ti FUS | 0.356 % | Ti XRF | 0.359 % | 0.069 | Means are Equal |
| Co 2A_MICP | 0.002 % | Co 4A_MICP | 0.015 % | Co FUS | 0.015 % | <0.001 | At least one mean is unequal |
| Ni 2A_MICP | 0.028 % | Ni 4A_MICP | 0.087 % | Ni FUS | 0.087% | <0.001 | At least one mean is unequal |

3. Intended Use

AMIS0770 is a matrix matched Certified Reference Material, fit for use as a control sample in routine assay laboratory quality control when inserted within runs of test samples and measured in parallel to test samples. This material can also be used for method development, use as independent calibration verification check standard (*i.e.* if not used as a calibration standard in an instrument calibration), or for validation of accuracy in a method validation exercise (see Appendix 3). The recommend procedure for the use of this CRM as a control standard in laboratory quality control is to develop a Shewhart chart, where a mean value and corresponding 1, 2 and 3 standard deviations are derived from replicate measurements of the CRM (see Appendix 6). This CRM can also be used to assess inter-laboratory or instrument bias and establish within-laboratory precision and within-laboratory reproducibility. The certified concentrations and expanded uncertainty for this material are property values based on an inter-laboratory measurement campaign and reflect consensus results from the laboratories that took part in the exercise

4. Abbreviations and Symbols

Abbreviations and symbols used in this document are shown in Table 5.

Table 5. Abbreviations, symbols and descriptions.

| Abbreviation/Symbol | Description |
|--------------------------|--|
| Alpha (α) | Significance level (denoted by alpha, ' α ') of 0.05 or 5% |
| ANOVA | Analysis of variance by statistical means |
| Bq | The becquerel is the SI derived unit of radioactivity. |
| BIF | Banded iron formation |
| CRM | Certified reference material |
| df | Degrees of freedom, typically, $n-1$, or $N-1$ |
| F_{calc} | Calculated F statistic from ANOVA or Fisher's test |
| F-critical or F_{crit} | F-critical value from F-distribution table |
| GOI | Gain on ignition |
| H_0 | Null hypothesis |
| H_1 | Alternate hypothesis |
| g/t | Grams per tonne |
| k | Coverage factor, e.g. $k=2$ for 95% level of confidence |
| LOC | Level of confidence or confidence level |
| LOD | Limit of detection |
| LOQ | Limit of quantitation |
| LOI | Loss on ignition |
| MS | Mean squares (ANOVA) |
| MSb | Mean squares between(ANOVA) |
| MSw | Mean squares within (ANOVA) |
| N | Number of labs |
| n | Number of replicates |
| μ | Property or certified value of a CRM |
| p | ' p -value' a measure of the strength of evidence against H_0 |
| P | Total number of data points in ANOVA |
| ppm | Parts per million. Equivalent to g/t |
| RSD | Relative standard deviation usually expressed as % at a 68% LOC |
| Replicates | Replication is the repetition of an experimental condition so that the variability associated with an analysis can be estimated (ASTM E1847) |
| s | Standard deviation |
| s_r | Within laboratory repeatability as derived from ANOVA |
| s_s | Between laboratory standard deviation as derived from ANOVA |
| SS | Sum of squares in ANOVA |
| SST | Total variation in ANOVA |
| SSB | Between group (laboratory) variance |
| SSW | Within group (laboratory) variance |

| Abbreviation/Symbol | Description |
|--------------------------|--|
| 2s | Two times standard deviation |
| SI | Standard International system of units |
| t_{calc} | Calculated t statistic from a one-sample, two-tailed t-test |
| t-critical or t_{crit} | t-critical value at given alpha and degrees of freedom |
| Tonne | A metric ton, is a unit of mass equaling 1000 kilograms |
| =TINV(5%, <i>df</i>) | MS Excel function for t-critical value at LOC 95% and <i>df</i> |
| <i>U</i> | Expanded uncertainty at a given k |
| <i>u</i> | Standard uncertainty at k=1 |
| u_c | Combined standard uncertainty at $k=1$ |
| μm | Micron, is an SI derived unit of length equaling 1×10^{-6} of a meter |

5. Uncertified Concentration Values

Appendix 1 gives uncertified concentrations for other elements present in the CRM.

6. Units

All results certified for major oxides are reported as oxides in percentages. All results certified for major elements analyses reported in percentages or ppm. Results for Au and the platinum group elements are reported in g/t or ppm. Specific gravity (SG) is the ratio of the density of a substance to the density of a reference substance, *i.e.* equivalently; it is the ratio of the mass of a substance to the mass of a reference substance for the same given volume. Since specific gravity is a ratio of densities its units are therefore dimensionless.

7. Analytical and Physical Methods

A complete list of analytical and physical methods as generic method codes with a brief description of the methods is available on the AMIS web site www.amis.co.za

8. Origin of Material

The material is Bushveld Complex UG2 ore supplied by Two Rivers Platinum Mine, located on the eastern Bushveld

9. Approximate Mineral and Chemical Composition

UG2 Reef at Two Rivers, consist of namely the 'normal reef' with a thick main chromitite layer; a 'split reef' characterised by an internal pyroxenite/norite lens within the main chromitite layer; and a 'multiple split reef' with numerous pyroxenite/norite lenses occurring within the main chromitite layer.

10. Quantitative Analysis by X-Ray Diffraction

Both natural and synthetic materials have a specific chemistry and atomic arrangement, known as phases. Phases can be identified and quantified using X-ray diffraction (XRD) which produces a plot of the intensity of X-rays scattered at different angles by crystalline phases in a material. Essentially, an X-ray diffraction pattern is the sum of the diffraction patterns produced by each phase. Simply put, an X-ray diffraction pattern is a fingerprint that allows the identification of what is in a target sample material. Knowledge of the mineral phase composition is useful in method development with techniques such as ICP-OES and XRF as potential matrix effects and spectral interferences can be recognized and accounted for. X-ray diffraction is effective in such a way that it allows the identification of different phases of compounds that are identical in chemistry, but have distinctly different atoms, e.g. quartz, cristobalite, and glass are all different phases of SiO₂. Where quantitative XRD results do not correspond to results of other analytical techniques, it should be borne in mind that even though the data are quantitative they are meant to be used for indicative purposes in development of other analytical methods. Mineral names may not reflect the actual compositions of minerals identified, but rather the mineral group.

The sample was micro milled for 10 minutes, with ethanol as the grinding liquid. The resultant sample was lightly pressed into a back-packed sample holder. These results unvalidated and not certified, they are to be used for informational purposes only. The results from a diffractogram of the material gave the following results:

Table 6. Mineral species identified and quantified two sub-samples of AMIS0770 using Rietveld Refinement.

1. Spinel group: Where A=Co, Cu, Fe²⁺, Ge, Mg, Mn²⁺, Ni, Ti, Zn; B=Al, Cr³⁺, Fe²⁺, Fe³⁺, Mg, Mn³⁺, Ti, V³⁺. Chromite was used in these QXRD refinements
2. Mica group: Where X=Ba, Ca, Cs, (H₃O), K, Na, (NH₄); Y=Al, Cr³⁺, Fe²⁺, Fe³⁺, Li, Mg, Mn²⁺, Mn³⁺, V³⁺, Zn; Z=Al, Be, Fe³⁺, Si. Phlogopite (Fe-rich) was used in these QXRD refinements.
3. Amphibole group: Where A=Ca, Na, K, Pb; B=Ca, Fe²⁺, Li, Mg, Mn²⁺, Na; Y=Al, Cr³⁺, Fe²⁺, Fe³⁺, Mg, Mn²⁺, Ti, Z=Al, Be, Si. Magnesiohornblende was used in these QXRD refinements.
4. Plagioclase: Anorthite, sodian was used in these QXRD refinements
5. Pyroxene group: Where A=Ca, Fe²⁺, Li, Mg, Mn²⁺, Na, Zn; B=Al, Cr³⁺, Fe²⁺, Fe³⁺, Mg, Mn²⁺, Sc, Ti, V³⁺; Z=Al, Si. Enstatite and Diopside were used in these QXRD refinements.
6. Due to the lack of other identifying peaks, a number of unidentified peaks/phases were observed. Further work is required to determine the mineral giving rise to this peak. Tr = Trace. i.e. the mineral that gives rise to this peak is assumed to be Presence in trace amounts.

Table 6. Results of XRD analysis.

| Mineral | Mineral Composition | Unit | Composite |
|------------------------------|--|------|-----------|
| ¹ Spinel group | AB ₂ O ₄ | % | 31 |
| ² Mica group | XY ₂₋₃ Z ₄ O ₁₀ (OH,F) ₂ | % | <1 |
| ³ Amphibole group | A ₀₋₁ B ₂ Y ₅ Z ₈ O ₂₂ (OH, F, Cl) ₂ | % | 5 |
| Quartz | SiO ₂ | % | <1 |
| ⁴ Plagioclase | (Na, Ca)Al(Al, Si)Si ₂ O ₈ | % | 24 |
| ⁵ Pyroxene group | ABZ ₂ O ₆ | % | 38 |
| ⁶ Unassigned peak | Unknown | % | Tr |
| Total | | | 100% |

For informational purposes only

11. Health and Safety

The material is a very fine 5YR 6/2 light brownish grey powder. Safety precautions for handling fine particulate matter are recommended, such as the use of safety glasses, breathing protection, gloves and a laboratory coat.

12. Method of Preparation

The particle size distribution for this material was shown to have a nominal top size of 54µm (95% passing 54µm). The procedure of preparation in brief is as follows: the material was crushed, dry-milled and air-classified to <54µm. It was then blended in a bi-conical mixer, systematically divided, and sealed into 1kg Laboratory Packs. Explorer Packs are then subdivided from the Laboratory Packs as required. Final packaged units were then selected on a random basis and submitted for analysis to an independent laboratory accredited with the ISO17025 standard of general requirements for the competence of testing and calibration laboratories. The results obtained from this laboratory are then evaluated statistically by AMIS for homogeneity.

13. Particle Size Determination

The sample has been analysed using a Malvern Mastersizer 2000. Particles are passed through a focused laser beam that scatter light at an angle inversely proportional to their size. The intensity of light is measured and converted to a volume in particle size distribution. The results for this standard are presented in Table 7.

Table 7. Particle Size Determination by laser diffraction.

| Size (µm) | Vol. Under % | Size (µm) | Vol. Under % | Size (µm) | Vol. Under % |
|-----------|--------------|-----------|--------------|-----------|--------------|
| 0.010 | 0.00 | 1.096 | 14.01 | 120.226 | 100.00 |
| 0.011 | 0.00 | 1.259 | 16.16 | 138.038 | 100.00 |
| 0.013 | 0.00 | 1.445 | 18.63 | 158.489 | 100.00 |
| 0.015 | 0.00 | 1.66 | 21.42 | 181.97 | 100.00 |
| 0.017 | 0.00 | 1.905 | 24.52 | 208.93 | 100.00 |
| 0.020 | 0.00 | 2.188 | 27.86 | 239.883 | 100.00 |
| 0.023 | 0.00 | 2.512 | 31.36 | 275.423 | 100.00 |
| 0.026 | 0.00 | 2.884 | 34.92 | 316.228 | 100.00 |
| 0.030 | 0.00 | 3.311 | 38.45 | 363.078 | 100.00 |
| 0.035 | 0.00 | 3.802 | 41.88 | 416.869 | 100.00 |
| 0.040 | 0.00 | 4.365 | 45.35 | 478.63 | 100.00 |
| 0.046 | 0.00 | 5.012 | 48.35 | 549.541 | 100.00 |
| 0.052 | 0.00 | 5.754 | 51.4 | 630.957 | 100.00 |
| 0.060 | 0.00 | 6.607 | 54.38 | 724.436 | 100.00 |
| 0.069 | 0.00 | 7.586 | 57.35 | 831.764 | 100.00 |
| 0.079 | 0.00 | 8.710 | 60.34 | 954.993 | 100.00 |
| 0.091 | 0.00 | 10.000 | 63.4 | 1096.478 | 100.00 |
| 0.105 | 0.00 | 11.482 | 66.57 | 1258.925 | 100.00 |
| 0.120 | 0.00 | 13.183 | 69.84 | 1445.44 | 100.00 |
| 0.138 | 0.00 | 15.136 | 73.21 | 1659.587 | 100.00 |
| 0.158 | 0.00 | 17.378 | 76.64 | 1905.461 | 100.00 |
| 0.182 | 0.00 | 19.953 | 80.07 | 2187.762 | 100.00 |
| 0.209 | 0.07 | 22.909 | 83.43 | 2511.886 | 100.00 |
| 0.240 | 0.44 | 26.303 | 86.64 | 2884.032 | 100.00 |
| 0.275 | 1.00 | 30.20 | 89.60 | 3311.311 | 100.00 |
| 0.316 | 1.76 | 34.674 | 92.25 | 3801.894 | 100.00 |
| 0.363 | 2.70 | 39.811 | 94.50 | 4365.158 | 100.00 |
| 0.417 | 3.78 | 45.709 | 96.34 | 5011.872 | 100.00 |
| 0.479 | 4.96 | 52.481 | 97.74 | 5754.399 | 100.00 |
| 0.550 | 6.22 | 60.256 | 98.75 | 6606.934 | 100.00 |
| 0.631 | 7.55 | 69.183 | 99.40 | 7585.776 | 100.00 |
| 0.724 | 8.95 | 79.433 | 99.79 | 8709.636 | 100.00 |
| 0.832 | 10.47 | 91.201 | 99.96 | 10000 | 100.00 |
| 0.955 | 12.13 | 104.713 | 100.00 | | |

For informational purposes only

14. Handling

The material is packaged in Laboratory Packs and Explorer Packs that must be shaken or otherwise agitated before use. The analyte concentrations are quoted on a dry basis; therefore, the user needs to determine the moisture content to convert any obtained assay values to an air-dry basis (see Appendix 7 for an example calculation).

15. Storage information

The material should be stored in a cool dry place, in such a way that it does not compromise the integrity of the CRM. The material should be stored in conditions which will ensure it does not absorb moisture.

16. Methods of Analysis Requested

The following methods of analysis were requested:

- a) Pt, Pd, Au, Rh, Ru, Ir – Nickel Sulphide with ICP-OES or ICP-MS or AAS finish
- b) Rh, Pt, Pd and Au-Pb collection with Ag as a co-collector, ICP-OES or ICP-MS or AAS finish
- c) Multi element scan to include especially Cu, Co, Ni and Cr-by 4-acid digestion by ICP-OES or ICP-MS or AAS finish
- d) Co, Cu and Ni: Aqua Regia digestion (2-acid), including by ICP-OES or ICP-MS
- e) Multi element scan to include especially Cu, Co, Ni and Cr-by peroxide fusion n by ICP-OES or ICP-MS or AAS finish
- f) Majors oxides (Al_2O_3 , CaO , Cr_2O_3 , Fe_2O_3 , K_2O , MgO , MnO , Na_2O , SiO_2 , TiO_2 , P_2O_5 and LOI) by XRF finish
- g) SG – gas pycnometer
- h) S and C Combustion/LECO
- i) Moisture

17. Information Requested of Participating Laboratories

The following information was requested of the participating laboratories for the development of this CRM:

- a) All aliquots used for all determinations to be stated.
- b) All results for major elements to be reported as oxides in percentage (%).
- c) All results for multi-element scans and fusion to be reported in parts per million (ppm).
- d) All results for 6E (Au, Pt, Pd, Rh, Ru and Ir) to be reported in grams per ton (g/t).
- e) All QC data to include replicates, blanks and any certified reference material used.
- f) All Round robin samples must be treated the same as routine test samples.
- g) All results must be reported to maximum decimal places i.e. dependent on laboratories capabilities.
- h) Moisture correction factor to be used to calculate and determine the moisture content. The moisture correction factor is to be stipulated when reporting results.
- i) Results to be recorded onto the PT Reporting template provided (Include uncertainty of measurements) together with a copy of the PDF.
- j) Ensure correct PPE is used i.e. gloves, dust masks and protective clothing when handling and analyzing the samples
- k) Storage and analysis shall be conducted under controlled environmental conditions.
- l) Ensure all safety requirements as per your laboratory requirements are adhered to.

18. Certification of Mean and Estimation of Measurement Uncertainty

The samples used in this certification process have been selected in such a way as to represent the entire batch of material and were taken from the final packaged units; therefore, all sources of uncertainty are included in the combined standard uncertainty determination. Initially the data submitted by all the laboratories are subjected to a z-score test, equation [1] to exclude outliers and the remaining data sets examined for their normality in distribution. This is followed by the exclusion of further outliers as defined by the IUPAC Harmonized Protocol of 1995 in which both Cochran and a Grubbs test are applied until all outliers are identified, equations [2] and [3]. A grand mean and standard deviation are re-calculated using all remaining data (Thompson, 2008; Carr, 2011) (see Appendix 2)

19. Two Standard Deviations

AMIS reports two-standard deviations (2s) with all certified values. Two -standard deviations are calculated using the expression:

$$\text{Two standard deviations} = 2 (u_c)$$

Where u_c is the standard combined uncertainty (see Appendix, equation [14]).

20. Confidence Interval

AMIS reports a confidence interval (CI) with all certified values. Confidence interval as used by AMIS is:

$$\text{Confidence Interval (CI)} = \frac{(t_{critical})s}{\sqrt{N}}$$

Where, N is the number of laboratories (accepted laboratory data), $t_{critical}$ is a two-tailed value for $N - 1$ degrees of freedom (df) and s , is the standard deviation of the accepted laboratory means. A two-tailed critical value is found for $N - 1$ degrees of freedom from either a t -distribution table (Appendix 9) or MS Excel as =TINV (5%, df).

21. Expanded Uncertainty

ANOVA gives an estimate of the repeatability and the reproducibility of the data accepted for certification of the candidate reference material (see equations, [15] and [16], in the Appendix). Therefore, random variables (e.g., subsampling, instrument effects, interferences, operators and measurement conditions) that occur during the analysis of the candidate reference material by the various laboratories is considered. This approach does not necessarily quantify each individual source of uncertainty; however, the combined effect of random uncertainties is assessed (Ramsey & Ellison, 2007). A combined standard uncertainty is calculated from equation [14], which when multiplied by the t -critical value for $N - 1$ laboratories, gives an *expanded uncertainty* at a 95% level of confidence. The expanded uncertainty is a measure of the doubt around the certified value at a level of 95% confidence. The expanded uncertainty is used in the validation of accuracy (see equation [18]).

22. Confidence Interval and Expanded Uncertainty

A combined standard uncertainty will be greater than a combined CI. This is because ANOVA considers the within-lab repeatability (that is repeatability within each lab group) as well as the repeatability between each lab data set. This attends to random variables that contribute to the measurement of uncertainty, during the analysis of the test sample at the participating laboratories. The within-lab repeatability and the between lab repeatability is combined as the square root of the sum of squares of these two values giving a combined standard uncertainty, at a 68% confidence level. Multiplying the

combined standard uncertainty by the t-critical value for $N-1$, gives the expanded uncertainty at 95% level of confidence. It is recommended that the procedure described in Appendix 6, "Using the CRM in Quality Control" be used, in setting the limits of the CRM. Table 8 below shows mean gold values obtained by fire assay lead collection, for nine different laboratories, the confidence interval, two-standard deviations and expanded uncertainty.

Table 8. Example of replicate assay data in which the *CI*, *2s* and *U* are shown.

| Lab No. | Mean Au (g/t) |
|---------|---------------|
| 1 | 0.268 |
| 2 | 0.273 |
| 3 | 0.270 |
| 4 | 0.288 |
| 5 | 0.274 |
| 6 | 0.256 |
| 7 | 0.263 |
| 8 | 0.258 |
| 9 | 0.288 |

| | |
|-----------|--------|
| CI | 0.0088 |
| 2s | 0.031 |
| U | 0.04 |

23. Participating Laboratories

The laboratories that provided timeous results are:

1. AGAT Calgary
2. ALS Geochemistry SA
3. Anglo American Technical Solutions
4. Antech Zimbabwe
5. Bureau Veritas Minerals Ultra Trace Pty Ltd
6. CM Solutions
7. Intertek Perth
8. Jordan Uranium Mining Company
9. Mintek (South Africa)
10. Professional Platinum Services
11. SGS Brits
12. SGS Rustenburg
13. Tharisa Mine

24. Accepted Assay Data

Data from the 13 laboratories used for certification are set out in Table 9.

Table 9. Data used to calculate the certified values after removal of outliers.

| NiS Au g/t | Pb Collection Au g/t | NiS Ir g/t | NiS Pd g/t | Pb Collection Pd g/t | NiS Pt g/t | Pb Collection Pt g/t | NiS Rh g/t | NiS Ru g/t |
|------------------|-------------------------------|------------------|------------------|-------------------------------|------------------|-------------------------------|------------------|------------------|
| 0.035 | 0.043 | 0.10 | 0.78 | 0.76 | 1.58 | 1.51 | 0.26 | 0.46 |
| 0.033 | 0.039 | 0.10 | 0.77 | 0.75 | 1.56 | 1.49 | 0.26 | 0.46 |
| 0.032 | 0.038 | 0.10 | 0.76 | 0.75 | 1.52 | 1.48 | 0.26 | 0.40 |
| 0.032 | 0.037 | 0.10 | 0.76 | 0.74 | 1.52 | 1.43 | 0.25 | 0.40 |
| 0.031 | 0.037 | 0.09 | 0.74 | 0.73 | 1.45 | 1.42 | 0.25 | 0.39 |
| 0.031 | 0.037 | 0.09 | 0.74 | 0.72 | 1.42 | 1.37 | 0.25 | 0.51 |
| 0.029 | 0.035 | 0.09 | 0.74 | 0.71 | 1.37 | 1.36 | 0.25 | 0.50 |
| 0.028 | 0.030 | 0.09 | 0.73 | 0.69 | 1.49 | 1.35 | 0.25 | 0.49 |
| 0.047 | 0.050 | 0.16 | 0.69 | 0.78 | 1.48 | 1.59 | 0.26 | 0.48 |
| 0.046 | 0.040 | 0.15 | 0.68 | 0.76 | 1.45 | 1.50 | 0.26 | 0.48 |
| 0.043 | 0.040 | 0.15 | 0.67 | 0.76 | 1.43 | 1.48 | 0.24 | 0.48 |
| 0.042 | 0.040 | 0.15 | 0.66 | 0.75 | 1.41 | 1.47 | 0.24 | 0.47 |
| 0.042 | 0.040 | 0.14 | 0.65 | 0.75 | 1.40 | 1.45 | 0.24 | 0.46 |
| 0.041 | 0.040 | 0.12 | 0.65 | 0.74 | 1.38 | 1.39 | 0.23 | 0.47 |
| 0.039 | 0.040 | 0.11 | 0.65 | 0.71 | 1.36 | 1.31 | 0.23 | 0.47 |
| 0.041 | 0.030 | 0.11 | 0.83 | 0.66 | 1.40 | 1.56 | 0.27 | 0.45 |
| 0.040 | 0.042 | 0.11 | 0.82 | 0.90 | 1.35 | 1.56 | 0.27 | 0.45 |
| 0.039 | 0.042 | 0.11 | 0.82 | 0.85 | 1.32 | 1.55 | 0.26 | 0.45 |
| 0.038 | 0.039 | 0.11 | 0.82 | 0.85 | 1.31 | 1.52 | 0.26 | 0.45 |
| 0.036 | 0.038 | 0.11 | 0.81 | 0.84 | 1.29 | 1.52 | 0.26 | 0.43 |
| 0.035 | 0.035 | 0.11 | 0.78 | 0.82 | 1.28 | 1.50 | 0.26 | 0.51 |
| 0.035 | 0.033 | 0.11 | 0.78 | 0.81 | 1.27 | 1.41 | 0.25 | 0.51 |
| 0.034 | 0.031 | 0.11 | 0.78 | 0.81 | 1.56 | 1.41 | 0.25 | 0.51 |
| 0.040 | 0.029 | 0.11 | 0.88 | | 1.55 | | 0.30 | 0.50 |
| 0.040 | 0.050 | 0.11 | 0.88 | | 1.54 | | 0.30 | 0.50 |
| 0.040 | 0.050 | 0.11 | 0.84 | | 1.50 | | 0.29 | 0.50 |
| 0.040 | 0.050 | 0.11 | 0.84 | | 1.49 | | 0.29 | 0.49 |
| 0.040 | 0.050 | 0.11 | 0.82 | | 1.48 | | 0.29 | 0.48 |
| 0.040 | 0.050 | 0.11 | 0.82 | | 1.46 | | 0.29 | 0.57 |
| 0.040 | 0.050 | | 0.81 | | 1.46 | | 0.28 | 0.56 |
| 0.040 | | | 0.82 | | 1.79 | | 0.27 | 0.54 |
| | | | 0.82 | | 1.79 | | | 0.54 |
| | | | 0.81 | | 1.78 | | | 0.54 |
| | | | 0.81 | | 1.73 | | | 0.53 |
| | | | 0.81 | | 1.71 | | | 0.52 |
| | | | 0.80 | | 1.69 | | | 0.51 |
| | | | 0.80 | | 1.68 | | | |
| | | | 0.80 | | 1.66 | | | |

Assay Data (Continued)

| 2A_MICP | 4A_MICP | FUS | 2A_MICP | 4A_MICP | 2A_MICP | 4A_MICP | FUS | 4A_MICP | XRF |
|---------|---------|--------|---------|---------|---------|---------|--------|---------|-------|
| Co | Co | Co | Cu | Cu | Ni | Ni | Ni | S | S |
| ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | % | % |
| 17.32 | 149.00 | 167.76 | 125.02 | 108.62 | 300.00 | 861.66 | 939.94 | 0.050 | 0.055 |
| 17.04 | 147.50 | 160.18 | 123.71 | 108.48 | 299.00 | 843.75 | 907.54 | 0.050 | 0.055 |
| 16.78 | 147.50 | 158.54 | 123.56 | 106.87 | 299.00 | 837.72 | 903.28 | 0.050 | 0.054 |
| 16.45 | 146.50 | 157.49 | 123.12 | 106.53 | 298.00 | 826.68 | 901.77 | 0.050 | 0.054 |
| 16.30 | 146.50 | 157.30 | 122.49 | 106.51 | 295.00 | 826.65 | 899.25 | 0.050 | 0.054 |
| 16.26 | 146.00 | 156.86 | 122.06 | 105.44 | 293.00 | 825.94 | 896.08 | 0.050 | 0.053 |
| 15.99 | 144.50 | 156.54 | 121.95 | 104.87 | 292.00 | 814.49 | 883.02 | 0.050 | 0.053 |
| 15.97 | 144.50 | 152.01 | 119.26 | 104.07 | 292.00 | 899.00 | 876.80 | 0.050 | 0.052 |
| 16.00 | 156.40 | 170.00 | 113.00 | 117.00 | 263.00 | 891.00 | 910.00 | 0.050 | 0.052 |
| 15.50 | 154.00 | 170.00 | 112.00 | 116.00 | 262.00 | 888.00 | 900.00 | 0.050 | 0.048 |
| 14.50 | 153.20 | 170.00 | 112.00 | 115.50 | 262.00 | 888.00 | 900.00 | 0.050 | 0.048 |
| 14.40 | 153.00 | 170.00 | 112.00 | 114.00 | 261.00 | 882.00 | 900.00 | 0.050 | 0.048 |
| 13.90 | 152.80 | 160.00 | 112.00 | 113.00 | 261.00 | 879.00 | 890.00 | 0.050 | 0.048 |
| 13.90 | 152.70 | 160.00 | 112.00 | 112.00 | 260.00 | 876.00 | 890.00 | 0.049 | 0.048 |
| 13.70 | 152.40 | 150.00 | 111.00 | 109.50 | 260.00 | 874.00 | 890.00 | 0.049 | 0.055 |
| 13.60 | 151.40 | 150.00 | 111.00 | 116.30 | 259.00 | 895.60 | 870.00 | 0.049 | 0.055 |
| 16.00 | | 155.00 | | 116.00 | | 884.00 | 918.00 | | 0.053 |
| 16.00 | | 153.00 | | 115.90 | | 883.30 | 916.00 | | 0.053 |
| 16.00 | | 153.00 | | 115.90 | | 879.60 | 913.00 | | 0.052 |
| 16.00 | | 153.00 | | 115.30 | | 878.10 | 911.00 | | 0.052 |
| 15.00 | | 152.00 | | 114.40 | | 867.70 | 908.00 | | 0.052 |
| 15.00 | | 152.00 | | 113.80 | | 866.50 | 907.00 | | 0.052 |
| 15.00 | | 150.00 | | 113.70 | | 865.70 | 905.00 | | |
| 15.00 | | 149.00 | | | | | 902.00 | | |
| | | 149.80 | | | | | 873.70 | | |
| | | 149.60 | | | | | 865.00 | | |
| | | 149.40 | | | | | 862.40 | | |
| | | 145.10 | | | | | 860.00 | | |
| | | 145.00 | | | | | 858.70 | | |
| | | 145.00 | | | | | 856.60 | | |
| | | 144.50 | | | | | 855.40 | | |
| | | 139.80 | | | | | 852.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 160.00 | | | | | 800.00 | | |
| | | 150.00 | | | | | 800.00 | | |
| | | 153.00 | | | | | | | |
| | | 151.00 | | | | | | | |
| | | 145.00 | | | | | | | |
| | | 142.00 | | | | | | | |
| | | 141.00 | | | | | | | |
| | | 141.00 | | | | | | | |

Assay Data (Continued)

| XRF | XRF | XRF | FUS | XRF | XRF | 4A_MICP | XRF | XRF | XRF |
|--------------------------------|------|--------------------------------|--------------------------------|--------------------------------|-------|-------------------|-------------------|------------------|-----------------|
| Al ₂ O ₃ | CaO | Cr ₂ O ₃ | Fe ₂ O ₃ | Fe ₂ O ₃ | MgO | Na ₂ O | Na ₂ O | SiO ₂ | SO ₃ |
| % | % | % | % | % | % | % | % | % | % |
| 14.56 | 5.85 | 15.15 | 16.73 | 16.74 | 12.80 | 0.88 | 0.80 | 33.42 | 0.14 |
| 14.52 | 5.85 | 15.13 | 16.58 | 16.74 | 12.79 | 0.86 | 0.79 | 33.40 | 0.14 |
| 14.52 | 5.85 | 15.00 | 16.44 | 16.74 | 12.78 | 0.86 | 0.79 | 33.39 | 0.14 |
| 14.52 | 5.85 | 14.95 | 16.30 | 16.72 | 12.77 | 0.85 | 0.78 | 33.39 | 0.14 |
| 14.52 | 5.84 | 14.94 | 16.30 | 16.71 | 12.76 | 0.85 | 0.78 | 33.37 | 0.13 |
| 14.51 | 5.84 | 14.92 | 16.23 | 16.70 | 12.76 | 0.83 | 0.78 | 33.33 | 0.13 |
| 14.51 | 5.84 | 14.90 | 16.08 | 16.70 | 12.75 | 0.83 | 0.78 | 33.33 | 0.13 |
| 14.51 | 5.84 | 15.01 | 15.73 | 16.67 | 12.73 | 0.82 | 0.78 | 33.32 | 0.13 |
| 14.61 | 5.84 | 15.00 | 17.03 | 16.58 | 12.82 | 0.77 | 0.84 | 33.56 | 0.13 |
| 14.58 | 5.84 | 14.95 | 17.03 | 16.53 | 12.79 | 0.75 | 0.83 | 33.47 | 0.12 |
| 14.57 | 5.83 | 14.89 | 16.94 | 16.49 | 12.78 | 0.73 | 0.83 | 33.44 | 0.12 |
| 14.51 | 5.83 | 14.88 | 16.93 | 16.48 | 12.78 | 0.73 | 0.82 | 33.44 | 0.12 |
| 14.51 | 5.83 | 14.76 | 16.90 | 16.47 | 12.78 | 0.82 | 0.81 | 33.42 | 0.12 |
| 14.50 | 5.82 | 14.74 | 16.83 | 16.47 | 12.76 | 0.81 | 0.80 | 33.37 | 0.12 |
| 14.61 | 5.82 | 14.66 | 16.80 | 16.47 | 12.74 | 0.81 | 0.80 | 33.36 | 0.14 |
| 14.58 | 5.81 | 14.94 | 16.77 | 16.46 | 12.72 | 0.81 | 0.80 | 33.28 | 0.14 |
| 14.57 | | 14.91 | 17.16 | | | 0.81 | 0.79 | | 0.13 |
| 14.56 | | 14.90 | 17.16 | | | 0.81 | 0.79 | | 0.13 |
| 14.56 | | 14.90 | 17.01 | | | 0.80 | 0.78 | | 0.13 |
| 14.54 | | 14.89 | 17.01 | | | 0.80 | 0.78 | | 0.13 |
| 14.51 | | 14.86 | 17.01 | | | 0.79 | 0.78 | | 0.13 |
| 14.49 | | 14.86 | 16.87 | | | 0.79 | 0.76 | | 0.13 |
| | | 14.85 | 16.87 | | | 0.78 | 0.76 | | |
| | | | 16.87 | | | 0.78 | | | |
| | | | 16.54 | | | 0.78 | | | |
| | | | 16.47 | | | 0.78 | | | |
| | | | 16.39 | | | 0.77 | | | |
| | | | 16.36 | | | 0.77 | | | |
| | | | 16.32 | | | 0.85 | | | |
| | | | 16.29 | | | 0.85 | | | |
| | | | 16.29 | | | 0.84 | | | |
| | | | 16.16 | | | 0.84 | | | |
| | | | | | | 0.84 | | | |
| | | | | | | 0.84 | | | |
| | | | | | | 0.83 | | | |
| | | | | | | 0.82 | | | |

Assay Data (Continued)

| 4A_MICP | FUS | XRF | 4A_MICP | 4A_MICP | 4A_MICP | 4A_MICP | 4A_MICP | FUS | XRF |
|----------|----------|----------|---------|---------|---------|---------|----------|----------|----------|
| Al | Al | Al | As | Ba | Be | Bi | Ca | Ca | Ca |
| ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm |
| 81995.65 | 76700.00 | 77080.00 | 1.99 | 55.42 | 0.50 | 0.14 | 40600.00 | 42918.55 | 41810.00 |
| 81294.40 | 76500.00 | 76870.00 | 1.96 | 54.66 | 0.48 | 0.13 | 40600.00 | 42434.84 | 41800.00 |
| 80335.26 | 75200.00 | 76850.00 | 1.92 | 54.43 | 0.47 | 0.13 | 40600.00 | 40592.22 | 41790.00 |
| 80067.25 | 74900.00 | 76850.00 | 1.86 | 53.78 | 0.47 | 0.13 | 40400.00 | 40359.25 | 41770.00 |
| 79714.38 | 74700.00 | 76850.00 | 1.83 | 53.69 | 0.46 | 0.13 | 40300.00 | 39560.79 | 41760.00 |
| 79699.62 | 73800.00 | 76810.00 | 1.82 | 53.62 | 0.44 | 0.13 | 40200.00 | 39393.98 | 41760.00 |
| 79400.00 | 73800.00 | 76790.00 | 1.77 | 53.13 | 0.44 | 0.12 | 40100.00 | 43500.00 | 41750.00 |
| 78800.00 | 72100.00 | 76770.00 | 1.76 | 52.76 | 0.43 | 0.12 | 40100.00 | 43100.00 | 41750.00 |
| 78600.00 | 77800.00 | 77433.00 | 2.40 | 60.00 | 0.51 | 0.12 | 41361.00 | 42600.00 | 41756.00 |
| 78400.00 | 77700.00 | 77274.00 | 2.20 | 60.00 | 0.50 | 0.12 | 41300.00 | 42300.00 | 41756.00 |
| 78200.00 | 77500.00 | 77221.00 | 2.10 | 60.00 | 0.49 | 0.11 | 40967.00 | 42000.00 | 41684.50 |
| 77900.00 | 77500.00 | 76903.00 | 2.10 | 60.00 | 0.48 | 0.11 | 40921.00 | 41500.00 | 41684.50 |
| 77800.00 | 77500.00 | 76903.00 | 2.00 | 60.00 | 0.46 | 0.11 | 40880.00 | 41300.00 | 41684.50 |
| 77100.00 | 77500.00 | 76850.00 | 2.00 | 60.00 | 0.45 | 0.11 | 40763.00 | 40500.00 | 41613.00 |
| 79478.00 | 76900.00 | 77433.00 | 2.00 | 60.00 | 0.45 | 0.09 | 40741.00 | 42000.00 | 41613.00 |
| 79195.00 | 76600.00 | 77274.00 | 2.00 | 60.00 | 0.44 | 0.11 | 40388.00 | 42000.00 | 41541.50 |
| 78513.00 | 76439.00 | 77221.00 | | 57.80 | 0.51 | 0.11 | 41391.18 | 42000.00 | |
| 78457.00 | 75985.00 | 77168.00 | | 57.10 | 0.51 | 0.10 | 41268.10 | 42000.00 | |
| 78279.00 | 75726.00 | 77168.00 | | 57.10 | 0.48 | 0.10 | 40497.77 | 42000.00 | |
| 78053.00 | 75704.00 | 77062.00 | | 57.00 | 0.48 | 0.10 | 40253.58 | 42000.00 | |
| 77824.00 | 75594.00 | 76903.00 | | 56.90 | 0.47 | 0.10 | 40197.11 | 41000.00 | |
| 77142.00 | 75291.00 | 76797.00 | | 56.50 | 0.47 | 0.10 | 39948.36 | 41000.00 | |
| | 75198.00 | | | 55.80 | 0.46 | 0.09 | 39322.26 | | |
| | 75080.00 | | | 55.70 | 0.46 | | 39306.51 | | |

| 4A_MICP | XRF | 4A_MICP | 4A_MICP | FUS | XRF | 4A_MICP | 4A_MICP | 4A_MICP | 4A_MICP |
|---------|-----------|---------|-----------|-----------|-----------|---------|---------|---------|---------|
| Ce | Cr | Cs | Fe | Fe | Fe | Ga | Hf | K | Li |
| ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm |
| 6.30 | 103620.00 | 0.29 | 114000.00 | 117000.00 | 117100.00 | 26.27 | 0.64 | 1613.53 | 7.13 |
| 6.19 | 103490.00 | 0.29 | 113000.00 | 116000.00 | 117050.00 | 25.98 | 0.64 | 1608.36 | 7.03 |
| 6.16 | 102640.00 | 0.28 | 113000.00 | 115000.00 | 117050.00 | 25.62 | 0.63 | 1578.55 | 6.95 |
| 6.16 | 102310.00 | 0.28 | 112500.00 | 114000.00 | 116960.00 | 25.60 | 0.63 | 1571.31 | 6.91 |
| 6.08 | 102190.00 | 0.28 | 112500.00 | 114000.00 | 116850.00 | 25.50 | 0.62 | 1550.51 | 6.72 |
| 6.06 | 102090.00 | 0.28 | 112000.00 | 113500.00 | 116810.00 | 25.45 | 0.62 | 1543.10 | 6.72 |
| 6.06 | 101970.00 | 0.27 | 112000.00 | 112500.00 | 116810.00 | 25.07 | 0.61 | 1496.22 | 6.71 |
| 6.02 | 102698.42 | 0.26 | 111000.00 | 110000.00 | 116610.00 | 25.03 | 0.60 | 1300.00 | 6.65 |
| 6.24 | 102630.00 | 0.35 | 119600.00 | 119100.00 | 115960.52 | 29.40 | 0.60 | 1300.00 | 7.70 |
| 6.12 | 102287.90 | 0.34 | 119400.00 | 119100.00 | 115610.82 | 29.30 | 0.60 | 1300.00 | 7.30 |
| 6.10 | 101877.38 | 0.34 | 118400.00 | 118500.00 | 115331.06 | 29.10 | 0.60 | 1300.00 | 7.30 |
| 6.08 | 101808.96 | 0.34 | 118300.00 | 118400.00 | 115261.12 | 28.80 | 0.60 | 1300.00 | 7.20 |
| 6.07 | 100987.92 | 0.33 | 118000.00 | 118200.00 | 115191.18 | 28.60 | 0.60 | 1300.00 | 7.20 |
| 6.04 | 100851.08 | 0.32 | 117900.00 | 117700.00 | 115191.18 | 28.00 | 0.60 | 1300.00 | 7.10 |
| 6.03 | 100303.72 | 0.32 | 117700.00 | 117500.00 | 115191.18 | 28.00 | 0.67 | 1300.00 | 7.10 |
| 6.03 | 102205.80 | 0.32 | 116600.00 | 117300.00 | 115121.24 | 27.80 | 0.65 | 1353.00 | 7.00 |
| | 101980.01 | 0.31 | 115405.61 | 103280.40 | | 29.15 | 0.63 | 1351.00 | 8.40 |
| | 101938.96 | 0.31 | 115260.59 | 103280.40 | | 29.06 | 0.63 | 1345.00 | 8.40 |
| | 101911.59 | 0.30 | 115171.62 | 100454.82 | | 29.03 | 0.63 | 1343.00 | 8.30 |
| | 101850.01 | 0.30 | 115098.01 | 97398.44 | | 28.68 | 0.62 | 1337.00 | 8.30 |
| | 101637.91 | 0.30 | 114888.00 | 120000.00 | | 28.67 | 0.62 | 1333.00 | 8.30 |
| | 101637.91 | 0.29 | 114774.84 | 120000.00 | | 28.58 | 0.61 | 1333.00 | 8.30 |
| | 101576.33 | 0.29 | 114623.05 | 119000.00 | | 28.36 | | 1327.00 | 8.20 |
| | | 0.29 | 114545.88 | 119000.00 | | 28.34 | | 1454.55 | 8.00 |
| | | | | 119000.00 | | | | 1453.69 | |
| | | | | 118000.00 | | | | 1441.72 | |
| | | | | 118000.00 | | | | 1417.16 | |
| | | | | 118000.00 | | | | 1402.08 | |
| | | | | 115722.00 | | | | 1387.89 | |
| | | | | 115222.00 | | | | 1355.97 | |
| | | | | 114632.00 | | | | 1321.20 | |
| | | | | 114452.00 | | | | | |
| | | | | 114162.00 | | | | | |
| | | | | 113971.00 | | | | | |
| | | | | 113967.00 | | | | | |
| | | | | 113008.00 | | | | | |

Assay Data (Continued)

| 4A_MICP | FUS | XRF | FUS | 4A_MICP | XRF | 4A_MICP | 4A_MICP | XRF |
|----------|----------|----------|---------|---------|---------|---------|---------|--------|
| Mg | Mg | Mg | Mn | Na | Na | Nb | P | P |
| ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm |
| 83360.49 | 78300.00 | 77160.00 | 1616.94 | 6557.03 | 5910.00 | 145.50 | 118.09 | 87.20 |
| 82284.75 | 78200.00 | 77100.00 | 1584.29 | 6387.08 | 5890.00 | 144.50 | 117.41 | 87.20 |
| 81900.67 | 76800.00 | 77050.00 | 1564.46 | 6384.72 | 5850.00 | 142.50 | 117.16 | 87.20 |
| 80838.87 | 76700.00 | 77020.00 | 1562.32 | 6329.13 | 5810.00 | 141.00 | 114.81 | 87.20 |
| 80833.81 | 76600.00 | 76970.00 | 1510.74 | 6283.60 | 5810.00 | 140.50 | 113.92 | 87.20 |
| 80298.43 | 75800.00 | 76950.00 | 1508.86 | 6176.97 | 5790.00 | 138.50 | 112.30 | 87.20 |
| 77300.00 | 75700.00 | 76860.00 | 1500.04 | 6123.15 | 5760.00 | 137.00 | 111.36 | 87.20 |
| 77000.00 | 73400.00 | 76780.00 | 1497.99 | 6078.31 | 5750.00 | 136.50 | 120.00 | 113.36 |
| 76700.00 | 77800.00 | 76920.00 | 1700.00 | 5700.00 | 6216.00 | 143.65 | 110.00 | 109.00 |
| 76700.00 | 77700.00 | 76740.00 | 1700.00 | 5600.00 | 6142.00 | 143.19 | 110.00 | 109.00 |
| 76500.00 | 77700.00 | 76680.00 | 1700.00 | 5400.00 | 6142.00 | 142.87 | 110.00 | 109.00 |
| 76400.00 | 77600.00 | 76680.00 | 1700.00 | 5400.00 | 6068.00 | 142.38 | 110.00 | 109.00 |
| 76200.00 | 77400.00 | 76680.00 | 1700.00 | 5300.00 | 5994.00 | 141.62 | 110.00 | 109.00 |
| 76100.00 | 77300.00 | 76560.00 | 1600.00 | 5200.00 | 5920.00 | 141.30 | 110.00 | 104.64 |
| 74984.00 | 76900.00 | 76440.00 | 1600.00 | 6100.00 | 5920.00 | 140.74 | 112.00 | 104.64 |
| 74806.00 | 76700.00 | 76320.00 | 1600.00 | 6000.00 | 5920.00 | 140.33 | 112.00 | |
| 74692.00 | 73568.00 | | 1623.00 | 6000.00 | 5846.00 | | 111.00 | |
| 74467.00 | 73519.00 | | 1609.00 | 6000.00 | 5846.00 | | 111.00 | |
| 74432.00 | 73188.00 | | 1607.00 | 6000.00 | 5772.00 | | 110.00 | |
| 74226.00 | 72997.00 | | 1606.00 | 6000.00 | 5772.00 | | 110.00 | |
| 74198.00 | 72949.00 | | 1605.00 | 5900.00 | 5772.00 | | 110.00 | |
| 73384.00 | 72920.00 | | 1597.00 | 5900.00 | 5624.00 | | 110.00 | |
| 75111.50 | 72830.00 | | 1586.00 | 5872.00 | 5624.00 | | | |
| 74961.89 | | | 1580.00 | 5859.00 | | | | |
| 74701.65 | | | 1556.00 | 5796.00 | | | | |
| 74366.77 | | | 1539.00 | 5793.00 | | | | |
| 73937.55 | | | 1537.00 | 5781.00 | | | | |
| 73591.88 | | | 1536.00 | 5754.00 | | | | |
| 73155.04 | | | 1534.00 | 5746.00 | | | | |
| 72650.70 | | | 1526.00 | 5721.00 | | | | |
| | | | 1526.00 | 6301.83 | | | | |
| | | | 1518.00 | 6300.24 | | | | |
| | | | | 6249.37 | | | | |
| | | | | 6241.56 | | | | |
| | | | | 6235.32 | | | | |
| | | | | 6223.18 | | | | |
| | | | | 6174.94 | | | | |
| | | | | 6104.10 | | | | |

Assay Data (Continued)

| 4A_MICP | 4A_MICP | 4A_MICP | 4A_MICP | XRF | 4A_MICP | 4A_MICP | 4A_MICP | 4A_MICP | FUS |
|---------|---------|---------|---------|-----------|---------|---------|---------|---------|---------|
| Pb | Rb | Sb | Sc | Si | Sr | Ta | Th | Ti | Ti |
| ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm | ppm |
| 7.97 | 4.89 | 1.77 | 17.84 | 156230.00 | 116.13 | 151.50 | 0.54 | 3290.00 | 3679.87 |
| 7.63 | 4.87 | 1.73 | 17.18 | 156120.00 | 115.98 | 151.50 | 0.53 | 3280.00 | 3547.80 |
| 7.58 | 4.86 | 1.71 | 17.00 | 156090.00 | 115.94 | 151.00 | 0.52 | 3260.00 | 3426.08 |
| 7.55 | 4.82 | 1.57 | 16.64 | 156060.00 | 115.76 | 150.50 | 0.52 | 3260.00 | 3406.51 |
| 7.49 | 4.80 | 1.57 | 16.48 | 155970.00 | 114.79 | 148.50 | 0.51 | 3240.00 | 3384.43 |
| 7.15 | 4.71 | 1.50 | 16.34 | 155810.00 | 114.09 | 148.50 | 0.51 | 3240.00 | 3371.51 |
| 7.15 | 4.70 | 1.47 | 16.18 | 155790.00 | 113.04 | 147.00 | 0.51 | 3230.00 | 3600.00 |
| 7.10 | 4.63 | 1.41 | 16.17 | 155730.00 | 110.87 | 145.00 | 0.51 | 3210.00 | 3600.00 |
| 7.40 | 5.48 | 1.98 | 20.00 | 157732.00 | 144.50 | 153.05 | 0.50 | 3359.00 | 3500.00 |
| 7.30 | 5.47 | 1.95 | 19.90 | 157309.00 | 144.00 | 152.50 | 0.49 | 3341.00 | 3500.00 |
| 7.30 | 5.44 | 1.75 | 19.60 | 157168.00 | 144.00 | 151.21 | 0.48 | 3315.00 | 3500.00 |
| 6.90 | 5.41 | 1.65 | 19.60 | 157168.00 | 143.00 | 150.96 | 0.47 | 3310.00 | 3500.00 |
| 6.80 | 5.35 | 1.64 | 19.60 | 157074.00 | 142.50 | 150.55 | 0.46 | 3300.00 | 3400.00 |
| 6.80 | 5.34 | 1.55 | 19.30 | 156839.00 | 142.00 | 149.45 | 0.46 | 3295.00 | 3300.00 |
| 6.80 | 5.29 | 1.51 | 18.80 | 156792.00 | 141.50 | 149.01 | 0.45 | 3291.00 | 3800.00 |
| 6.70 | 5.25 | 1.82 | 18.70 | 156416.00 | 140.00 | 148.51 | 0.55 | 3271.00 | 3800.00 |
| | | 1.80 | 19.90 | | 134.21 | | 0.52 | | 3700.00 |
| | | 1.73 | 19.50 | | 133.69 | | 0.51 | | 3700.00 |
| | | 1.67 | 19.40 | | 132.20 | | 0.51 | | 3700.00 |
| | | 1.55 | 19.30 | | 132.08 | | 0.50 | | 3700.00 |
| | | 1.54 | 19.30 | | 131.04 | | 0.49 | | 3700.00 |
| | | 1.54 | 19.20 | | 130.32 | | 0.48 | | 3700.00 |
| | | 1.54 | 19.00 | | 130.06 | | 0.48 | | |
| | | | 18.80 | | 129.66 | | | | |

| 4A_MICP | 4A_MICP | 4A_MICP | 4A_MICP | 4A_MICP | FUS | 4A_MICP |
|---------|---------|---------|---------|---------|--------|---------|
| U | V | W | Y | Zn | Zn | Zr |
| ppm | ppm | ppm | ppm | ppm | ppm | ppm |
| 0.50 | 848.00 | 2.22 | 4.52 | 308.00 | 343.28 | 14.47 |
| 0.50 | 845.00 | 2.17 | 4.46 | 306.00 | 330.50 | 14.24 |
| 0.50 | 844.00 | 2.14 | 4.44 | 306.00 | 329.30 | 14.20 |
| 0.50 | 842.00 | 2.12 | 4.39 | 305.00 | 326.63 | 14.12 |
| 0.50 | 838.00 | 2.06 | 4.35 | 302.00 | 325.94 | 13.90 |
| 0.50 | 838.00 | 2.06 | 4.32 | 302.00 | 324.50 | 13.76 |
| 0.50 | 835.00 | 2.00 | 4.30 | 302.00 | 319.06 | 13.67 |
| 0.50 | 833.00 | 2.00 | 4.21 | 302.00 | 317.57 | 13.56 |
| 0.58 | 909.00 | 2.00 | 5.30 | 309.00 | 337.00 | 16.50 |
| 0.58 | 907.00 | 2.00 | 5.20 | 307.00 | 336.00 | 16.50 |
| 0.58 | 900.00 | 2.00 | 5.20 | 306.00 | 336.00 | 16.40 |
| 0.58 | 897.00 | 2.00 | 5.20 | 305.00 | 332.00 | 16.40 |
| 0.57 | 896.00 | 2.00 | 5.10 | 305.00 | 332.00 | 16.40 |
| 0.56 | 893.00 | 2.00 | 5.00 | 305.00 | 328.00 | 16.40 |
| 0.56 | 889.00 | 2.10 | 5.00 | 304.00 | 328.00 | 16.30 |
| 0.55 | 883.00 | 2.10 | 4.90 | 302.00 | 328.00 | 16.20 |
| | 837.21 | 2.00 | 5.03 | | | 14.81 |
| | 835.43 | 2.00 | 5.00 | | | 14.67 |
| | 835.43 | 2.00 | 4.99 | | | 14.59 |
| | 835.39 | 2.00 | 4.98 | | | 14.48 |
| | 834.72 | 2.00 | 4.97 | | | 14.41 |
| | 834.01 | 2.00 | 4.96 | | | 14.40 |
| | 833.89 | | 4.93 | | | 14.40 |
| | 832.28 | | 4.90 | | | 14.30 |

25. Reported Values

The certified values listed in this certificate fulfil the AMIS statistical criteria (see section 18) regarding agreement for certification and have been independently validated by Allan Fraser.

26. Validation of Accuracy (Trueness)

This CRM can be used to validate accuracy (trueness) as required in method validation as stated in the ISO17025 standard. See Appendix 3 for an example on the validation of accuracy using replicate data derived from the analysis of a CRM.

27. Limit of Detection and Limit of Quantitation in Gravimetric Fire Assay

In the determination of limit of detection (LOD) and limit of quantitation (LOQ) in gravimetric fire assay (*i.e.* lead collection and weighing of resulting gold prill), the minimum mass that an assay microbalance is capable of weighing and the original test sample mass determines the LOD and the LOQ in the assay (Fraser, 2015), (see Appendix 8 for an example of the calculation LOD and LOQ and Table 14 for a recommend reporting scheme for LOD and LOQ values).

28. Metrological Traceability

The values quoted herein are based on the consensus values derived from statistical analysis of the data from an inter-laboratory measurement program. Traceability to SI units is via the standards used by the individual laboratories the majority of which are accredited to the ISO17025 general requirements for the competence of testing and calibration laboratories and who have maintained measurement traceability during the analytical process.

29. Period of Validity

The certified values are valid for this product, while still sealed in its original packaging, until notification to the contrary. The stability of the material will be subject to continuous testing for the duration of the inventory. Should product stability become an issue, all customers will be notified and notification to that effect will be placed on the www.amis.co.za website.

30. Minimum Sample Size

Most of the laboratories reporting, used a 0.5g sample size for the ICP-OES and a 30g sample size for the fire assay. These are the recommended minimum sample sizes for the use of this material.

31. Availability

This product is available in Laboratory Packs containing 1kg of material and Explorer Packs containing custom weights (from 30g to 250g) of material. The Laboratory Packs are sealed bottles delivered in sealed foil pouches. The Explorer Packs contain material in standard geochem envelopes, nitrogen flushed, and vacuum sealed in foil pouches.

32. Recommended use in Quality Control

Users should set their own limits *i.e.* 1, 2 and 3 standard deviations from an obtained mean value based on at least 10 replicate analyses using this CRM (see Appendix 6 for detail on the use of this CRM in quality control).

33. Legal Notice

This certificate and the reference material described in it have been prepared with due care and attention. However, AMIS, Melesha Gopi Mungaroo; Daishnee Padayachee and Allan Fraser; accept no liability for any decisions or actions taken following the use of the reference material.

Date of Version 000: 07 July 2021

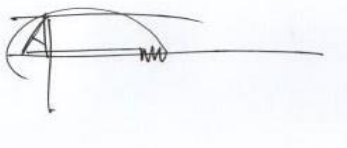
Version: 000

Approving Officer:

African Mineral Standards: _____

Melesha Gopi Mungaroo (Senior Quality Specialist)

Certifying Officer:

A handwritten signature in black ink, appearing to be 'A. Fraser', written over a horizontal line.

Geochemist: _____

Allan Fraser

M.Sc. (Geology), N.D. (Analytical Chem.),
Pr.Sci.Nat. Pr.Chem.SA

References

- Abzalov, M. (2011). Sampling Errors and Control of Assay Data Quality in Exploration and Mining Geology, Applications and Experiences of Quality Control. Ivanov, O., (Ed.), InTech, DOI: 10.5772/14965. Available from: <http://www.intechopen.com/books/applications-and-experiences-of-quality-control/sampling-errors-and-control-of-assay-data-quality-in-exploration-and-mining-geology>. Accessed 25th September, 2016.
- ASTM E122-09e1, Standard Practice for Calculating Sample Size to Estimate, With Specified Precision, the Average for a Characteristic of a Lot or Process, ASTM International, West Conshohocken, PA, 2011, www.astm.org
- Barwick V.J., Pritchard E. (Eds). (2011). Eurachem Guide; Terminology in Analytical Measurement – Introduction to VIM 3 (2011). ISBN 978-0-948926-29-7.
- Carr R.H. (2011). Estimating errors using graphs and taking good data. California State University Los Angeles. March 2011. 14. <http://web.calstatela.edu/faculty/kaniol/Err-Gph-Meas-IBooklet.pdf>. Retrieved: 29 August, 2016.
- Cochran, W.G. (1950). The Comparison of Percentages in Matched Samples. *Biometrika*, 37, 256-66.
- Ellison, S., Barwick, V., Duguid Farrant, T. (2009). *Practical statistics for the analytical scientist, a bench guide. 2nd Edition*. RSC Publishing. 25-172.
- ERM Application Note 1. (2005). Comparison of a measurement value with a certified value. European Commission – Joint Research Centre Institute for Reference Materials and Measurements (IRMM). 1-2.
- EURACHEM / CITAC Guide CG 4. (2012). Quantifying uncertainty in analytical measurement. 3rd Edition. Editors: Ellison, S. L. R., Williams, A. 4-121.
- EuroLab Technical Report No.1/2007. (2007). Measurement uncertainty revisited: Alternative approaches to uncertainty evaluation. European Federation of National Associations of Measurement, Testing and Laboratories. 43-46.
- Fraser, A.W. (2015). Minimising uncertainty in measurement and improving limit of detection in gold bearing ores from concentrations predicted by linear regression in atomic absorption spectrometry. M.Sc. thesis (unpublished), University of Johannesburg: 199 pages.
- Grubbs, F.E. (1969). Procedures for detecting outlying observations in samples. *Technometrics*, 11, 1969. 1-21.
- Horwitz, W., Albert, R. (2006). The Horwitz Ratio (HorRat): A useful index of method performance with respect to precision, *Journal of Association of Official Analytical Chemists International*, 89: 1095-1109.
- ISO/IEC 17025:2005(E).(2005). *General requirements for the competence of testing and calibration laboratories*. 2nd Edition: 14-59.
- ISO Guide 35 (2017). Certification of reference materials — General and Statistical principles, 4th edition. ISO/REMCO WG 1.
- ISO 5725-2:1994. Accuracy (trueness and precision) of measurement methods and results -- Part 2: Basic method for the determination of repeatability and reproducibility of a standard measurement method.
- Long, J., Winefordner, J. (1983). Limit of detection - a closer look at the IUPAC definition. *Analytical Chemistry* 55: 712A - 724A.
- Miller, J., Miller, J. (2010). *Statistics for analytical chemistry*. 6th Edition. New York: Ellis Horwood.36-126.
- Nelsen T.C., Wehling P. (2008). Collaborative studies for quantitative chemical analytical methods. AACC International Report. *Cereal Foods World*. September – October 2008, Vol. 53, No. 5. 285-288.
- Skoog, D., West, D. (1982). *Fundamentals of analytical chemistry*. 4th Edition. CBS College Publishing. Holt Saunders International Editions:39-73.
- Thompson, M. (Ed.) (2008). Test for 'sufficient homogeneity' in a reference material. Analytical Methods Committee, AMCTB 17A, ISSN 1757-5958.
- Thompson, M. (Ed.) (2010). Internal quality control in routine analysis. AMC Technical Brief. Analytical Methods Committee. AMCTB No.46. 2010. Issn 1757-5958.
- Thompson, M., Lowthian, P. (2011). *Notes on statistics and data quality for analytical chemists*. Imperial College Press: 15-115.

Appendices

Appendix 1: Uncertified Element Statistics

Uncertified element statistics are shown in **Table 10**.

Table 10. Uncertified element concentrations statistics (Following the removal of outliers).

| Element | Generic Method | n | Mean | SD | RSD % | Unit |
|--------------------------------|-------------------------------|----|----------|--------|-------|------|
| 4E | Pb Collection | 7 | 2.6 | 0.04 | 1 | g/t |
| Ag | 2A_MICP | 8 | 0.06 | * | * | ppm |
| Ag | 4A_MICP | 4 | 0.07 | 0.01 | 14 | ppm |
| Ag | FUS | 7 | 2.3 | 0.04 | 2 | ppm |
| Al | 2A_MICP | 15 | 42017.1 | 1646.1 | 4 | ppm |
| As | 2A_MICP | 15 | 1.7 | 0.08 | 5 | ppm |
| As | FUS | 4 | 100.0 | * | * | ppm |
| Au | 2A_MICP | 8 | 0.03 | 0.009 | 27 | g/t |
| Ba | ¹² 2A_ICPES | 8 | 49.7 | 3.0 | 6 | ppm |
| Ba | 2A_MICP | 15 | 42.9 | 2.9 | 7 | ppm |
| Ba | ¹³ 4A_ICPES | 8 | 56.4 | 1.9 | 3 | ppm |
| Ba | FS_ICPES | 16 | 61.2 | 5.8 | 9 | ppm |
| Ba | FUS | 15 | 62.3 | 2.2 | 4 | ppm |
| Ba | XRF | 11 | 72.6 | 18.4 | 25 | ppm |
| BaO | XRF | 11 | 0.008 | 0.002 | 22 | % |
| Be | 2A_ICPES | 8 | 0.3 | 0.01 | 4 | ppm |
| Be | 2A_MICP | 16 | 0.3 | 0.02 | 5 | ppm |
| Be | 4A_ICPES | 8 | 0.4 | 0.03 | 6 | ppm |
| Be | FUS | 8 | 0.5 | 0.08 | 15 | ppm |
| Bi | 2A_MICP | 16 | 0.1 | 0.01 | 13 | ppm |
| Bi | FUS | 16 | 0.2 | 0.06 | 35 | ppm |
| C | ¹⁴ Combustion/LECO | 20 | 427.3 | 132.8 | 31 | ppm |
| Ca | 2A_MICP | 16 | 29037.6 | 851.5 | 3 | ppm |
| Cd | 2A_MICP | 7 | 0.06 | 0.02 | 40 | ppm |
| Cd | 4A_MICP | 11 | 0.1 | 0.1 | 88 | ppm |
| Cd | FUS | 1 | 10.0 | * | * | ppm |
| Ce | 2A_ICPES | 8 | 5.2 | 0.4 | 8 | ppm |
| Ce | 2A_MICP | 15 | 4.2 | 0.1 | 4 | ppm |
| Ce | FUS | 8 | 6.4 | 0.1 | 2 | ppm |
| Co | 2A_ICPES | 8 | 22.7 | 1.8 | 8 | ppm |
| Co | 4A_ICPES | 8 | 146.3 | 4.7 | 3 | ppm |
| Co | FS_ICPES | 15 | 148.3 | 3.3 | 2 | ppm |
| Cr | 2A_ICPES | 8 | 1096.0 | 47.3 | 4 | ppm |
| Cr | 2A_MICP | 16 | 1216.9 | 144.2 | 12 | ppm |
| Cr | 4A_ICPES | 7 | 66130.0 | 2831.7 | 4 | ppm |
| Cr | 4A_MICP | 8 | 61765.0 | 6989.9 | 11 | ppm |
| Cr | FS_ICPES | 8 | 90841.3 | 3693.1 | 4 | ppm |
| Cr | FUS | 40 | 101201.5 | 6192.4 | 6 | ppm |
| Cr ₂ O ₃ | FUS | 40 | 14.8 | 0.9 | 6 | % |
| Cs | 2A_MICP | 7 | 0.3 | 0.008 | 3 | ppm |
| Cs | FUS | 16 | 0.3 | 0.04 | 12 | ppm |
| Cu | 2A_ICPES | 8 | 120.5 | 5.5 | 5 | ppm |
| Cu | 4A_ICPES | 8 | 117.9 | 3.9 | 3 | ppm |
| Cu | FS_ICPES | 8 | 116.5 | 4.8 | 4 | ppm |
| Cu | FUS | 32 | 145.8 | 39.6 | 27 | ppm |
| Cu | XRF | 8 | 50.0 | * | * | ppm |
| CuO | XRF | 7 | 0.006 | * | * | % |
| Dy | 4A_MICP | 8 | 0.9 | 0.02 | 2 | ppm |
| Dy | FUS | 8 | 0.9 | 0.05 | 6 | ppm |
| Er | 4A_MICP | 8 | 0.6 | 0.007 | 1 | ppm |
| Er | FUS | 8 | 0.6 | 0.02 | 3 | ppm |
| Eu | 4A_MICP | 8 | 0.3 | 0.01 | 4 | ppm |

| Element | Generic Method | n | Mean | SD | RSD % | Unit |
|-------------------------------|--------------------------|----|----------|--------|-------|---------------|
| Nb | 2A_MICP | 12 | 0.6 | 0.3 | 57 | ppm |
| Nb | FUS | 8 | 158.9 | 5.1 | 3 | ppm |
| Nd | 4A_MICP | 8 | 3.2 | 0.07 | 2 | ppm |
| Nd | FUS | 8 | 3.2 | 0.05 | 2 | ppm |
| Ni | 2A_ICPES | 8 | 287.0 | 11.0 | 4 | ppm |
| Ni | 4A_ICPES | 8 | 964.6 | 41.1 | 4 | ppm |
| Ni | FS_ICPES | 8 | 808.5 | 31.7 | 4 | ppm |
| Ni | XRF | 7 | 882.9 | 4.9 | 0.6 | ppm |
| NiO | XRF | 7 | 0.1 | 0.0005 | 0.5 | % |
| Os | NiS | 16 | 0.06 | 0.005 | 10 | ppm |
| P | 2A_MICP | 16 | 120.7 | 12.0 | 10 | ppm |
| P | FUS | 16 | 131.3 | 33.7 | 26 | ppm |
| P ₂ O ₅ | XRF | 15 | 0.02 | 0.003 | 15 | % |
| Pb | 2A_ICPES | 8 | 6.6 | 1.2 | 18 | ppm |
| Pb | 2A_MICP | 15 | 6.7 | 0.3 | 5 | ppm |
| Pb | FUS | 15 | 51.1 | 47.4 | 93 | ppm |
| Pd | ¹⁸ NiS_ICPOES | 8 | 0.7 | 0.07 | 10 | g/t |
| Pr | 4A_MICP | 8 | 0.8 | 0.009 | 1 | ppm |
| Pr | FUS | 8 | 0.8 | 0.04 | 5 | ppm |
| Pt | NiS_ICPOES | 8 | 1.3 | 0.1 | 10 | g/t |
| Rb | 2A_MICP | 15 | 3.8 | 0.2 | 5 | ppm |
| Rb | FUS | 15 | 5.6 | 0.2 | 3 | ppm |
| Re | 2A_MICP | 3 | 0.001 | * | * | ppm |
| Re | 4A_MICP | 2 | 0.002 | * | * | ppm |
| Rh | NiS_ICPOES | 8 | 0.2 | 0.02 | 11 | g/t |
| Rh | Pb Collection | 16 | 0.2 | 0.02 | 9 | g/t |
| Ru | NiS_ICPOES | 8 | 0.4 | 0.04 | 10 | g/t |
| S | 2A_MICP | 16 | 0.06 | 0.009 | 15 | % |
| S | Combustion/LECO | 24 | 0.05 | 0.009 | 17 | % |
| S | FUS | 16 | 0.08 | 0.01 | 16 | % |
| Sb | 2A_MICP | 16 | 1.0 | 0.3 | 28 | ppm |
| Sb | FUS | 16 | 1.7 | 0.2 | 9 | ppm |
| Sc | 2A_ICPES | 8 | 3.7 | 0.2 | 5 | ppm |
| Sc | 2A_MICP | 16 | 3.3 | 0.3 | 8 | ppm |
| Sc | 4A_ICPES | 8 | 17.5 | 0.7 | 4 | ppm |
| Sc | FS_ICPES | 8 | 17.1 | 0.9 | 5 | ppm |
| Sc | FUS | 16 | 18.5 | 3.6 | 19 | ppm |
| Se | 2A_MICP | 7 | 0.2 | 0.05 | 22 | ppm |
| Se | 4A_MICP | 8 | 4.0 | 1.2 | 30 | ppm |
| SG | ¹⁹ SG | 15 | 3.5 | 0.2 | 4 | Dimensionless |
| Si | FUS | 16 | 163437.5 | 3501.8 | 2 | ppm |
| Sm | 4A_MICP | 8 | 0.7 | 0.02 | 2 | ppm |
| Sm | FUS | 8 | 0.8 | 0.03 | 4 | ppm |
| Sn | 2A_MICP | 8 | 6.4 | 0.4 | 7 | ppm |
| Sn | 4A_MICP | 7 | 25.2 | 2.3 | 9 | ppm |
| Sr | 2A_ICPES | 8 | 127.0 | 4.8 | 4 | ppm |
| Sr | 2A_MICP | 16 | 119.8 | 6.6 | 5 | ppm |
| Sr | 4A_ICPES | 8 | 130.8 | 5.3 | 4 | ppm |
| Sr | FS_ICPES | 7 | 134.7 | 2.0 | 1 | ppm |
| Sr | FUS | 15 | 142.1 | 4.3 | 3 | ppm |
| Sr | XRF | 8 | 556.3 | 39.3 | 7 | ppm |
| SrO | XRF | 8 | 0.07 | 0.005 | 7 | % |

| Element | Generic Method | n | Mean | SD | RSD % | Unit |
|------------------|------------------------|----|---------|-------|-------|------|
| Eu | FUS | 8 | 0.3 | 0.01 | 3 | ppm |
| Fe | 2A_MICP | 16 | 13602.5 | 442.7 | 3 | ppm |
| Ga | 2A_MICP | 15 | 5.6 | 0.2 | 4 | ppm |
| Ga | FUS | 15 | 30.1 | 0.7 | 2 | ppm |
| Gd | 4A_MICP | 8 | 0.8 | 0.02 | 3 | ppm |
| Gd | FUS | 8 | 0.9 | 0.02 | 3 | ppm |
| Ge | 4A_MICP | 16 | 0.5 | 0.3 | 64 | ppm |
| Ge | FUS | 1 | 1.0 | * | * | ppm |
| GOI | ¹⁵ GOI | 23 | 0.2 | 0.04 | 15 | % |
| Hf | 2A_MICP | 16 | 0.08 | 0.02 | 32 | ppm |
| Hf | FUS | 8 | 1.1 | 0.05 | 4 | ppm |
| Hg | 2A_MICP | 8 | 0.02 | 0.006 | 34 | ppm |
| Ho | 4A_MICP | 8 | 0.2 | 0.005 | 2 | ppm |
| Ho | FUS | 8 | 0.2 | 0.006 | 3 | ppm |
| In | 2A_MICP | 7 | 0.008 | 0.001 | 12 | ppm |
| In | 4A_MICP | 7 | 0.05 | 0.03 | 63 | ppm |
| K | 2A_MICP | 15 | 856.9 | 53.5 | 6 | ppm |
| K | FUS | 23 | 1827.3 | 282.8 | 15 | ppm |
| K | XRF | 16 | 1238.8 | 41.0 | 3 | ppm |
| K ₂ O | FUS | 16 | 0.2 | 0.03 | 15 | % |
| K ₂ O | XRF | 16 | 0.1 | 0.005 | 3 | % |
| La | 2A_ICPES | 8 | 3.2 | 0.2 | 5 | ppm |
| La | 2A_MICP | 15 | 2.1 | 0.09 | 4 | ppm |
| La | 4A_ICPES | 8 | 1.9 | 0.2 | 10 | ppm |
| La | 4A_MICP | 22 | 2.9 | 0.06 | 2 | ppm |
| La | ¹⁶ FS_ICPES | 8 | 2.8 | 0.2 | 8 | ppm |
| La | FUS | 8 | 3.1 | 0.08 | 3 | ppm |
| Li | 2A_ICPES | 8 | 8.7 | 0.7 | 8 | ppm |
| Li | 2A_MICP | 16 | 4.8 | 0.4 | 8 | ppm |
| Li | 4A_ICPES | 8 | 5.7 | 0.8 | 13 | ppm |
| Li | FUS | 16 | 8.0 | 2.1 | 26 | ppm |
| Lu | 4A_MICP | 8 | 0.09 | * | * | ppm |
| Mg | 2A_MICP | 16 | 12978.0 | 565.1 | 4 | ppm |
| Mn | 2A_ICPES | 8 | 322.3 | 14.6 | 5 | ppm |
| Mn | 2A_MICP | 15 | 296.7 | 8.5 | 3 | ppm |
| Mn | 4A_ICPES | 8 | 1591.7 | 66.4 | 4 | ppm |
| Mn | FS_ICPES | 8 | 1428.5 | 50.4 | 4 | ppm |
| Mn | XRF | 16 | 1615.5 | 189.5 | 12 | ppm |
| MnO | XRF | 16 | 0.2 | 0.02 | 12 | % |
| Mo | 2A_MICP | 7 | 0.2 | 0.02 | 9 | ppm |
| Mo | 4A_MICP | 15 | 0.4 | 0.05 | 13 | ppm |
| Moisture | ¹⁷ Moisture | 15 | 0.2 | 0.04 | 24 | % |
| Na | 2A_MICP | 16 | 4823.3 | 239.4 | 5 | ppm |

| Element | Generic Method | n | Mean | SD | RSD % | Unit |
|-------------------------------|----------------|----|--------|-------|-------|------|
| Ta | 2A_MICP | 16 | 4.2 | 4.2 | 101 | ppm |
| Ta | FUS | 16 | 179.2 | 6.2 | 3 | ppm |
| Tb | 4A_MICP | 7 | 0.1 | * | * | ppm |
| Tb | FUS | 8 | 0.2 | 0.005 | 3 | ppm |
| Te | 2A_MICP | 8 | 0.04 | 0.009 | 26 | ppm |
| Te | 4A_MICP | 8 | 0.09 | 0.02 | 19 | ppm |
| Th | 2A_ICPES | 8 | 1.0 | 0.2 | 23 | ppm |
| Th | 2A_MICP | 16 | 0.4 | 0.02 | 6 | ppm |
| Th | FUS | 15 | 0.5 | 0.01 | 3 | ppm |
| Ti | 2A_MICP | 16 | 741.5 | 29.8 | 4 | ppm |
| Ti | XRF | 21 | 3580.0 | 22.8 | 0.6 | ppm |
| TiO ₂ | XRF | 21 | 0.6 | 0.004 | 0.6 | % |
| Tl | 2A_MICP | 8 | 0.09 | 0.008 | 9 | ppm |
| Tl | 4A_MICP | 15 | 0.1 | 0.005 | 5 | ppm |
| Tm | 4A_MICP | 8 | 0.09 | 0.005 | 5 | ppm |
| Tm | FUS | 8 | 0.09 | 0.005 | 6 | ppm |
| U | 2A_MICP | 14 | 0.3 | 0.01 | 4 | ppm |
| U | FUS | 8 | 0.6 | 0.02 | 4 | ppm |
| V | 2A_ICPES | 8 | 61.2 | 3.3 | 5 | ppm |
| V | 2A_MICP | 16 | 48.7 | 1.7 | 4 | ppm |
| V | 4A_ICPES | 8 | 768.0 | 86.3 | 11 | ppm |
| V | FS_ICPES | 8 | 839.9 | 33.5 | 4 | ppm |
| V | FUS | 15 | 914.7 | 18.6 | 2 | ppm |
| V | XRF | 16 | 908.4 | 12.4 | 1 | ppm |
| V ₂ O ₅ | XRF | 16 | 0.2 | 0.002 | 1 | % |
| W | 2A_MICP | 16 | 0.1 | 0.03 | 19 | ppm |
| W | FUS | 8 | 2.6 | 0.2 | 6 | ppm |
| Y | 2A_ICPES | 8 | 2.3 | 0.3 | 15 | ppm |
| Y | 2A_MICP | 15 | 1.7 | 0.07 | 4 | ppm |
| Y | 4A_ICPES | 8 | 4.4 | 0.2 | 4 | ppm |
| Y | FS_ICPES | 8 | 4.9 | 0.2 | 5 | ppm |
| Y | FUS | 7 | 5.4 | 0.06 | 1 | ppm |
| Yb | 4A_MICP | 8 | 0.6 | 0.02 | 3 | ppm |
| Yb | FUS | 8 | 0.6 | 0.02 | 4 | ppm |
| Zn | 2A_ICPES | 8 | 23.9 | 2.5 | 11 | ppm |
| Zn | 2A_MICP | 16 | 18.9 | 1.0 | 5 | ppm |
| Zn | 4A_ICPES | 8 | 294.8 | 12.4 | 4 | ppm |
| Zn | FS_ICPES | 8 | 288.5 | 14.4 | 5 | ppm |
| Zn | XRF | 7 | 300.0 | * | * | ppm |
| ZnO | XRF | 7 | 0.04 | * | * | % |
| Zr | 2A_ICPES | 8 | 2.8 | 0.07 | 3 | ppm |
| Zr | 2A_MICP | 16 | 2.4 | 0.6 | 25 | ppm |
| Zr | FUS | 8 | 23.6 | 0.7 | 3 | ppm |

* denotes that the results were the same and SD and RSD% could not be calculated

12. 2A_ICPES is Aqua Regia Digest, with ICP-OES finish.
13. 4A_ICPES is 4 Acid Digestion (HNO₃-HClO₄-HF-HCl), ICP-OES finish.
14. Combustion/ LECO
15. GOI is gain on ignition.
16. FS_ICPES is Fusion, ICP-OES finish.
17. Moisture is residual moisture content.
18. NiS_ICPOES is Nickel Sulphide with ICPOES finish
19. SG is specific gravity.

Appendix 2 through 9, prepared by Allan Fraser.

Appendix 2. Certification of Reference Material and Estimation of Measurement Uncertainty

In the establishment of a consensus value for the CRM, outlier tests are carried out followed by performance statistics and the estimation of the measurement uncertainty. In practice, it is highly likely that data generated by multiple laboratories as an inter-laboratory comparison of material for certification, will contain erroneous as well as extreme measurements (outliers). The influence of outliers on summary statistics needs to be minimised by the application of procedures for outlier identification on raw data. The use of z-scoring, Cochran's test for suspect repeatability variances, along with Grubbs test for suspect measurement values allows for the detection of outliers (IUPAC, 2006). Method performance in terms of precision as relative standard deviation is judged by the application of the Horwitz ratio, which gives an indication of whether the observed relative standard deviation at the concentration levels of analyte determined are acceptable (Horwitz & Albert, 2006).

In the absence of an extensive uncertainty budget, measurement uncertainty is estimated from the reproducibility standard deviation from inter-laboratory data and reported as an expanded uncertainty at a level of confidence of 95% (Miller & Miller, 2010).

The steps below give detail on the establishment of a consensus value through the elimination of outliers, method performance and estimation of measurement uncertainty using standard uncertainties and the analysis of variance.

Z-Score

A z-score is calculated using equation [1]:

$$z = \frac{x - x_a}{s_p} \quad [1]$$

Where, x is the result of a submitted sample, x_a is the mean and s_p is the standard deviation of the submitted results from all the participating laboratories. Z-Scores are interpreted as follows:

$|z| \leq 2$ satisfactory performance
 $2 < |z| \leq 3$ questionable performance
 $|z| > 3$ unsatisfactory performance

(Thompson & Lowthian, 2011)

Data with z-scores exceeding two are discarded and are not included for further assessment.

Cochran's Test

The test of Cochran (1950) as shown in equation [2] is applied to any suspect repeatability variances:

$$C_{calc} = \frac{s_{max}^2}{\sum_{i=1}^l s_i^2} \quad [2]$$

Where, C_{calc} , s_{max}^2 and $\sum_{i=1}^l s_i^2$, are the calculated values for Cochran's test, data set with the maximum variance and the sum of the variances of all of the participating, l laboratory datasets. The C_{calc} value is compared with a critical value, C_{crit} at a level of confidence of 95% and an alpha of 0.05% (see Ellison, *et al.*, 2009, Appendix A, Table A.3a, page 209 for a table of critical values for the test of Cochran at LOC 95%).

According to ISO 5725-2 (1999), results from a laboratory with a suspect repeatability variance can be excluded if it is shown by the Cochran test to be an outlier. Therefore, if $C_{calc} > C_{crit}$, the laboratory with the maximum variance is removed. The data found to be excluded should not be $>2/9$, or 22% of the total data.

Grubbs Test

The test of Grubbs (1969) calculates a test statistic, G_{calc} and in the detection of a single outlier, G_1 is found by using

$$G_{1\,calc} = \frac{|Suspect\ value - \bar{x}|}{s} \quad [3]$$

Where, the sample mean and standard deviation, \bar{x} and s , are calculated with the suspect value included. The $G_{1\,calc}$ statistic is compared to a critical value for N measurements. See Ellison, *et al.*, 2009, Appendix A, Table A.2, page 208 for a table of critical values for the test of Cochran at LOC 95%.

Method Performance

The Horwitz function is used to assess the performance of the data under consideration, with respect to precision (Horwitz & Albert, 2006). A calculated %RSD is found using the Horwitz expression

$$\%RSD = \pm 2^{(1-0.5\log C)} \quad [4]$$

where, C is the analyte concentration in percent divided by 100 and \log is the natural logarithm. The observed %RSD is calculated as

$$Observed\ \%RSD = \frac{s}{Mean} \times 100 \quad [5]$$

where s is the standard deviation of n replicates.

The ratio of the observed %RSD and the calculated %RSD gives the Horwitz ratio (HorRat):

$$HorRat = \frac{\%RSD\ Observed}{\%RSD\ Calculated} \quad [6]$$

A HorRat <2 indicates that the method is of adequate precision. Should the HorRat be >2 the overall data are discarded, and the candidate material considered not suitable for certification as the precision is excessive for the concentration of the analyte being determined (Nelsen & Wehling, 2008).

Grand Mean

The grand mean ($\bar{\bar{x}}$) *i.e.* the certified value of a dataset is the total of all the data values divided by the total sample size (n):

$$\bar{\bar{x}} = \sum \frac{x}{n} \quad [7]$$

Certified Value

From ANOVA as per the description in section 18, an 'appropriate precision' as shown in [8] is calculated for sufficient homogeneity (Thompson, 2008):

$$s_r \leq 0.3u_c \quad [8]$$

Where, s_r is the within laboratory repeatability, as determined from [14]. Once [8] is satisfied, a grand mean [7] is calculated and this is taken to be the certified value.

Total Variation (SST)

The total variation (not the variance) comprises the sum of the squares of the differences of each mean with the grand mean.

$$SST = \sum (x - \bar{x})^2 \quad [9]$$

Between Group Variation (SSB)

The *variation* due to the between the laboratories is denoted SSB or Sum of Squares Between laboratories and given by [10]. If the laboratory means are close to each other (and therefore the Grand Mean) SSB will be a small value. There are P samples involved with one datum value for each sample (the sample mean), so there are P-1 degrees of freedom.

$$SSB = \sum n(\bar{x} - \bar{\bar{x}})^2 \quad [10]$$

The *variance* due to the interaction between the laboratories is denoted MSB for Mean Square Between groups and is the SSB divided by its degrees of freedom.

$$MS = \frac{SSB}{n - 1} \quad [11]$$

Within Group Variation (SSW)

The variation due to differences within individual samples is denoted SSW for Sum of Squares Within laboratories. The degrees of freedom are equal to the sum of the individual degrees of freedom for each sample. Since each sample has degrees of freedom (*df*) equal to one less than their sample sizes, and there are *k* samples, the total degrees of freedom is P less than the total sample size: $df = n - P$.

$$SSW = \sum df \cdot s^2 \quad [12]$$

The variance due to the differences within individual samples is denoted MSW for Mean Square Within groups. This is the within group variation divided by its degrees of freedom:

$$MSW = \frac{SSW}{P - n} \quad [13]$$

From equations [9] through [13], the ANOVA table as shown in Table 11 is developed.

Table 11. A single-factor ANOVA table showing key elements. Where P is the total number of groups, or laboratories. P-1 is 1 less than number of laboratories, P (n-1) is the number of data values minus number of groups (equals degrees of freedom for each group added together), and P-1 + P(n-1) is 1 less than the number of data points. MS is the mean squares of between laboratories and within laboratories. After Ellison *et al.*, (2009), Table 6.2, page 61.

| Source | Sum of Squares | df | Mean Sum of Squares | F | p | F _{crit} |
|----------------------|----------------|--------------|---------------------|---------|-----------------|-------------------|
| Between Laboratories | SSB | P-1 | MSB=SSB/df | MSB/MSW | =FDIST(x,df,df) | F-table |
| Within Laboratories | SSW | P(n-1) | MSW=SSW/df | – | – | – |
| Total | SSB+SSW | P-1 + P(n-1) | – | – | – | – |

Combined Standard Uncertainty

The combined standard uncertainty (u_c) represents the effects of random events such as days, instruments, and analysts on the precision of the analytical procedures of all accepted data of the participating laboratories. Using the output from ANOVA, the combined standard uncertainty (u_c) is determined from the square root of the sum of squares of the variances of the within laboratory repeatability, s_r , and the between laboratory precision, s_s :

$$u_c = \sqrt{s_r^2 + s_s^2} \quad [14]$$

Within laboratory repeatability is determined as

$$s_r = \sqrt{MSB} \quad [15]$$

and, the between laboratory precision as

$$s_s = \sqrt{\frac{(MSW - MSB)}{n}} \quad [16]$$

where MSW is the mean squares of the within laboratory variance, MSB is the mean squares for the between laboratories and n in this case, is the number of replicates in a group of the accepted data (Thompson & Lowthian, 2011).

Expanded Uncertainty

The expanded uncertainty (U) at a confidence level of 95% is determined by multiplication of the combined uncertainty (u_c) by a coverage factor (k) found from $N-1$ degrees of freedom (df), where N is the number of laboratory means accepted in the establishment of the certified value. The t-critical value for 5% significance can be found in a t-critical table (see Appendix 9, or from MS Excel as =TINV (5%, df)).

Uncertainty Statement

Typically, an uncertainty statement is presented as follows: Au =0.77±0.04 g/t, where the number following the symbol ± is the numerical value of an expanded uncertainty, $U = ku_c$, with U determined from a combined standard uncertainty multiplied by a coverage factor $k = 2$ or, a t-critical value for $N-1$ accepted laboratories. Since it can be assumed that the possible estimated values of the standard are approximately normally distributed with standard uncertainty, u_c , the certified value of the CRM is believed to lie in the interval defined by U with a level of confidence of approximately 95 %, e.g. a mean value of 0.77±0.04g/t will have intervals of: 0.73≤0.77≤0.81 g/t.

Appendix 3. Example: Comparison of Mean and Certified Value for Validation of Accuracy

According to ERM (2005); Eurolab (2007); Abzalov (2011) and Carr (2011), the validation of accuracy for a given mean and certified value requires the inclusion of the measurement uncertainty of the CRM in a t-test for statistical significance. The classical Student's t-test as shown in [17], does not consider the measurement uncertainty of the CRM. To compensate for this, Eurolab Technical Report No.1/2007 recommends equation [18] for the validation of CRMs with stated measurement uncertainties.

$$t_{calc} = \frac{|\bar{x} - \mu|}{\frac{s}{\sqrt{n}}} \quad [17]$$

$$t_{calc} = \frac{|\bar{x} - \mu|}{\sqrt{(u_{\mu})^2 + \frac{s^2}{n}}} \quad [18]$$

Where, t_{calc} is the calculated t-statistic, \bar{x} the mean of n replicates with a standard deviation of s for a CRM of μ certified value. The standard uncertainty u is the stated expanded uncertainty (U) of the CRM divided by the coverage factor (k) as stated on the certificate of analysis. Note that the $| \quad |$ bars indicate that the absolute value between the mean and the certified value is to be used, *i.e.* ignore the sign.

An example in which [18] is used for validation of accuracy is given below.

Example

A CRM is independently replicated nine times for Al_2O_3 concentration by XRF analysis, *i.e.* 9 individual fused glass beads were prepared. The observed mean and standard deviation of the replicate data are shown with the certified value and expanded uncertainty in Table 12. In validation of accuracy, the hypothesis question is: Is the difference between the observed mean and the certified value statistically significant at a level of confidence of 95%? Alternatively put, is there sufficient evidence to conclude that the data *i.e.* replicates generated, are inaccurate?

The relevant hypotheses are:

Null hypothesis: H_0 : Mean = Certified value of CRM with stated measurement uncertainty. The acceptance of H_0 means that accuracy is demonstrated; *i.e.* insufficient evidence to reject H_0 ;

Alternate hypothesis: H_1 : Mean \neq Certified value of CRM with stated measurement uncertainty. The acceptance of H_1 means that accuracy is not demonstrated, *i.e.* there is sufficient evidence to accept H_1 ;

Table 12. CRM certified value, quoted expanded uncertainty U , the coverage factor for the CRM, $k=2.25$ and mean for $n=9$ replicates and corresponding standard deviation for the replicate data.

| CRM Certified Value | Expanded % (U) | Coverage Factor (k) | Mean ($n=9$) | n | Standard Deviation (s) |
|---------------------|--------------------|-------------------------|----------------|-----|----------------------------|
| 4.62% | 0.08% | 2.25 | 4.59% | 9 | 0.01015 |

The standard uncertainty (u) is found by dividing the expanded uncertainty by the coverage factor:

$$u = \frac{0.08}{2.25} = 0.0356 \%$$

Using the observed mean for the replicate data ($n=9$) obtained for the CRM and substituting into [18]:

$$t_{calc} = \frac{|\bar{x} - \mu|}{\sqrt{0.0356^2 + \frac{0.01015^2}{9}}} = \frac{|4.59 - 4.62|}{\sqrt{0.00126 + 0.00001145}} = 0.84$$

Therefore, $t_{calc} = 0.84$ and $t_{crit}(5\%, 8) = 2.31$ (df is 8, therefore, $t_{crit}=2.31$, see Appendix 9, page 36) which is >0.84 . Similarly, the p -value= 0.43 which is >0.05 . This is strong evidence in favour of accepting the null hypothesis that there is no significant statistical difference between the certified value and the observed mean. Therefore, under the conditions that the uncertainty associated with the certified value is known, the accuracy is validated for the CRM tested. If the null hypothesis is accepted that the mean obtained is not statistically different from the certified value, then the principle of traceability has been conformed to.

Appendix 4. Two-standard Deviations

Two-standard deviations are calculated using the expression:

$$\text{Two standard deviations} = 2 (u_c) \quad [19]$$

Where, u_c is the standard combined uncertainty (equation [14]).

Appendix 5. Confidence Interval

Confidence interval is calculated as:

$$\text{Confidence Interval (CI)} = \frac{(t_{critical})s}{\sqrt{N}} \quad [20]$$

Where, N is the number of laboratories (accepted laboratory data), $t_{critical}$ is a two-tailed value for $N - 1$ degrees of freedom (df) and s , is the standard deviation of the accepted laboratory means. A two-tailed critical value is found for $N - 1$ degrees of freedom from either a t -distribution table (Appendix 9) or MS Excel as =TINV (5%, df).

Appendix 6. Using the CRM in Quality Control

QC chart control limits should not be determined by the certified value and stated measurement uncertainty of the certified reference material used. These parameters although “certified” will never be known; it is only the corresponding statistical estimates, *i.e.* standard deviation and the mean calculated from replicated results that are known and these should be used in quality control charts. However, should the laboratory choose to use the certified value as the mean then the quoted $2s$, or CI value for the CRM can be used in the quality control chart.

It is recommended that a Shewhart chart be developed for the use if this CRM is to be used as a control sample in laboratory quality control. A Shewhart chart is a plot of sequential assay results obtained from quality control material such as an AMIS CRM. The warning and control limits are based on the standard deviation obtained from the mean of the replicates of a CRM (Ellison, *et al.*, 2009; Thompson, 2010). The procedure in preparing a Shewhart chart is as follows:

1. Analyse 10 to 15 replicates or more of the AMIS CRM.
2. Apply the Grubbs test for outliers.
3. Determine the mean of the replicates after application of the Grubbs test.
4. Determine the standard deviation, using equation [21], of the replicates following Grubbs test.
5. Calculate the standard deviation, s from:

$$s = \sqrt{\frac{\sum(x_i - \bar{x})^2}{n - 1}} \quad [21]$$

where, x_i is an individual measurement in the data set, \bar{x} is the mean of the data set at $n-1$ degrees of freedom (df) and n is the number of replicates. The sample standard deviation can be found using the MS Excel formula “=stdev.s (number1;)”.

6. Verify accuracy of the mean value using equation [18].
7. Once accuracy is verified, calculate $\pm 2s$ and $\pm 3s$, where s is the standard deviation calculated from [21].
8. Construct the Shewhart control chart around the mean of n replicates.
9. Use $\pm 2s$ as the warning limits.
10. Use $\pm 3s$ as the control limits.
11. It is recommended that if 2 to 3 points are outside the warning limits analyse another sample and if it is then within warning limits, continue. If it is outside the warning limits, stop and troubleshoot.
12. It is recommended that if any point is outside control limits, analyse another portion (sample) of the CRM. If it is within control limits, continue. If it is outside control limits, stop and troubleshoot.
13. For reference purposes, the CRM certified value can be plotted on the Shewhart chart alongside the mean value.

On a regular basis the accuracy of the replicates of the CRM should be assessed in terms of the certified value of the CRM using equation [18].

Appendix 7. Conversion to Air-dry Basis (Prepared by Allan Fraser)

Since AMIS certified analyte values are reported on a dry-basis, the user laboratory is required to dry a portion (accurately weigh out 1.0 grams in duplicate) of the CRM material in air at 105°C in a drying oven to constant mass to determine the moisture content. Use a crucible with a flat inner surface with a surface area not smaller than 10 cm² with the CRM material spread evenly over same; this represents a 0.1 gram spread per cm². In correcting the certified value for moisture content, a moisture correction factor is calculated:

$$\text{Moisture correction factor (MCF)} = \frac{100 - \% \text{Moisture at } 105^{\circ}C}{100} \quad [22]$$

$$\text{Air dry basis concentration} = \text{MCF} \times \text{certified value on a dry basis} \quad [23]$$

Example

The moisture content determined at 105°C on a CRM is 0.500%. The certified analyte concentration for the CRM is 12.62±0.52% (dry basis). Calculating the moisture correction factor using [22] gives:

$$\text{Moisture correction factor} = \frac{100 - 0.500}{100} = 0.995$$

Multiplying the factor of 0.995 by the certified value as stated on the certificate of analysis on a dry basis (as in [23]) gives the analyte concentration on an air-dry basis:

$$0.995 \times 12.62\% = 12.56\%$$

The stated measurement uncertainty also needs to be corrected using [22] and [23], e.g. $0.995 \times 0.52 = 0.51_{(7)}$, rounded to 0.52%. The air-dry basis concentration *i.e.* $12.56 \pm 0.52\%$ is to be used as the certified value with its corresponding measurement of uncertainty.

Appendix 8. Example of Determination of LOD and LOQ in Fire Assay

The limit of detection (LOD) is the minimum detectable quantity of the analyte of interest (Skoog & West, 1985). To determine the LOD in fire assay by lead collection, the minimum mass that an assay microbalance is capable of weighing (m in micrograms, and the original test sample mass, $Mass_{assay}$ in grams) determines the LOD. The smallest prill mass most assay microbalances can measure is $1\mu\text{g}$ or 0.001mg . Even with a microscope it may be difficult to locate and pick up a prill weighing ten times that amount (*i.e.* 0.01mg or $10\mu\text{g}$) and weigh it. If an analyst can weigh a prill of $1\mu\text{g}$ then the LOD becomes $1\mu\text{g}$. However, the concentration factor would be 50 times for a 50-gram assay sample and therefore the LOD in g/t becomes $1\mu\text{g}$ divided by the original mass of the sample in grams taken for fire assay [24]. Therefore, the LOD in fire assay is computed as:

$$LOD = \frac{m (\mu\text{g})}{Mass_{assay} (g)} (\text{g/t}) \quad [24]$$

The limit of quantitation (LOQ), is simply the LOD multiplied by 10 (Long & Winefordner, 1983):

$$LOQ = 10 \cdot \frac{m (\mu\text{g})}{Mass_{assay} (g)} (\text{g/t}) \quad [25]$$

Therefore, with a sample mass of 50g taken for fire assay, the limit of detection would be 0.02g/t . *i.e.* $1\mu\text{g} = 1\text{g/t}$, therefore $1\mu\text{g}/50\text{g} = 0.02\text{g/t}$. If no prill was found, then the LOD result would be $<0.02\text{g/t}$ or "not detected". Using a larger assay sample mass improves the LOD and LOQ (Table 12). Table 14 gives a recommended reporting scheme for LOD and LOQ.

Table 13. Mass of assay sample and corresponding limit of detection and limit of quantitation for an assay microbalance capability of smallest prill mass of $1\mu\text{g}$ or 0.001mg .

| Mass Assay Sample (g) | LOD (g/t) | LOQ (g/t) |
|-----------------------|-----------|-----------|
| 30 | 0.03 | 0.3 |
| 50 | 0.02 | 0.2 |
| 100 | 0.01 | 0.1 |

Table 14. Recommended reporting scheme for LOD and LOQ in fire assay.

| Data | Report as |
|------------|---------------------|
| <LOD | Not detected |
| <LOQ | Detected |
| \geq LOQ | Report assay result |

Appendix 9. T-distribution table

Table 15. T-distribution table for t-critical values (t crit.) for a two-tailed t-test at a 95% level of confidence.

| df | Two-tailed | df | Two-tailed |
|-----------|-------------------|-----------|-------------------|
| 1 | 12.71 | 23 | 2.06 |
| 2 | 4.30 | 24 | 2.06 |
| 3 | 3.18 | 25 | 2.06 |
| 4 | 2.78 | 26 | 2.05 |
| 5 | 2.57 | 27 | 2.05 |
| 6 | 2.44 | 28 | 2.04 |
| 7 | 2.36 | 29 | 2.04 |
| 8 | 2.30 | 30 | 2.04 |
| 9 | 2.26 | 35 | 2.03 |
| 10 | 2.22 | 40 | 2.02 |
| 11 | 2.20 | 45 | 2.01 |
| 12 | 2.17 | 50 | 2.00 |
| 13 | 2.16 | 55 | 2.00 |
| 14 | 2.14 | 60 | 2.00 |
| 15 | 2.13 | 70 | 1.99 |
| 16 | 2.12 | 80 | 1.98 |
| 17 | 2.11 | 90 | 1.98 |
| 18 | 2.10 | 100 | 1.98 |
| 19 | 2.09 | 120 | 1.98 |
| 20 | 2.08 | Infinity | 1.96 |
| 21 | 2.08 | | |
| 22 | 2.07 | | |

End of certificate